


ERLAB'S LIST OF
APPROVED CHEMICALS

GreenFumeHood® 3



TO PROTECT

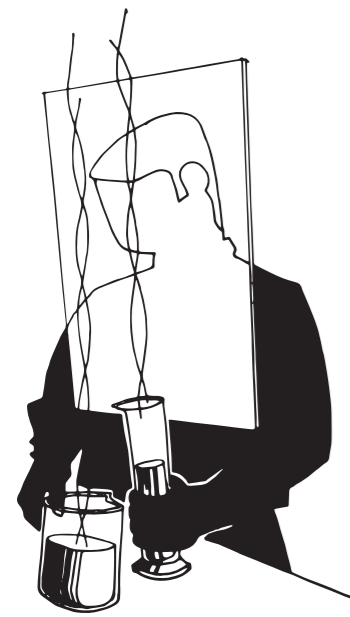

Le 11 avril 1983

Erlab is the inventor and world leader in laboratory-grade chemical filtration products.

This edition of the **Chemical listing booklet** has been developed by **Erlab's R&D laboratory**. It is the **result of 50 years of research and development into filtration technologies** and demonstrates the expertise of Erlab's R&D laboratory in the field of molecular and particulate filtration.

In compliance with **AFNOR NFX 15-211: 2009** Standard, this booklet is supplied with every Erlab ductless filtering fume hood and includes a full list of chemicals certified by Erlab for handling uses.

You may contact Erlab at any time:

- for information regarding chemicals not listed in this booklet
- to ensure you have the latest copy of this Chemical listing
- if you require information related to the handling of your chemicals

BEFORE YOU START

Caution	4
AFNOR NFX 15-211: 2009 Standard	5
Erlab molecular filters	7
Erlab filter performance tests	8
Efficiency test certificate	10
The Erlab Difference	11
List of Chemicals not recommended	12
Definitions of the column headings	13

CHEMICAL AGENTS

By alphabetical order	14
By formula	94
By CAS number	106



Caution

The retention capacities given in this Chemical listing only apply to Neutrodine filtration technology, manufactured and marketed by Erlab.

Occupational Exposure Limits (OEL) or Threshold Limit Values (TLV) are specific to each chemical and under no circumstances may be applied to chemical agents when used in a combination.

Therefore, the values given in this Chemical listing may alter in accordance with the latest findings regarding the chemicals included in it:

- OEL or TLV values may be reviewed and corrected by government bodies responsible for establishing them.
- Regulations regarding these chemicals (classification, storage, etc.) are updated in accordance with research and health monitoring authorities.
- Improvement in the performance of carbon filters manufactured by Erlab has a direct effect on retention capacities given in this Chemical listing.

The LP - Pay Back Survey has been designed to generate data to accurately forecast the life-cycle of the filtration media.

Based on this calculated life-cycle, Erlab forecast the savings of using Neutrodine Technology versus the energy cost associated with the classic approach of ducting to atmosphere.

Handling CMR Agents (Carcinogenic Mutagenic Reprotoxic)

Strict regulations to protect people exposed to CMR agents at their workstation are laid out by the French Labor Code. These regulations only apply to Category 1A and 1B CMR agents, as Category 2 CMR agents are not governed by these regulations.

Safety regulations are extremely demanding but at the same time pragmatic, in that they suggest alternative solutions allowing in particular the use of filtering hoods, on condition that they reduce exposure to the product to the lowest level possible.

Labor safety code requires that the following measures be taken if CMR agents are present (categories 1A or 1B):

- Regular risk assessment represented by each CMR agent.
- If it is not possible to replace a CMR agent with a non-CMR agent, the CMR agent must be handled in a closed system.

- If it is impossible to use a closed system, exposure to CMR agents must be reduced to the lowest possible levels.
- Work must be stopped when exposure to a CMR agent reaches the occupational exposure limit or TLV.

AFNOR NFX 15-211 sets minimum rules to ensure that a fume hood used by an operator performs at a high level (release into the room at a concentration less than 1% of the occupational exposure limit or TLV).

To avoid interfering with the current laws and regulations, AFNOR NFX 15-211 does not apply to CMR agents Category 1A and 1B. Category 2 CMR agents are governed by AFNOR NFX 15-211.

The purpose of the Labor Code and AFNOR NFX 15-211 is to reduce exposure to the lowest possible level.

AFNOR NFX 15-211: 2009 Standard

The AFNOR NFX 15-211: 2009 Standard was established by the Union de Normalisation de la Mécanique (UNM), composed of a team of experts (INRS, government bodies and professional unions), mandated by AFNOR. This standard applies to filtering fume hoods (also known as recirculating fume cupboards or ETRAF) designed for research work, analysis, teaching, etc. for all laboratories in which chemicals subject to occupational exposure limits (OEL or TLV-TMA) are handled. This text requires performance criteria relating to:

- **Filtration efficiency**
- **Containment efficiency**
- **Air face velocity**

Therefore, a list of approved chemicals for our filters and **a specific user manual must be provided** with filtering fume hoods.

The classes established by the standard

Class 1	Class 2 ⁽¹⁾
Filtering fume hood with safety reserve	Filtering fume hood without safety reserve
A main filtering level and a safety filtering level	One level of filtration

Classification according to the type of filtration

	Filtration type according to AFNOR NFX 15-211:2009	Equivalent Erlab filtration type
Particle filtration*	Type P	HEPA
Vapor filtration**	Type V	neutrodine Unisorb
Particle and vapor filtration**	Type PV	HEPA - neutrodine Unisorb

⁽¹⁾: GFH only available with class 1

*: the particle filter must be at least type H14 in accordance with standard NF EN 1822-1

** : vapor filters must undergo two successive tests using cyclohexane and isopropanol for filters designed to capture Volatile Organic Compounds (VOC). Another test designed for acid vapors is performed with hydrochloric acid.

Filtration efficiency (see description of the test method on page 8)

Filtration efficiency is defined by the capacity of the filter to capture dangerous molecules handled inside the enclosure and determines the quality of air filtered downstream of the filter.

	Class 1	Class 2 ⁽¹⁾
Normal operating phase	Normal operating phase during which the concentration downstream of the filters must be less than 1% of the OEL	
Detection phase	Detection phase during which the concentration downstream of the filters must be less than 1% of the OEL and failing filtration indication	Detection phase during which the concentration downstream of the filters must be less than 50% of the OEL
Safety phase	Safety phase during which the concentration downstream of the filters must be less than 50% of the OEL, and that must not be less than 1/12 of the duration of the normal operating phase	

⁽¹⁾: GFH only available with class 1

For class 2, the quantity of chemicals handled in the hood cannot exceed 1/8 of the filter retention capacity of chemicals.

Enclosure containment efficiency

Containment efficiency is determined by the ability of the hood to keep the fumes or particles inside the enclosure preventing their release into the laboratory environment.

To prove this efficiency, a test is performed according to the protocol described in the AFNOR NFX 15-211.

Test procedure: SF tracer gas (sulfur hexafluoride) is released into the enclosure and a grid holding sensors is placed opposite the handling ports. Samples are then taken on the grid.

On the basis of gas concentrations released and samples taken that will allow the average exposure of an operator to this tracer gas to be calculated, it is possible to establish a level of containment performance for the fume hood. The containment threshold set by standard AFNOR NFX 15-211: 2009 imposes a maximum SF concentration level of 0.1 ppm at the measuring points on the grid.

Air face velocity

Air face velocity is described as the ability of the hood to create a dynamic barrier between the operator and the handling.

For fume hoods with fixed front panels, the air face velocity at all of the openings must be between 0.4 and 0.6 m/s. Therefore, they must be equipped with a device to continuously monitor ventilation which is also an indicator of good containment.

Erlab Molecular filters

Quality design

Erlab guarantees that users of our filtering fume hoods and storage cabinets are working with a high quality product composed of the following:

- Each filter is designed so that the density of the carbon within it remains constant over time (US patent 4946480).
- Total control of negative wall effects (US patent 4946480).
- A stable, uniform adsorption column.
- Each filter is packaged in a sealed plastic bag identified by a serial number and manufacture date to guarantee it is in perfect working condition upon delivery.

Compliance with filtration standards

The design of our filters guarantees total protection. Inspired in part by military-type gas mask technology, our filters have been subject to rigorous testing and meet all safety requirements set forth by the following filtration standards:

- ASTM (American Standard Test Method) standard: This standard pertains to carbon as a raw material used in the design of filtration cartridges. It is an evaluation standard that pertains only to the quality of the raw material. It thus allows Erlab to select high-performance carbons.
- AFNOR (French standardization organization) standard NF X 15-211: This standard guarantees the filtration performance of the filters used in our units. It sets forth requirements regarding the air quality downstream from the filter. This level of quality is proven through tests carried out by an independent laboratory, these tests demonstrate the filtration quality of our filters.

Types of filter and optional prefilters available for GFH

The type of filter recommended is specific to the type of chemicals being handled and allows users to benefit from high retention capacities for the following:

Filter / Optional prefilters	For Use With
Neutrodine Unisorb	Universal filter gases and vapors
HP	HEPA H14 for powders
Particular prefilter	Dust
Optional ammonia prefilter	Ammonia vapors + dust
Optional acid prefilter	Heated acid vapors + dust

Neutrodine Unisorb Filter performance tests procedure

The Neutrodine Unisorb filters are subject to performance tests conducted in accordance to the requirements of the AFNOR NFX 15-211: 2009 Standard. The results of the tests given in this list of approved chemicals demonstrates the technological performance developed by Erlab.

The fume hood that the tests were performed in was fitted with new filters and installed in a closed space. The chemical used for the test was evaporated in the fume hood to give a constant concentration during all operating phases stated in the standard.

The three chemicals selected for the Erlab filter performance tests were:

- Isopropanol
- Cyclohexane
- Hydrochloric acid

The concentration of the chemical downstream of the filtration system was checked at least three times per hour during all the filtering fume hood operating phases and was expressed in ppm by volume.

The maximum values of the reference chemicals are given in the list of approved chemicals provided with every Erlab filtered fume hood.

The test was performed in 8 hour sequences, 16 hours apart.

The Analyzers

Whatever the chemical being tested, the analysis procedure was adapted so as to obtain a detection threshold of less than 1% of the occupational exposure level or TLV.

The procedure can be one of those described below or any other equivalent method:

- The concentrations of hydrochloric acid in the air were sampled by capturing a known volume of air on a cartridge impregnated with a buffer solution of Na CO /NaHCO . The samples prepared in this way were analyzed by ion chromatography (IC).
- The concentrations of organic gas were sampled by capturing a known volume of air on a cooled cartridge of adsorbent Tenax and active carbon. The samples prepared in this way were analyzed by gas chromatography (GC-FID) after thermal desorption.

The sample prepared in this way was then desorbed by a solution of carbon disulphide (CS) before being analyzed using a gas chromatograph (GC) equipped with a suitable detector (FID).

Test procedure

The tests were carried out at (20 ± 2) °C with a relative humidity between 40% and 70%.

The filtering fume hood being tested was placed in a closed test enclosure with an interior volume between 10 and 50 times the internal volume of the filtering fume hood.

The difference between the temperature inside the filtering fume hood and the temperature of the test enclosure must not exceed 5°C.

The chemical used for the test was introduced using a peristaltic pump, drop by drop, into a heated recipient in the center of the worktop in the filtering fume hood being tested. The system was set so as to produce the desired concentration to more or less 10% in the filtering fume hood for the whole duration of the test.

When necessary, the recipient was heated to slightly more than the boiling point of the test chemical in order to ensure instant evaporation.

Diagram of the test assembly (Evaporation and air sampling principle)

The air is sampled in three zones according to a procedure to be adapted according to the measurement protocol adopted:

- Zone 1: During the whole test, air is regularly sampled 30 cm downstream of the filtering system to check the purifying performance of the filtering fume hood being tested.
- Zone 2: As soon as the test begins, (when the evaporation concentration is stable) the air is sampled inside the fume hood, 30 cm upstream of the filtering system to check that the concentration released before the filters has evaporated.
- Zone 3: A few minutes after the beginning of the test, the air is sampled in the breathing zone to check that the concentration is less than 1% of TLV.

All of the necessary precautions must be taken during the test to avoid anything affecting the air samples between the sampling zone and the analyzer. Sampling must be carried out so as to provide a measurement result that is representative of the air analyzed (e.g. by using multipoint sampling grids).

Normal operating phase

The concentration of the chemical used for the test downstream of the filter must not exceed 1% of the authorized Occupational Exposure Limit or TLV. The evaporation must last for

the entire quantity of the chemical considered (given in the list of approved chemicals supplied by Erlab)

Detection phase

The concentration of the chemical used for the test downstream of the filter must not exceed 1% of the authorized Occupational Exposure Limit for class 1 filtering fume hoods and 50% of the authorized occupational exposure limit for class 2 filtering fume hoods.

Safety phase (for class 1 apparatus)

The concentration of the chemical used for the test downstream of the filter must not exceed 50 % of the authorized Occupational Exposure Limit. This phase must not be less than 1/12 of the duration of the normal operating phase.



Filtration test reports

The test report for each test performed must indicate:

- the reference of the test (name of the laboratory performing the test and date completed)
- the volume of the test enclosure that the filtering fume hood is placed in, the type and the reference of the filtering fume hood being tested
- the type and reference of the filtering fume hood being tested
- the type and reference of the filter(s) in the filtering fume hood being tested
- the nature of the chemical used for the test

The report for each phase of the test must indicate:

- Duration in hours
- Weight of the chemical(s)
- The concentration of the chemical used for the test in ppm by volume of air extracted.

Efficiency test certificate



The retention capacities recorded during the tests demonstrates the technological performance developed by Erlab. These results ensure the highest protection level to the GFH users.



Example of test carried out on a GFH with: 6 modular columns equipped with Neutrodine filtration technology.

Isopropanol	Cyclohexane	HCl (37%)
8526 g	10812 g	13104 g

The Erlab Difference - Our commitment to your safety

The air exhausted from our filters is so pure that we guarantee your safety. We make this guarantee because the following ingredients are built into every product - without all of them, your health and safety are compromised.

~ Stéphane Hauville
President & CEO

We need to be the best to keep you safe.

We invest in a state-of-the-art lab with top-notch engineers and chemists that are experts in molecular filtration with access to sophisticated analytical equipment that allows us to determine and optimize the retention capacity of our filters.

We continually work to ensure that our products are safe, improve the quality and technology of our products, and develop new products to offer greater protection for laboratory personnel.

Chemical Assessment for your chemical handlings

Assurance that a ductless hood is safe for you. We won't sell a hood without verifying that it is right for your chemical handlings. Your specific chemical handlings are analyzed by our in-house test laboratory. We determine if your chemical handlings can be done safely using our hood and filters. If yes, we let you know the most efficient filter type that is needed, the filter lifetime, and the best method of failing filtration indication.

We must be liable for your safety.

Certificate of validation

It is important to know what chemicals your hood can handle. A certificate that gives the precise details of the chemicals to be used, the filter type, and an estimation of the filter life expectancy is provided. This informs the user which chemicals are approved for use in the hood.

We must put it in writing.

Erlab Services

It is important that a safety specialist from the manufacturer regularly follows up with you to be sure the hood is in proper working order, check if any chemicals used in the hood have changed, and inform you when your filters need replacement.

We guarantee your safety for life.

List of chemicals not recommended*

GreenFumeHood® protects users with an innovating filter technology called Neutrodine Unisorb. This breakthrough filtering technology captures and neutralizes a large spectrum of chemicals such as acids, bases or volatile organic solvents on the same media. The few rare chemicals that follow are however not recommended for use under GreenFumeHood®.

Gaseous chemicals in normal temperature and pressure conditions with a very low boiling point (25°C, 1 Atmosphere).

Examples:

- Hydrogen (H₂)
- Helium and noble gases
- Methane
- Ethane
- Ethylene oxide
- Carbon monoxide
- Carbon dioxide
- Nitrogen monoxide
- Propyne
- Acetylene
- Propylene

Organophosphoric compounds :

Because of their very high toxicity (compounds that can be used as Chemical weapons).

Mercury :

In spite of the fact that this product is very well retained by Neutrodine®, it remains extremely toxic (TLV = 0,05 ppm) and very hard to detect.

Hydrogene Cyanide :

Immediately lethal

Our technical and sales departments remain fully at your service to study any specific project requirements. Do not hesitate to contact our local offices for more information.

*Non exhaustive list

Definitions of the column headings

The content of the various lists of chemicals given in this Chemical listing may differ from one table to another according to the relevance of the information linked to a specific chemical.

Chemical name: standard name or brand name of the chemical. For chemical names followed by an ®, the brand names have been registered by their owners

Formula: empirical chemical formula

CAS number (Chemical Abstract Standard Number): unique registration number of a chemical given by the American Chemical Society.

Suitable filter: the type of Erlab filter suitable for handling the chemical and/or providing the highest retention capacity:

- Neutrodine Unisorb
- HEPA: Powders 0.1 micron or higher
- PF: Prefilter to protect Hepa and/or Molecular Filters
- PF Ammonia: Prefilter to protect Hepa and/or Molecular Filters + Ammonia vapors

Retention capacity: filter retention capacity for the chemical, expressed in grams, during the normal operating phase described in AFNOR NFX 15-211: 2009, class 1

VP (Vapor Pressure): saturation vapor pressure at room temperature. Temperature is given in the corresponding box if data is not given at room temperature

MM: molar mass

Boiling point: The highest temperature a substance can reach before evaporating freely, expressed in °C at a pressure of 1 atmosphere

NIOSH 8h: American average limit values established by the National Institute for Occupational Safety and Health

France 8h: French average limit values established by the French Labour Ministry and published by the French INRS

AGS 8h: German average limit values established by the Deutsche Forschungsgemeinschaft

DFG 8h: German average limit values established by the Ausschuss für Gefahrstoffe

Japan 8h: Japanese average limit values established by the Japan Society for Occupational Health

China 8h: Chinese average limit values established by the GBZ 2.1-2007 – Occupational exposure limits for hazardous agent in the workplace

UK 8h: United Kingdom limit values established by the Health and Safety Executive.

European union 8h: European average limit values established by the Scientific Committee for Occupational Exposure Limits to Chemical Agents.

Detection:

- Sensor: Automatic method of detection. This equipment is required by class 1 of AFNOR NFX 15-211. The box is checked if the alarm can detect the relevant chemical. 3 versions of sensors are available : A for acids, F for formaldehyde and S for VOC (Volatile Organic Compounds).

Chemical agents by alphabetical order

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
1, 4-DIOXANE	C4H8O2	123-91-1	Neutrodine Unisorb	1518	4.95 kPa	88	101	-	20 ppm	20 ppm	20 ppm	10 ppm	70 mg/m ³	25 ppm	-	S
1,1,1-TRICHLOROETHANE	C2H3Cl3	71-55-6	Neutrodine Unisorb	1080	16.5 kPa	133	74	-	100 ppm	200 ppm	200 ppm	200 ppm	900 mg/m ³	-	-	
1,1,2,2-TETRABROMOETHANE	C2H2Br4	79-27-6	Neutrodine Unisorb	4680	0.003 kPa	346	239	-	1 ppm	-	-	-	-	0,5 ppm	-	
1,1,2,2-TETRACHLOROETHANE	C2H2Cl4	79-34-5	Neutrodine Unisorb	3108	0.622 kPa	168	146	1 ppm	1 ppm	1 ppm	1 ppm	1 ppm	-	-	-	
1,1'-BIPHENYL-4,4'-DIAMINE	C8H16N2	92-87-5	Neutrodine Unisorb	768	Low	184,3	400	-	0,001 ppm	-	-	-	-	-	50 ppm	S
1,1-DICHLOROETHANE	C2H4Cl2	75-34-3	Neutrodine Unisorb	456	30.5 kPa	98	57	100 ppm	100 ppm	100 ppm	100 ppm	100 ppm	-	-	100 ppm	S
1,2-DIBROMOETHANE	C2H4Br2	106-93-4	Neutrodine Unisorb	4680	1.55 kPa	188	131	-	-	-	-	-	-	0,5 ppm	-	S
1,2-DICHLOROBENZENE	C6H4Cl2	95-50-1	Neutrodine Unisorb	3060	0.18 kPa (125°C)	147	180	-	20 ppm	10 ppm	10 ppm	25 ppm	50 mg l/m ³	25 ppm	20 ppm	
1,2-DICHLOROETHANE	C2H4Cl2	107-06-2	Neutrodine Unisorb	1056	10.6 kPa	99	83	1 ppm	10 ppm	-	-	10ppm	7 mg/m ³	5 ppm	-	S
1,2-DICHLOROETHYLENE	C2H2Cl2	540-59-0	Neutrodine Unisorb	630	44.2 kPa	96	59	200 ppm	-	200 ppm	200 ppm	150 ppm	800 mg/m ³	-	-	S
1,2-EPOXY-3-ISOPROPOXYPROPANE	C6H12O2	4016-14-2	Neutrodine Unisorb	1788	1.2 kPa	116	127	-	50 ppm	-	-	-	-	50 ppm	-	
1,2-ETHANEDIOL	C2H6O2	107-21-1	Neutrodine Unisorb	1116	0.010 kPa	66	198	-	20 ppm	10 ppm	10 ppm	-	20 mg/m ³	20 ppm	-	S
1,3-BUTADIENE	C4H6	106-99-0	Neutrodine Unisorb	60	120 kPa (0°C)	54	-4,5	0,19 ppm	-	2 ppm	-	-	5 mg/m ³	10 ppm	-	S
1,3-CYCLOPENTADIENE	C5H6	542-92-7	Neutrodine Unisorb	840	58.5 kPa	66	42	75 ppm	75 ppm	-	-	-	-	-	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
1,3-DICHLOROPROPENE	C3H4Cl2	542-75-6	Neutrodine Unisorb	1578	3.73 kPa	110	103	1 ppm	-	-	-	-	4 mg/m ³	-	-	S
1,3-DICHLOROPROPYLENE	C3H4Cl2	542-75-6	Neutrodine Unisorb	1578	3.73 kPa	110	103	1 ppm	-	-	-	-	4 mg/m ³	-	-	S
1,3-DIOXOLANE	C3H6O2	646-06-0	Neutrodine Unisorb	984	14.6 kPa	74	78	-	-	100 ppm	100 ppm	-	-	-	-	-
1,3-DIVINYLBENZENE	C10H10	1321-74-0	Neutrodine Unisorb	1500	0.087 kPa	130	200	10 ppm	-	-	-	-	50 mg/m ³	10 ppm	-	S
1-AMINOBUTANE	C4H9NH2	109-73-9	Neutrodine Unisorb	456	12.2 kPa	73	78	-	-	-	2 ppm	-	-	-	-	S
1-AMINOPROPANE	C3H9NO	107-10-8	Neutrodine Unisorb	384	42.1 kPa	75	180	-	-	-	-	-	-	-	-	-
1-BUTANETHIOL	C4H10S	109-79-5	Neutrodine Unisorb	156	6.07 kPa	90,19	97	-	0,5 ppm	0,5 ppm	0,5 ppm	-	2 mg/m ³	-	-	S
1-BUTANOL	C4H10O	71-36-3	Neutrodine Unisorb	1760	0.86 kPa	74	117,5	-	-	100 ppm	100 ppm	25 ppm	100 mg/m ³	-	-	S
1-CHLORO BUTANE	C4H9Cl	109-69-3	Neutrodine Unisorb	1308	13.7 kPa	92	78,5	-	-	25 ppm	-	-	-	-	-	-
1-CHLORO-2,3-EPOXYPROPANE	C3H5ClO	106-89-8	Neutrodine Unisorb	1248	2.20 kPa	93	115	-	-	2 ppm	-	-	1 mg/m ³	0,5 ppm	-	-
1-HEXANOL	C6H14O	111-27-3	Neutrodine Unisorb	2755	0,1 kPa	102,2	157	-	-	50 ppm	-	-	-	-	-	-
1-MERCAPTOBUTANE	C4H10S	109-79-5	Neutrodine Unisorb	156	6.07 kPa	90	97	-	0,5 ppm	0,5 ppm	0,5 ppm	-	2 mg/m ³	-	-	S
1-METHYL-2-PYRROLIDINONE	C5H9NO	872-50-4	Neutrodine Unisorb	1500	0.04 kPa	99	202	-	10 ppm	20 ppm	20 ppm	1 ppm	-	10 ppm	10 ppm	-
1-PROPANETHIOL	C3H8S	107-03-9	Neutrodine Unisorb	78	20.6 kPa	76,2	67	0,3 ppm	-	-	-	-	-	-	-	-
1-PROPANOL	C3H8O	71-23-8	Neutrodine Unisorb	834	2.76 kPa	60	97	200 ppm	200 ppm	-	-	-	200 mg/m ³	-	-	S
2, 2'-DICHLORODIETHYL ETHER	C4H8OCl2	111-44-4	Neutrodine Unisorb	492	0.143 kPa	143	179	5 ppm	5 ppm	10 ppm	10 ppm	-	-	-	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
2, 4-DIMETHYL PENTANE	C7H16	108-08-7	Neutrodine Unisorb	918	36.6 kPa	100	80,5	-	-	-	-	-	-	-	-	-
2,4-DIMETHYL-3-PENTANONE	C7H14O	565-80-0	Neutrodine Unisorb	1458	6.87 kPa	114	124	-	-	-	-	-	-	-	-	-
2,6-DIMETHYL-4-HEPTANONE	C9H18O	108-83-8	Neutrodine Unisorb	1458	0.23 kPa	142	166	25 ppm	25 ppm	-	-	-	145 mg/m ³	25 ppm	-	S
2-AMINO 1-PROPANOL	C3H9NO	35320-23-1	Neutrodine Unisorb	384	0.1 kPa	75	180	-	-	-	-	-	-	-	-	-
2-AMINO BUTANE	C4H9NH2	13952-84-6	Neutrodine Unisorb	420	23 kPa	73	63	-	-	-	2 ppm	-	-	-	-	S
2-AMINO PYRIDINE	C5H6N2	504-29-0	Neutrodine Unisorb	1680	0.11 kPa	94	211	0,5 ppm	0,5 ppm	-	-	-	2 mg/m ³	-	-	-
2-AMINOETHANOL	C2H7NO	141-43-5	Neutrodine Unisorb	360	0.050 kPa	61	171	3 ppm	1 ppm	2 ppm	2 ppm	-	8 mg/m ³	1 ppm	-	-
2-AMINOPROPANE	C3H9N	75-31-0	Neutrodine Unisorb	234	78 kPa	59	34	-	5 ppm	5 ppm	5 ppm	-	12 mg/m ³	-	-	-
2-BUTANOL	C4H10O	78-92-2	Neutrodine Unisorb	1392	2.32 kPa	74	99,5	100 ppm	100 ppm	-	-	100 ppm	-	100 ppm	-	S
2-BUTANONE	C4H8O	78-93-3	Neutrodine Unisorb	984	12.6 kPa	72	80	200 ppm	200 ppm	200 ppm	200 ppm	200 ppm	300 mg/m ³	-	200 ppm	S
2-BUTENAL	C4H6O	4170-30-3	Neutrodine Unisorb	990	4.92 kPa	70	102	-	2 ppm	-	-	-	-	-	-	S
2-BUTOXYETHANOL	C6H14O2	111-76-2	Neutrodine Unisorb	2142	16.5 kPa	118	164	5 ppm	10 ppm	10 ppm	10 ppm	25 ppm	-	25 ppm	-	S
2-CHLOROACETALDEHYDE	C2H3OCl	107-20-0	Neutrodine Unisorb	744	13.3 kPa	78	90	-	-	-	-	-	-	-	-	S
2-CHLOROETHANAL	C2H3OCl	107-20-0	Neutrodine Unisorb	744	13.3 kPa	78	90	-	-	-	-	-	-	-	-	S
2-CHLOROETHANOL	C2H5OCl	107-07-3	Neutrodine Unisorb	1440	4.45 kPa (50°C)	81	129	-	-	1 ppm	1 ppm	-	-	-	-	-
2-CHLOROETHYL ALCOHOL	C2H5OCl	107-07-3	Neutrodine Unisorb	1440	4.45 kPa (50°C)	81	129	-	-	1 ppm	1 ppm	-	-	-	-	-

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
2-CHLOROPROPYLENE OXIDE	C3H5OCl	106-89-8	Neutrodine Unisorb	1248	2.20 kPa	93	115	-	-	2 ppm	-	-	1 mg/m ³	0,5 ppm	-	
2-ETHOXY ACETATE	C6H12O3	111-15-9	Neutrodine Unisorb	1518	0.24 kPa	132	157	0,5 ppm	2 ppm	2 ppm	2 ppm	5 ppm	30 mg/m ³	10 ppm	-	S
2-ETHOXYETHANOL	C4H10O2	110-80-5	Neutrodine Unisorb	3182	0.71 kPa	90	135	0,5 ppm	2 ppm	2 ppm	2 ppm	5 ppm	18 mg/m ³	10 ppm	-	S
2-ETHYL-1-HEXANOL	C8H18O	104-76-7	Neutrodine Unisorb	1500	0.019 kPa	130	190	-	-	10 ppm	10 ppm	-	-	-	-	
2-FURYL METHANOL	C5H6O2	98-00-0	Neutrodine Unisorb	1674	0.097 kPa	98	170	10 ppm	10 ppm	-	-	5 ppm	40 mg/m ³	-	-	S
2-HEPTANONE	C7H14O	110-43-0	Neutrodine Unisorb	1620	0.49 kPa	114,9	149	100 ppm	50 ppm	238 mg/m ³	-	-	-	50 ppm	-	S
2-HEXANONE	C6H12O	591-78-6	Neutrodine Unisorb	1488	1.54 kPa	100	127	1 ppm	5 ppm	5 ppm	5 ppm	5 ppm	20 mg/m ³	-	-	S
2-HYDROXYMETHYLFURAN	C5H6O2	98-00-0	Neutrodine Unisorb	1674	0.097 kPa	98	170	10 ppm	10 ppm	-	-	5 ppm	40 mg/m ³	-	-	S
2-METHOXY-2-METHYLPROPANE	C5H12O	1634-04-4	Neutrodine Unisorb	1869	27 kPa	88,2	55,05	-	50 ppm	50 ppm	50 ppm	-	-	25 ppm	50 ppm	
2-METHYL-1,3-BUTADIENE	C5H8	78-79-5	Neutrodine Unisorb	768	73.4 kPa	68	34	-	-	3 ppm	3 ppm	-	-	-	-	S
2-METHYL-1-PROPANOL	C4H10O	78-83-1	Neutrodine Unisorb	1542	1.39 kPa	74	108	50 ppm	50 ppm	100 ppm	100 ppm	50 ppm	-	-	-	S
2-METHYLBUTANE	C5H12	78-78-4	Neutrodine Unisorb	672	91.7 kPa	72	28	120 ppm	-	1000 ppm	1000 ppm	-	500 mg/m ³	-	1000 ppm	S
2-METHYLPROPYL ACETATE	C6H12O2	110-19-0	Neutrodine Unisorb	1740	2.39 kPa	116	117	150 ppm	150 ppm	62 ppm	100 ppm	150 ppm	-	150 ppm	-	S
2-METHYLPROPYL ESTER OF ACETIC ACID	C6H12O2	110-19-0	Neutrodine Unisorb	1740	2.39 kPa	116	117	150 ppm	150 ppm	62 ppm	100 ppm	150 ppm	-	150 ppm	-	S
2-PENTANONE	C5H10O	107-87-9	Neutrodine Unisorb	1500	4.97 kPa	86,13	101	150 ppm	200 ppm	-	-	-	-	200 ppm	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
2-PHENYL PROPANE	C9H12	98-82-8	Neutrodine Unisorb	1776	0.61 kPa	120	152	50 ppm	20 ppm	10 ppm	10 ppm	-	-	25 ppm	20 ppm	S
2-PROPANOL	C3H8O	67-63-0	Neutrodine Unisorb	1421	6.02 kPa	60	83	400 ppm	-	200 ppm	200 ppm	400 ppm	350 mg/m ³	-	-	S
2-PROPANONE	C3H6O	67-64-1	Neutrodine Unisorb	697	30.8 kPa	58	56,5	250 ppm	500 ppm	500 ppm	500 ppm	750 ppm	300 mg/m ³	-	-	S
2-PROPEN-1-OL	C3H6O	107-18-6	Neutrodine Unisorb	678	3.14 kPa	58	97	2 ppm	0,2 ppm	2 ppm	-	-	2 mg/m ³	2 ppm	-	S
2-PROPENAL	C3H4O	107-02-8	Neutrodine Unisorb	240	36.2 kPa	56	53	0,1 ppm	-	0,09 ppm	-	-	-	0,1 ppm	-	
2-PROPENAMIDE	C3H5NO	79-06-1	HEPA	-	0.014 kPa (75°C)	71,1	125	0,03 mg/m ³	0,1 ppm	0,07 mg/m ³	-	0,1 mg/m ³	0,3 mg/m ³	0,3 mg/m ³	-	
2-PROPENENITRILE	C3H3N	107-13-1	Neutrodine Unisorb	486	11.3 kPa	53	77	1 ppm	2 ppm	1,2 ppm	-	2 ppm	1 mg/m ³	2 ppm	-	
2-PROPENOIC ACID	C3H4O2	79-10-7	Neutrodine Unisorb	2040	2.45 kPa (50°C)	72	142	2 ppm	2 ppm	10 ppm	10 ppm	-	6 mg/m ³	-	10 ppm	
2-PROPENOL	C3H6O	107-18-6	Neutrodine Unisorb	678	3.14 kPa	58	97	2 ppm	0,2 ppm	2 ppm	-	-	2 mg/m ³	2 ppm	-	S
2-PROPYL ACETATE	C5H10O2	108-21-4	Neutrodine Unisorb	1572	5.59 kPa	102	88	-	250 ppm	-	100 ppm	250 ppm	-	-	-	S
2-PROPYLAMINE	C3H9N	75-31-0	Neutrodine Unisorb	234	78 kPa	59	34	-	5 ppm	5 ppm	5 ppm	-	12 mg/m ³	-	-	
2-PROPYN-1-OL	C3H4O	107-19-7	Neutrodine Unisorb	750	1.59 kPa	56	113	1 ppm	1 ppm	2 ppm	2 ppm	-	-	1 ppm	-	
2-PROPYNYL ALCOHOL	C3H4O	107-19-7	Neutrodine Unisorb	750	1.59 kPa	56	113	1 ppm	1 ppm	2 ppm	2 ppm	-	-	1 ppm	-	
3-AMINO-1-PROPANOL	C3H9NO	156-87-6	Neutrodine Unisorb	384	0.04 kPa	75	184	-	-	-	-	-	-	-	-	
3-CHLORO-1-PROPENE	C3H5Cl	107-05-1	Neutrodine Unisorb	462	76.5 kPa	76	45	1 ppm	1 ppm	-	-	-	2 mg/m ³	-	-	S
3-CRESOL	C7H8O	108-39-4	Neutrodine Unisorb	1578	0.019 kPa	108,14	203	2,3 ppm	-	-	-	-	-	5 ppm	-	
3-HYDROXYTOLUENE	C7H8O	108-39-4	Neutrodine Unisorb	1578	0.019 kPa	108,14	203	2,3 ppm	-	-	-	-	-	5 ppm	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
3-METHOXY-3-METHYL-1-BUTANOL	C6H14O2	56539-66-3	Neutrodine Unisorb	1514	0.125 kPa	118	173	-	-	-	-	-	-	-	-	S
3-METHYL PHENOL	C7H8O	108-39-4	Neutrodine Unisorb	1578	0.019 kPa	108,14	203	2,3 ppm	-	-	-	-	-	5 ppm	-	
3-METHYL-1-BUTANOL	C5H12O	123-51-3	Neutrodine Unisorb	1542	0.315 kPa	88	132	100 ppm	100 ppm	-	20 ppm	100 ppm	-	100 ppm	-	S
3-METHYL-3-PENTEN-2-ONE	C6H10O	565-62-8	Neutrodine Unisorb	2178	2 kPa	98	130	-	-	-	-	-	-	-	-	S
3-OCTANONE	C8H16O	106-68-3	Neutrodine Unisorb	1116	0.286 kPa	128,21	167	-	-	-	-	-	-	-	-	S
3-PENTANONE	C5H10O	96-22-0	Neutrodine Unisorb	1506	4.72 kPa	86	102	200 ppm	200 ppm	-	-	-	700 mg/m ³	200 ppm	5 ppm	S
4,4'-BIANILINE	C8H16N2	92-87-5	Neutrodine Unisorb	768	Low	184,3	400	-	0,001 ppm	-	-	-	-	-	50 ppm	S
4,4'-BIPHENYLDIAMINE	C8H16N2	92-87-5	Neutrodine Unisorb	768	Low	184,3	400	-	0,001 ppm	-	-	-	-	-	50 ppm	S
4,4'-DIAMINOBIHENYL	C8H16N2	92-87-5	Neutrodine Unisorb	768	Low	184,3	400	-	0,001 ppm	-	-	-	-	-	50 ppm	S
4-AMINOTOLUENE	C7H9N	106-49-0	Neutrodine Unisorb	465	1.74 kPa	107,2	200	-	-	-	-	-	-	-	-	
4-CRESOL	C7H8O	106-44-5	Neutrodine Unisorb	1578	0.017 kPa	108,14	202	2,3 ppm	-	-	-	-	-	5 ppm	-	
4-HYDROXYTOLUENE	C7H8O	106-44-5	Neutrodine Unisorb	1578	0.017 kPa	108,14	202	2,3 ppm	-	-	-	-	-	5 ppm	-	
4-METHYL 2-PENTANONE	C6H12O	108-10-1	Neutrodine Unisorb	1572	2.64 kPa	100	116	50 ppm	20 ppm	20 ppm	20 ppm	50ppm	-	50 ppm	-	S
4-METHYLANILINE	C7H9N	106-49-0	Neutrodine Unisorb	465	1.74 kPa	107,2	200	-	-	-	-	-	-	-	-	
4-TERT-BUTYL TOLUENE	C11H16	98-51-1	Neutrodine Unisorb	1560	0.090 kPa	148,24	193	10 ppm	10 ppm	-	-	-	6 mg/m ³	-	-	S
5-METHYL-3-HEPTANONE	C2H6O2	541-85-5	Neutrodine Unisorb	1116	0.27 kPa	138	157	25 ppm	10 ppm	10 ppm	10 ppm	-	130 mg/m ³	-	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
ABSOLUTE ALCOHOL	C2H6O	64-17-5	Neutrodine Unisorb	319	7.87 kPa	46	78	1000 ppm	1000 ppm	500 ppm	500 ppm	-	-	-	-	S
ACETALDEHYDE	C2H4O	75-07-0	Neutrodine Unisorb	33	120 kPa	44	20	18 ppm	100 ppm	50 ppm	50 ppm	50 ppm	-	-	-	S
ACETIC ACID	C2H4O2	64-19-7	Neutrodine Unisorb	5100	2.07 kPa	60	118	10 ppm	-	10 ppm	10 ppm	-	10 mg/m ³	-	-	S
ACETIC ANHYDRE	C4H6O3	108-24-7	Neutrodine Unisorb	2280	0.68 kPa	102	140	-	-	5 ppm	5 ppm	-	16 mg/m ³	0,5 ppm	-	S
ACETIC OXIDE	C4H6O3	108-24-7	Neutrodine Unisorb	2280	0.68 kPa	102	140	-	-	5 ppm	5 ppm	-	16 mg/m ³	0,5 ppm	-	S
ACETONE	C3H6O	67-64-1	Neutrodine Unisorb	697	30.8 kPa	58	56,5	250 ppm	500 ppm	500 ppm	500 ppm	750 ppm	300 mg/m ³	-	-	S
ACETONITRILE	C2H3N	75-05-8	Neutrodine Unisorb	*	11.9 kPa	41	82	20 ppm	40 ppm	20 ppm	20 ppm	-	30 mg/m ³	-	40 ppm	S
ACETYLENE	C2H2	74-86-2	Neutrodine Unisorb	48	4400 kPa	26	-84	2500 ppm	-	-	-	-	-	-	-	S
ACETYLENE DICHLORIDE	C2H2Cl2	540-59-0	Neutrodine Unisorb	630	44.2 kPa	96	59	200 ppm	-	200 ppm	200 ppm	150 ppm	800 mg/m ³	-	-	S
ACETYLENE TETRABROMIDE	C2H2Br4	79-27-6	Neutrodine Unisorb	4680	0.003 kPa	346	151	-	1 ppm	-	-	-	-	0,5 ppm	-	S
ACETYLENE TETRACHLORIDE	C2H2Cl4	79-34-5	Neutrodine Unisorb	3108	0.622 kPa	168	146	1 ppm	1 ppm	1 ppm	1 ppm	1 ppm	-	-	-	S
ACETYLSALICYCLIC ACID	C9H8O4	50-78-2	HEPA	-	-	180,2	-	5 mg/m ³ inhalable aerosol	-	-	-	-	5 mg/m ³	-	-	S
A-CHLOROTOLUENE	C7H7Cl	100-44-7	Neutrodine Unisorb	2544	0.164 kPa	127	179	-	1 ppm	-	-	-	-	0,5 ppm	-	S
ACROLEIC ACID	C3H4O2	79-10-7	Neutrodine Unisorb	2040	2.45 kPa (50°C)	72	142	2 ppm	2 ppm	10 ppm	10 ppm	-	6 mg/m ³	-	10 ppm	S
ACROLEIN	C3H4O	107-02-8	Neutrodine Unisorb	240	36.2 kPa	56	53	0,1 ppm	-	0,09 ppm	-	-	-	0,1 ppm	-	S
ACRYLAMIDE	C3H5NO	79-06-1	HEPA	-	0.014 kPa (75°C)	71,1	125	0,03 mg/m ³	0,1 ppm	0,07 mg/m ³	-	0,1 mg/m ³	0,3 mg/m ³	0,3 mg/m ³	-	S
ACRYLIC ACID	C3H4O2	79-10-7	Neutrodine Unisorb	2040	2.45 kPa (50°C)	72	142	2 ppm	2 ppm	10 ppm	10 ppm	-	6 mg/m ³	-	10 ppm	S

* More accurate estimation through eValiquet service

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapour pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
ACRYLIC ALDEHYDE	C3H4O	107-02-8	Neutrodine Unisorb	240	36.2 kPa	56	53	0,1 ppm	-	0,09 ppm	-	-	-	0,1 ppm	-	
ACRYLONITRILE	C3H3N	107-13-1	Neutrodine Unisorb	486	11.3 kPa	53	77	1 ppm	2 ppm	1,2 ppm	-	2 ppm	1 mg/m ³	2 ppm	-	
ALCOHOL	C2H6O	64-17-5	Neutrodine Unisorb	319	7.87 kPa	46	78	1000 ppm	1000 ppm	500 ppm	500 ppm	-	-	-	-	S
ALLYL ALCOHOL	C3H6O	107-18-6	Neutrodine Unisorb	678	3.14 kPa	58	97	2 ppm	0,2 ppm	2 ppm	-	-	2 mg/m ³	2 ppm	-	S
ALLYL ALDEHYDE	C3H4O	107-02-8	Neutrodine Unisorb	240	36.2 kPa	56	53	0,1 ppm	-	0,09 ppm	-	-	-	0,1 ppm	-	
ALLYL CHLORIDE	C3H5Cl	107-05-1	Neutrodine Unisorb	462	76.5 kPa	77	44,5	1 ppm	1 ppm	-	-	-	2 mg/m ³	-	-	S
ALLYLENE	C3H4	74-99-7	Neutrodine Unisorb	24	145 kPa (-25°C)	40	-23	1000 ppm	1000 ppm	-	-	-	-	-	-	
ALLYLGLYCIDYLETHER	C6H10O2	106-92-3	Neutrodine Unisorb	2292	1.77 kPa (50°C)	114	154	5 ppm	5 ppm	-	-	-	-	5 ppm	-	
ALLYLIC ALCOHOL	C3H6O	107-18-6	Neutrodine Unisorb	678	3.14 kPa	58	97	2 ppm	0,2 ppm	2 ppm	-	-	2 mg/m ³	2 ppm	-	S
ALUMINA	Al2O3	1344-28-1	PF + HEPA or Neutrodine Unisorb	-	1 Pa (1209°C)	101,96	2980	-	10 mg/m ³ respirable aerosol	-	"4 mg/m ³ inhalable aerosol 1.5 mg/m ³ inhalable respirable aerosol"	-	3 mg/m ³	"10 mg/m ³ inhalable aerosol 4 mg/m ³ respirable aerosol"	-	
ALUMINIUM	Al	7429-90-5	PF + HEPA or Neutrodine Unisorb	-	-	27	2327	10 mg/m ³ total dust	10 mg/m ³ inhalable aerosol	-	4 mg/m ³ inhalable aerosol	-	3 mg/m ³	-	-	
ALUMINUM OXIDE	Al2O3	1344-28-1	PF + HEPA or Neutrodine Unisorb	-	1 Pa @ 1209°C	101,96	2980	-	10 mg/m ³ respirable aerosol	-	"4 mg/m ³ inhalable aerosol 1.5 mg/m ³ inhalable respirable aerosol"	-	3 mg/m ³	"10 mg/m ³ inhalable aerosol 4 mg/m ³ respirable aerosol"	-	
ALUMINUM TRIOXIDE	Al2O3	1344-28-1	PF + HEPA or Neutrodine Unisorb	-	1 Pa (1209°C)	101,96	2980	-	10 mg/m ³ respirable aerosol	-	"4 mg/m ³ inhalable aerosol 1.5 mg/m ³ inhalable respirable aerosol"	-	3 mg/m ³	"10 mg/m ³ inhalable aerosol 4 mg/m ³ respirable aerosol"	-	
AMINO-BENZENE	C6H5NH2	62-53-3	Neutrodine Unisorb	1944	0.09 kPa	93	184	-	2 ppm	2 ppm	2 ppm	-	3 mg/m ³	-	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
AMINOCYCLOHEXANE	C6H11NH2	108-91-8	Neutrodine Unisorb	630	1.2 kPa	99	135	10 ppm	10 ppm	2 ppm	2 ppm	-	10 mg/m ³	10 ppm	-	S
AMINOETHANE	C2H7N	75-04-7	Neutrodine Unisorb	114	116 kPa	45	17	10 ppm	5 ppm	5 ppm	5 ppm	-	9 mg/m ³	-	-	S
AMINOMETHANE	CH5N	74-89-5	Neutrodine Unisorb	36	353 kPa	31	-7	10 ppm	-	10 ppm	10 ppm	-	5 mg/m ³	-	-	
AMMONIA (30% solution)	NH3	7664-41-7	PF Ammonia + Neutrodine Unisorb	384	1003 kPa	17	-33	25 ppm	10 ppm	20 ppm	20 ppm	25 ppm	20 mg/m ³	25 ppm	20 ppm	S
AMMONIUM CHLORIDE	NH4Cl	12125-02-9	PF + HEPA or Neutrodine Unisorb	-	0.13 kPa	53,49	-	10 mg/m ³	10 mg/m ³	-	-	-	10 mg/m ³	10 mg/m ³	-	
AMMONIUM CHLORIDE FUME	NH4Cl	12125-02-9	PF + HEPA or Neutrodine Unisorb	-	0.13 kPa	53,49	-	10 mg/m ³	10 mg/m ³	-	-	-	10 mg/m ³	10 mg/m ³	-	
AMYL ALCOHOL N	C5H12O	71-41-0	Neutrodine Unisorb	3416	0.259 kPa	88	138	-	-	-	20 ppm	-	100 mg/m ³	-	-	S
ANHYDROUS HYDROGEN BROMIDE	HBr	10035-10-6	Neutrodine Unisorb	1626	0.15 kPa	80,91	-66	-	-	6,7 mg/m ³	2 ppm	-	10mg/m ³	3 ppm	-	A
ANILINE	C6H5NH2	62-53-3	Neutrodine Unisorb	1944	0.09 kPa	93	184	-	2 ppm	2 ppm	2 ppm	-	3 mg/m ³	-	-	
AQUA FORTIS	HNO3	7697-37-2	Neutrodine Unisorb	1608	6.39 kPa	63	120	2 ppm	-	-	-	2 ppm	-	1 ppm	-	A
AQUA REGIA	HCL+HNO3	-	Neutrodine Unisorb	2040	-	44	120	-	-	-	-	-	-	-	-	A
AQUEOUS HYDROGEN BROMIDE (I.E.)	HBr	10035-10-6	Neutrodine Unisorb	1626	0.15 kPa	80,91	-66	-	-	6,7 mg/m ³	2 ppm	-	10mg/m ³	3 ppm	-	A
AQUEOUS HYDROGEN CHLORIDE (I.E.)	HCl aq. sol.	7647-01-0	Neutrodine Unisorb	2184	4103 kPa	37	120	-	-	2 ppm	2 ppm	-	-	1 ppm	5 ppm	A
ARSENIC (INORGANIC COMPOUNDS, AS AS)	As	7440-38-2	HEPA	-	1 Pa (280°C)	74,92	614	-	-	0,0083 mg/m ³	-	0,003 mg/m ³	0,01 mg/m ³	-	-	
ASBESTOS	Hydrated mineral silicates	1332-21-4	HEPA	-	-	-	-	0,1 fibers/cm3	0,01 fibres per cm3	0,01 fibres/cm3	-	0,15 fibers/cm3	0,8 mg/m ³ inhalable fraction	0,1 fibres per cm ³	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
ASPIRIN	C9H8O4	50-78-2	HEPA	-	-	180,2	-	5 mg/m ³ inhalable aerosol	-	-	-	-	5 mg/m ³	-	-	
ATRAZINE	C8H14ClN5	1912-24-9	HEPA	-	4 10-5 Pa	215,7	-	5 mg/m ³	5 mg/m ³	1 mg/m ³ inhalable aerosol	1 mg/m ³ inhalable aerosol	-	-	-	-	
AZINE	C5H5N	110-86-1	Neutrodine Unisorb	960	2.13 kPa	79	115	5 ppm	5 ppm	-	-	-	4 mg/m ³	5 ppm	-	S
BARIUM CHLORIDE	BaCl2.2H2O	10326-38-9	PF + HEPA or Neutrodine Unisorb	-	-	244,26	-	-	-	-	-	-	-	-	-	
BENZENAMINE	C6H5NH2	62-53-3	Neutrodine Unisorb	1944	0.09 kPa	93	184	-	2 ppm	2 ppm	2 ppm	-	3 mg/m ³	-	-	
BENZENE	C6H6	71-43-2	Neutrodine Unisorb	948	12.7 kPa	78	80	0.1 ppm	1 ppm	1 ppm	-	10 ppm	6 mg/m ³	-	-	S
BENZENE CHLORIDE	C6H5Cl	108-90-7	Neutrodine Unisorb	2364	1.6 kPa	113	133	-	5 ppm	10 ppm	10 ppm	10 ppm	50 mg/m ³	1 ppm	-	S
BENZINE 35 80	C8H16N2	92-87-5	Neutrodine Unisorb	768	Low	184,3	400	-	0,001 ppm	-	-	-	-	-	50 ppm	S
BENZYL ALCOHOL	C6H5CH2OH	100-51-6	Neutrodine Unisorb	1347	0.015 kPa	108,14	203	-	-	-	-	-	-	-	-	
BENZYL CHLORIDE	C7H7Cl	100-44-7	Neutrodine Unisorb	2544	0.164 kPa	127	179	-	1 ppm	-	-	-	-	0,5 ppm	-	S
BERYLLIUM COMPOUNDS (AS BE)	Be	7440-41-7	HEPA	-	-	9,01	2471	0,0005 mg/m ³	0,002 mg/m ³	-	-	0,001 mg/m ³	0,0005 mg/m ³	-	-	
BET	C21H20N3Br	1239-45-8	Neutrodine Unisorb	1560	-	314	238	-	-	-	-	-	-	-	-	
BETA-AMINOETHYL ALCOHOL	C2H7NO	141-43-5	Neutrodine Unisorb	360	0.050 kPa	61	171	3 ppm	1 ppm	2 ppm	2 ppm	-	8 mg/m ³	1 ppm	-	
BETA-CHLOROPRENE	C4H4Cl	126-99-8	Neutrodine Unisorb	456	29.5 kPa	87	60	-	10 ppm	-	-	-	4 mg/m ³	10 ppm	-	S
BETA-METHYL ACROLEIN	C4H6O	4170-30-3	Neutrodine Unisorb	990	4.92 kPa	70	102	-	2 ppm	-	-	-	-	-	-	S
BETA-METHYLPROPYL ETHANOATE	C6H12O2	110-19-0	Neutrodine Unisorb	1740	2.39 kPa	116	117	150 ppm	150 ppm	62 ppm	100 ppm	150 ppm	-	150 ppm	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapour pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
BICYCLOPENTADIENE	C10H12	77-73-6	Neutrodine Unisorb	1524	-	132	167	5 ppm	5 ppm	0,5 ppm	0,5 ppm	-	25 mg/m ³	5 ppm	-	S
BIOTITE	K(Mg, Fe) 3AlSi ₃ O ₁₀ (F, OH) ₂	12001-26-2	PF + HEPA or Neutrodine Unisorb	-	-	797	-	3 mg/m ³ respirable fraction	-	-	-	-	"2 mg/m ³ inhalable fraction 1.5 mg/m ³ respirable fraction"	"10mg/m ³ inhalable fraction 0.8 mg/m ³ respirable fraction"	-	
BLEU DE THYMOL	C27H39O5S	76-61-9	PF + HEPA or Neutrodine Unisorb	-	-	466,59	-	-	-	-	-	-	-	-	-	
BORAX	Na ₂ B ₄ O ₇ • 10H ₂ O	1303-96-4	PF + HEPA or Neutrodine Unisorb	-	-	381,4	-	5 mg/m ³	5 mg/m ³	-	0.75 mg/m ³ inhalable aerosol	-	-	5 mg/m ³	-	
BORON OXIDE	B ₂ O ₃	1303-86-2	PF + HEPA or Neutrodine Unisorb	-	-	69,62	-	10 mg/m ³ total dust	10 mg/m ³	-	-	-	-	-	-	
BORON TRIFLUORIDE	BF ₃	02-07-7637	No filtration	-	101 kPa	67,81	-	-	-	0,35 ppm	-	0,3 ppm	-	-	-	
BROMINE	Br ₂	7726-95-6	Neutrodine Unisorb	858	28.7 kPa	160	59	0,1 ppm	-	0,7 mg/m ³	-	0,1 ppm	0,6 mg/m ³	0,1 ppm	0,1 ppm	
BROMOCHLOROMETHANE	CH ₂ BrCl	74-97-5	Neutrodine Unisorb	1992	19,5 kPa	129	68	200 ppm	200 ppm	-	-	-	-	-	-	
BROMOETHANE	C ₂ H ₅ Br	74-96-4	Neutrodine Unisorb	1080	62.5 kPa	113	38,5	-	200 ppm	-	-	-	-	-	-	S
BROMOETHENE	C ₂ H ₃ Br	593-60-2	Neutrodine Unisorb	48	141 kPa	107	16	-	-	-	-	-	-	-	-	
BROMOETHYLENE	C ₂ H ₃ Br	593-60-2	Neutrodine Unisorb	48	141 kPa	107	16	-	-	-	-	-	-	-	-	
BROMOFORM	CHBr ₃	75-25-2	Neutrodine Unisorb	900	0.726 kPa	253	149,5	0,5 ppm	0,5 ppm	-	-	1 ppm	-	-	-	S
BUTANOIC ACID	C ₄ H ₈ O ₂	107-92-6	Neutrodine Unisorb	2280	0.221 kPa	88	163,5	-	-	-	-	-	-	-	-	
BUTYL ACRYLATE	C ₇ H ₁₂ O ₂	141-32-2	Neutrodine Unisorb	2064	0.731 kPa	128	146	-	2 ppm	2 ppm	2 ppm	-	25 mg/m ³	1 ppm	-	S
BUTYL ALCOHOL	C ₄ H ₁₀ O	71-36-3	Neutrodine Unisorb	1760	0.86 kPa	74	117,5	-	-	100 ppm	100 ppm	25 ppm	100 mg/m ³	-	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
BUTYL ALCOHOL SEC	C4H10O	78-92-2	Neutrodine Unisorb	1392	2.32 kPa	74	99,5	100 ppm	100 ppm	-	-	101 ppm	-	100 ppm	-	S
BUTYL ALCOHOL TER	C4H10O	75-65-0	Neutrodine Unisorb	1170	5.42 kPa	74	83	100 ppm	100 ppm	20 ppm	20 ppm	50 ppm	-	100 ppm	-	S
BUTYL CARBINOL	C5H12O	71-41-0	Neutrodine Unisorb	3416	0.259 kPa	88	138	-	-	-	20 ppm	-	100 mg/m ³	-	-	S
BUTYL CELLOSOLVE®	C6H14O2	111-76-2	Neutrodine Unisorb	2142	16.5 kPa	118	164	5 ppm	10 ppm	10 ppm	10 ppm	25 ppm	-	25 ppm	-	S
BUTYL ETHER	C8H18O	142-96-1	Neutrodine Unisorb	1302	0.898 kPa	130	142	-	-	-	-	-	-	-	-	S
BUTYL GLYCIDYL ETHER	C7H14O2	06-08-2426	Neutrodine Unisorb	1998	426 Pa	130,18	164	-	25 ppm	-	-	-	60 mg/m ³	25 ppm	-	S
BUTYL GLYCOL	C6H14O2	111-76-2	Neutrodine Unisorb	2142	16.5 kPa	118	164	5 ppm	10 ppm	10 ppm	10 ppm	25 ppm	-	25 ppm	-	S
BUTYL LACTATE	C7H14O3	138-22-7	Neutrodine Unisorb	2298	0.053 kPa	146	188	5 ppm	5 ppm	-	-	-	25 mg/m ³	5 ppm	-	S
BUTYL METACRYLATE	C18H14O2	97-88-1	Neutrodine Unisorb	2250	266 Pa	142	164	-	10 ppm	-	-	-	-	-	-	S
BUTYL VINYL ETHER	C6H12O	111-34-2	Neutrodine Unisorb	1044	6.65 kPa	100	94	-	-	-	-	-	-	-	-	S
BUTYLENE HYDRATE	C4H10O	78-92-2	Neutrodine Unisorb	1392	2.32 kPa	74	99,5	100 ppm	100 ppm	-	-	102 ppm	-	100 ppm	-	S
BUTYRIC ACID	C4H8O2	107-92-6	Neutrodine Unisorb	2280	0.221 kPa	88	163,5	-	-	-	-	-	-	-	-	S
BVE	C6H12O	111-34-2	Neutrodine Unisorb	1044	6.65 kPa	100	94	-	-	-	-	-	-	-	-	S
CADMIUM DUST (AS CD)	Cd	7440-43-9	HEPA	-	-	112,41	765	0.01 mg/m ³	-	-	-	-	0.01 mg/m ³	-	-	S
CADMIUM FUME (AS CD)	Cd	7440-43-9	HEPA	-	-	112,4	767	-	0.05 mg/m ³	-	-	0,05 mg/m ³	-	-	-	S
CALCIUM CARBONATE	CaCO3	1317-65-3	PF + HEPA or Neutrodine Unisorb	-	-	100,09	-	"10 mg/m ³ total dust 5 mg/m ³ respirable aerosol"	-	-	-	-	"8 mg/m ³ inhalable fraction 4 mg/m ³ respirable fraction"	10 mg/m ³ total dust	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
CALCIUM HYDRATE	Ca(OH)2	1305-62-0	PF + HEPA or Neutrodine Unisorb	-	-	74,09	-	5 mg/m ³	5 mg/m ³	1 mg/m ³ inhalable fraction	1 mg/m ³ inhalable fraction	-	-	5 mg/m ³	-	
CALCIUM HYDROXIDE	Ca(OH)2	1305-62-0	PF + HEPA or Neutrodine Unisorb	-	-	74,09	-	5 mg/m ³	5 mg/m ³	1 mg/m ³ inhalable fraction	1 mg/m ³ inhalable fraction	-	-	5 mg/m ³	-	
CALCIUM OXIDE	CaO	1305-78-8	PF + HEPA or Neutrodine Unisorb	-	-	56,08	-	2 mg/m ³	2 mg/m ³	1 mg/m ³ inhalable fraction	1 mg/m ³ inhalable fraction	-	2 mg/m ³	2 mg/m ³	-	
CALCIUM SULFATE	Ca(SO4). 2H2O	7778-18-9	PF + HEPA or Neutrodine Unisorb	-	-	172,17	-	5 mg/m ³ respirable aerosol	-	6 mg/m ³ respirable aerosol	4 mg/m ³ inhalable aerosol	-	-	-	-	
CARBON BLACK	C	1333-86-4	HEPA	-	-	-	-	3,5 mg/m ³	3,5 mg/m ³	-	-	"1 mg/m ³ respirable dust 4 mg/m ³ total dust"	4 mg/m ³ inhalable fraction	3.5 mg/m ³	-	
CARBON BROMIDE	CBr4	558-13-4	Neutrodine Unisorb	2700	5.3 kPa	332	189,5	0,1 ppm	0,1 ppm	-	-	-	1,5 mg/m ³	-	-	
CARBON DIOXIDE	CO2	124-38-9	No filtration	-	3483 kPa (0°C)	44	-	5000 ppm	5000 ppm	5000 ppm	5000 ppm	5000 ppm	9000 mg/m ³	5000 ppm	5000 ppm	
CARBON DISULFIDE	CS2	75-15-0	Neutrodine Unisorb	294	48.2 kPa	76	46	1 ppm	5 ppm	10 ppm	5 ppm	1 ppm	5 mg/m ³	-	-	S
CARBON MONOXIDE	CO	630-08-0	No filtration	-	> 3545 kPa	28	-191,5	35 ppm	50 ppm	30 ppm	30 ppm	50 ppm	20 mg/m ³	30 ppm	20 ppm	
CARBON TETRABROMIDE	CBr4	558-13-4	Neutrodine Unisorb	2700	5.3 kPa	332	189,5	0,1 ppm	0,1 ppm	-	-	-	1,5 mg/m ³	-	-	
CARBON TETRACHLORIDE	CCl4	56-23-5	Neutrodine Unisorb	2790	15.2 kPa	154	77	-	2 ppm	0,5 ppm	0,5 ppm	-	15 mg/m ³	-	-	
CAUSTIC POTASH	KOH	1310-58-3	PF + HEPA or Neutrodine Unisorb	-	0.13 kPa (1044°C)	56,11	-	-	-	-	-	-	-	-	-	
CAUSTIC SODA	NaOH	1310-73-2	PF + HEPA or Neutrodine Unisorb	-	-	40	1390	-	2 mg/m ³	-	-	-	-	-	-	
CELLOSOLVE «ACETATE»	C6H12O3	111-15-9	Neutrodine Unisorb	1518	0.24 kPa	132	157	0,5 ppm	2 ppm	2 ppm	2 ppm	5 ppm	30 mg/m ³	10 ppm	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapour pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
CELLOSOLVE®	C4H10O2	110-80-5	Neutrodine Unisorb	3182	0.71 kPa	90	135	0,5 ppm	2 ppm	2 ppm	2 ppm	5 ppm	18 mg/m ³	10 ppm	-	S
CELLULOSE	(C6H10O5) _n	9004-34-6	PF + HEPA or Neutrodine Unisorb	-	-	160,000-560,000	-	10 mg/m ³ total dust	10 mg/m ³ inhalable aerosol	-	-	-	10 mg/m ³	10 mg/m ³ inhalable aerosol	-	-
CHLORINE	Cl ₂	7782-50-5	Neutrodine Unisorb	858	780 kPa (50°C)	35	-34,5	-	-	0,5 ppm	0,5 ppm	0,5 ppm	-	-	-	-
CHLORINE DIOXIDE	ClO ₂	10049-04-4	Neutrodine Unisorb	270	101 kPa	68	10	0,1 ppm	0,1 ppm	0,1 ppm	0,1 ppm	-	0,3mg/m ³	0,1 ppm	-	-
CHLORINE OXIDE	ClO ₂	10049-04-4	Neutrodine Unisorb	270	101 kPa	68	10	0,1 ppm	0,1 ppm	0,1 ppm	0,1 ppm	-	0,3mg/m ³	0,1 ppm	-	-
CHLORO-1-NITROPROPANE 1	C3H6ClNO2	600-25-9	Neutrodine Unisorb	2250	0.79 kPa	123	171	2 ppm	2 ppm	-	-	-	-	-	-	S
CHLOROBENZENE	C6H5Cl	108-90-7	Neutrodine Unisorb	2364	1.6 kPa	113	133	-	5 ppm	10 ppm	10 ppm	10 ppm	50 mg/m ³	1 ppm	-	S
CHLOROBROMOMETHANE	CH2BrCl	74-97-5	Neutrodine Unisorb	1992	19.5 kPa	129	68	200 ppm	200 ppm	-	-	-	-	-	-	-
CHLOROBUTADIENE	C4H4Cl	126-99-8	Neutrodine Unisorb	456	29.5 kPa	87	60	-	10 ppm	-	-	-	4 mg/m ³	10 ppm	-	S
CHLOROETHANE	C2H5Cl	75-00-3	Neutrodine Unisorb	48	160 kPa	65	12	-	100 ppm	40 ppm	-	-	-	-	-	S
CHLOROETHENE	C2H3Cl	75-01-4	Neutrodine Unisorb	48	355 kPa	61	-14	-	1 ppm	3 ppm	-	2 ppm	10 mg/m ³	-	-	S
CHLOROETHYLENE	C2H3Cl	75-01-4	Neutrodine Unisorb	48	355 kPa	61	-14	-	1 ppm	3 ppm	-	2 ppm	10 mg/m ³	-	-	S
CHLOROFORM	CHCl ₃	67-66-3	Neutrodine Unisorb	780	26.2 kPa	119	61	-	2 ppm	0,5 ppm	0,5 ppm	3 ppm	20 mg/m ³	-	-	-
CHLOROMETHANE	CH3Cl	74-87-3	Neutrodine Unisorb	18	574 kPa	51	-24	-	50 ppm	50 ppm	50 ppm	-	60 mg/m ³	-	-	S
CHLOROPRENE	C4H4Cl	126-99-8	Neutrodine Unisorb	456	29.5 kPa	87	60	-	10 ppm	-	-	-	4 mg/m ³	10 ppm	-	S
CHLOROTHENE	C2H3Cl ₃	71-55-6	Neutrodine Unisorb	1080	16.5 kPa	133	74	-	100 ppm	200 ppm	200 ppm	200 ppm	900 mg/m ³	-	-	-

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapour pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
CHLOROTOLUENE (ORTHO)	C7H7Cl	106-43-4	Neutrodine Unisorb	2544	0.482 kPa	126	159	-	-	-	-	-	-	-	-	S
CHLORURE DE SODIUM	NaCl	7647-14-5	PF + HEPA or Neutrodine Unisorb	-	-	58,44	-	-	-	-	-	-	-	-	-	-
CHROMIC ACID	CrO3	1333-82-0	HEPA + Neutrodine Unisorb	270	1.87 kPa	151,99	4000	-	0,05 mg/m ³	-	-	-	-	-	-	-
CHROMIC OXIDE	CrO3	1333-82-0	HEPA + Neutrodine Unisorb	270	1.87 kPa	151,99	4000	-	0,05 mg/m ³	-	-	-	-	-	-	-
CHROMIUM(VI) OXIDE (1:3)	CrO3	1333-82-0	HEPA + Neutrodine Unisorb	270	1.87 kPa	151,99	4000	-	0,05 mg/m ³	-	-	-	-	-	-	-
CINAMENE	C8H8	100-42-5	Neutrodine Unisorb	1260	0.81 kPa	104	146	50 ppm	50 ppm	20 ppm	20 ppm	50 ppm	50 mg/m ³	100 ppm	-	S
CLAY	Al2Si2O5(OH)4	1332-58-7	PF + HEPA or Neutrodine Unisorb	-	-	-	-	"10 mg/m ³ total dust 5 mg/m ³ respirable fraction"	10 mg/m ³ respirable aerosol	-	-	-	-	2 mg/m ³ respirable aerosol	-	-
COPPER (DUSTS AND MISTS, AS CU)	Cu	7440-50-8	PF + HEPA or Neutrodine Unisorb	-	-	63,55	2562	-	1 mg/m ³	-	0,01 mg/m ³	-	1 mg/m ³	-	-	-
COPPER(II) OXIDE FUME	CuO	1317-38-0	PF + HEPA or Neutrodine Unisorb	-	-	79,55	-	-	-	-	-	-	-	-	-	-
CRESOL ALL ISOMERS	C7H8O	1319-77-3	Neutrodine Unisorb	1578	0.018 kPa	108	191	2,3 ppm	5 ppm	-	-	5 ppm	10 mg/m ³	5 ppm	-	-
CROTONALDEHYDE	C4H6O	4170-30-3	Neutrodine Unisorb	990	4.92 kPa	70	102	-	2 ppm	-	-	-	-	-	-	S
CUMENE	C9H12	98-82-8	Neutrodine Unisorb	1776	0.61 kPa	120	152	50 ppm	20 ppm	10 ppm	10 ppm	-	-	25 ppm	20 ppm	S
CUMOL	C9H12	98-82-8	Neutrodine Unisorb	1776	0.61 kPa	120	152	50 ppm	20 ppm	10 ppm	10 ppm	-	-	25 ppm	20 ppm	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
CYANOMETHANE	C2H3N	75-05-8	Neutrodine Unisorb	*	11.9 kPa	41	82	20 ppm	40 ppm	20 ppm	20 ppm	-	30 mg/m ³	-	40 ppm	S
CYCLOHEXANE	C6H12	110-82-7	Neutrodine Unisorb	1802	13 kPa	84	81	300 ppm	200 ppm	200 ppm	200 ppm	-	250 mg/m ³	100 ppm	-	S
CYCLOHEXANOL	C6H12O	108-93-0	Neutrodine Unisorb	1830	0.1 kPa	100	162	50 ppm	50 ppm	50 ppm	50 ppm	25 ppm	100 mg/m ³	50 ppm	-	S
CYCLOHEXANONE	C6H10O	108-94-1	Neutrodine Unisorb	2094	0.53 kPa	98	157	25 ppm	10 ppm	20 ppm	-	20 ppm	50 mg/m ³	10 ppm	-	S
CYCLOHEXENE	C6H10	110-83-8	Neutrodine Unisorb	1344	11.8 kPa	82	83	300 ppm	300 ppm	-	-	-	-	300 ppm	-	S
CYCLOHEXYL ALCOHOL	C6H12O	108-93-0	Neutrodine Unisorb	1830	0.1 kPa	100	162	50 ppm	50 ppm	50 ppm	50 ppm	25 ppm	100 mg/m ³	50 ppm	-	S
CYCLOHEXYL KETONE	C6H10O	108-94-1	Neutrodine Unisorb	2094	0.53 kPa	98	157	25 ppm	10 ppm	20 ppm	-	20 ppm	50 mg/m ³	10 ppm	-	S
CYCLOHEXYLAMINE	C6H11NH2	108-91-8	Neutrodine Unisorb	630	1.2 kPa	99	135	10 ppm	10 ppm	2 ppm	2 ppm	-	10 mg/m ³	10 ppm	-	S
CYCLOPENTANE	C5H10	287-92-3	Neutrodine Unisorb	852	42.3 kPa	70	49	600 ppm	600 ppm	-	-	-	-	619 ppm	-	S
DEA	C4H11NO2	111-42-2	Neutrodine Unisorb	840	< 1 Pa	105	217	3 ppm	3 ppm	-	1 mg/m ³ inhalable aerosol	-	-	3 ppm	-	S
DECANE	C10H22	124-18-5	Neutrodine Unisorb	1908	190 Pa	142	174	-	-	-	-	-	-	-	-	S
DIACETONE	C6H12O2	123-42-2	Neutrodine Unisorb	1536	0.224 kPa	116	168	50 ppm	50 ppm	20 ppm	20 ppm	-	240 mg/m ³	50 ppm	-	S
DIACETONE ALCOHOL	C6H12O2	123-42-2	Neutrodine Unisorb	1536	0.224 kPa	116	168	50 ppm	50 ppm	20 ppm	20 ppm	-	240 mg/m ³	50 ppm	-	S
DIAMINE	N2H4	302-01-2	Neutrodine Unisorb	768	1.3 kPa	32,05	-	-	0,1 ppm	0,017 ppm	-	-	0,06 mg/m ³	0,02 ppm	-	S
DIBUTYL ETHER	C8H18O	142-96-1	Neutrodine Unisorb	1302	0.898 kPa	130	142	-	-	-	-	-	-	-	-	S
DICHLOROMETHANE	CH2Cl2	75-09-2	Neutrodine Unisorb	604	58.2 kPa	85	40	-	50 ppm	50 ppm	50 ppm	100 ppm	200 mg/m ³	-	100 ppm	S
DICHLOROPROPANE 1, 2	C3H6Cl2	78-87-5	Neutrodine Unisorb	1350	6.62 kPa	113	97	-	75 ppm	-	-	1 ppm	350 mg/m ³	-	-	S

* More accurate estimation through eValiquet service

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
DICYCLOPENTADIENE	C10H12	77-73-6	Neutrodine Unisorb	1524	-	132	167	5 ppm	5 ppm	0,5 ppm	0,5 ppm	-	25 mg/m ³	-	-	S
DIETHAMINE	C4H11N	109-89-7	Neutrodine Unisorb	504	30.1 kPa	73	56	10 ppm	5 ppm	5 ppm	5 ppm	-	-	5 ppm	-	S
DIETHANOLAMINE	C4H11NO2	111-42-2	Neutrodine Unisorb	840	< 1 Pa	105	217	3 ppm	3 ppm	-	1 mg/m ³ inhalable aerosol	-	-	3 ppm	-	
DIETHYL ETHER	C4H10O	60-29-7	Neutrodine Unisorb	822	71.7 kPa	74	35	-	100 ppm	400 ppm	400 ppm	400 ppm	300 mg/m ³	-	-	S
DIETHYL KETONE	C5H10O	96-22-0	Neutrodine Unisorb	1506	4.72 kPa	86	102	200 ppm	200 ppm	-	-	-	700 mg/m ³	200 ppm	5 ppm	S
DIETHYL OXIDE	C4H10O	60-29-7	Neutrodine Unisorb	822	71.7 kPa	74	35	-	100 ppm	400 ppm	400 ppm	400 ppm	300 mg/m ³	-	-	S
DIETHYLAMINE	C4H11N	109-89-7	Neutrodine Unisorb	384	30.1 kPa	73	56	10 ppm	5 ppm	5 ppm	5 ppm	-	-	5 ppm	-	S
DIETHYLAMINOETHANOL-2	C6H15NO	100-37-8	Neutrodine Unisorb	1218	30 kPa @ 125°C	117	162	10ppm	10 ppm	5 ppm	5 ppm	-	50mg/m ³	10 ppm	-	S
DIETHYLENE DIOXIDE	C4H8O2	123-91-1	Neutrodine Unisorb	1518	4.95 kPa	88	101	-	20 ppm	20 ppm	20 ppm	10 ppm	70 mg/m ³	25 ppm	-	S
DIETHYLENE GLYCOL MONOBUTYL ETHER	C8H18O3	112-34-5	Neutrodine Unisorb	1899	0.029 kPa	162,23	224	-	10 ppm	10 ppm	10 ppm	-	-	12 ppm	-	
DIETHYLENE OXIDE	C4H8O	109-99-9	Neutrodine Unisorb	1044	21.6 kPa	72	65	200 ppm	50 ppm	50 ppm	50 ppm	50 ppm	300 mg/m ³	50 ppm	-	S
DIETHYLENE TRIAMINE	C4H13N3	111-40-0	Neutrodine Unisorb	570	0.053 kPa	103	207	1 ppm	1 ppm	-	-	-	4 mg/m ³	1 ppm	-	
DIISOBUTYL KETONE	C9H18O	108-83-8	Neutrodine Unisorb	1458	0.23 kPa	142	166	25 ppm	25 ppm	-	-	-	145 mg/m ³	25 ppm	-	S
DIISOPROPYL ETHER	C6H14O	108-20-3	Neutrodine Unisorb	2486	8.35 kPa	102	69	500 ppm	250 ppm	200 ppm	200 ppm	-	-	250 ppm	-	
DIISOPROPYL KETONE	C7H14O	565-80-0	Neutrodine Unisorb	1458	6.87 kPa	114	124	-	-	-	-	-	-	-	-	
DIISOPROPYL OXIDE	C6H14O	108-20-3	Neutrodine Unisorb	2486	8.35 kPa	102	69	-	250 ppm	200 ppm	200 ppm	-	-	250 ppm	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
DIISOPROPYLAMINE	C6H15N	108-18-9	Neutrodine Unisorb	318	10.7 kPa	101	83	-	5 ppm	-	-	-	-	5 ppm	-	S
DIMETHOXYMETHANE	C3H8O2	109-87-5	Neutrodine Unisorb	744	53.1 kPa	76	43	1000 ppm	1000 ppm	1000 ppm	1000 ppm	-	-	1000 ppm	-	S
DIMETHYL AMINE	C2H7N	124-40-3	Neutrodine Unisorb	90	203 kPa	45	7	10 ppm	1 ppm	2 ppm	2 ppm	-	5 mg/m ³	2 ppm	-	S
DIMETHYL BENZENE (AND ISOMERS)	C8H10	95-47-6	Neutrodine Unisorb	1920	0.88 kPa	106	138	100 ppm	50 ppm	100 ppm	100 ppm	50 ppm	50 mg/m ³	50 ppm	50 ppm	S
DIMETHYL CARBINOL	C3H8O	67-63-0	Neutrodine Unisorb	1421	6.02 kPa	60	83	400 ppm	-	200 ppm	200 ppm	400 ppm	350 mg/m ³	-	-	S
DIMETHYL ETHER	C2H6O	115-10-6	Neutrodine Unisorb	18	273 kPa (0°C)	46	-25	-	1000 ppm	1000 ppm	1000 ppm	-	-	400 ppm	-	S
DIMETHYL KETONE	C3H6O	67-64-1	Neutrodine Unisorb	697	30.8 kPa	58	56,5	250 ppm	500 ppm	500 ppm	500 ppm	750 ppm	300 mg/m ³	-	-	S
DIMETHYL SULFOXIDE	C2H6SO	67-68-5	Neutrodine Unisorb	1170	0.1 kPa	78	189	-	-	-	50 ppm	-	-	-	-	S
DIMETHYLACETONE	C5H10O	96-22-0	Neutrodine Unisorb	1506	4.72 kPa	86	102	200 ppm	200 ppm	-	-	-	700 mg/m ³	200 ppm	5 ppm	S
DIMETHYLFORMAMIDE	C3H7NO	68-12-2	Neutrodine Unisorb	1944	0.439 kPa	73	153	10 ppm	5 ppm	5 ppm	5 ppm	10 ppm	20 mg/m ³	-	-	S
DIMETHYLMETHANE	C3H8	74-98-6	Neutrodine Unisorb	12	939 kPa	44	-42	1000 ppm	-	1000 ppm	1000 ppm	-	-	-	-	S
DINITROGEN MONOXIDE	N2O	10024-97-2	No filtration	-	5197 kPa	44,01	-91	25 ppm	-	100 ppm	100 ppm	-	-	100 ppm	-	S
DIOXIN	C12H4Cl4O2	1746-01-6	No filtration	-	1.733191e-005 (25°C)	322	-	-	-	-	0.00000001 g/m3 inhalable aerosol	-	-	-	-	S
DIOXINE	C12H4Cl4O2	1746-01-6	No filtration	-	1.733191e-005 (25°C)	322	-	-	-	-	0.00000001 g/m3 inhalable aerosol	-	-	-	-	S
DIPHENYL OXIDE	C12H10O	101-84-8	Neutrodine Unisorb	2304	< 1 hPa	170	259	1 ppm	1 ppm	1 ppm	1 ppm	-	7mg/m ³	1 ppm	-	S
DIPROPYL KETONE	C7H14O	123-19-3	Neutrodine Unisorb	1530	0.164 kPa	114	173	50 ppm	50 ppm	-	-	-	-	-	-	S
DIPROPYLMETHANE	C7H16	142-82-5	Neutrodine Unisorb	1482	6.09 kPa	100	99	85 ppm	400 ppm	500 ppm	500 ppm	-	500 mg/m ³	500 ppm	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
DIURON	C9H10Cl2N2O	330-54-1	HEPA	-	2.6 10 ⁻⁷ Pa	233,1	-	10 mg/m ³ inhalable aerosol	10 mg/m ³	-	-	-	10 mg/m ³	-	-	
DIVINYL	C4H6	106-99-0	Neutrodine Unisorb	60	120 kPa (0°C)	54	-4,5	0,19 ppm	-	2 ppm	-	-	5 mg/m ³	10 ppm	-	S
DMA	C2H7N	124-40-3	Neutrodine Unisorb	75	203 kPa	45	7	10 ppm	1 ppm	2 ppm	2 ppm	-	5 mg/m ³	2 ppm	-	S
DMF	C3H7NO	68-12-2	Neutrodine Unisorb	1944	0.439 kPa	73	153	10 ppm	5 ppm	5 ppm	5 ppm	10 ppm	20 mg/m ³	-	-	
DMSO	C2H6SO	67-68-5	Neutrodine Unisorb	1170	0.1 kPa	78	189	-	-	-	50 ppm	-	-	-	-	
EAK	C8H16O	106-68-3	Neutrodine Unisorb	1116	0.286 kPa	128,21	167	-	-	-	-	-	-	-	-	S
EDTA	C10H16N2O8	60-00-4	HEPA	-	-	292,25	-	-	-	-	-	-	-	-	-	
EPICHLORHYDRINE	C3H5ClO	106-89-8	Neutrodine Unisorb	1248	2.20 kPa	93	115	-	-	2 ppm	-	-	1 mg/m ³	0,5 ppm	-	
EPOXY-2,3-PROPANOL-1	C3H6O2	556-52-5	Neutrodine Unisorb	2190	0.12 kPa	74	167	25 ppm	25 ppm	-	-	-	-	-	-	
ERYTHRENE	C4H6	106-99-0	Neutrodine Unisorb	60	120 kPa (0°C)	54	-4,5	0,19 ppm	-	2 ppm	-	-	5mg/m ³	10 ppm	-	S
ESSENCE OF MIRBANE	C6H5NO2	98-95-3	Neutrodine Unisorb	1284	0.04 kPa	123	210	1 ppm	0,2 ppm	1 mg/m ³	0,1 ppm	1 ppm	2 mg/m ³	0,2 ppm	0,2 ppm	
ETHANAL	C2H4O	75-07-0	Neutrodine Unisorb	66	120 kPa	44	20	18 ppm	100 ppm	50 ppm	50 ppm	50 ppm	-	-	-	S
ETHANEDIOIC ACID	C2H2O4	144-62-7	PF + HEPA or Neutrodine Unisorb	-	0.13 Pa	90,03	-	1 mg/m ³	1 mg/m ³	1 mg/m ³ inhalable aerosol	-	-	1 mg/m ³	1 mg/m ³	-	
ETHANOIC ACID	C2H4O2	64-19-7	Neutrodine Unisorb	5100	2.07 kPa	60	118	10 ppm	-	10 ppm	10 ppm	-	10 mg/m ³	-	-	S
ETHANOL	C2H6O	64-17-5	Neutrodine Unisorb	319	7.87 kPa	46	78	1000 ppm	1000 ppm	500 ppm	500 ppm	-	-	-	-	S
ETHANOLAMINE	C2H7NO	141-43-5	Neutrodine Unisorb	360	0.050 kPa	61	171	3 ppm	1 ppm	2 ppm	2 ppm	-	8 mg/m ³	1 ppm	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h		France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
ETHER	C4H10O	60-29-7	Neutrodine Unisorb	822	71.7 kPa	74	35	-		100 ppm	400 ppm	400 ppm	400 ppm	300 mg/m ³	-	-	S
ETHIDIUM BROMIDE	C21H20N3Br	1239-45-8	Neutrodine Unisorb	1560	-	314	238	-		-	-	-	-	-	-	-	-
ETHYL ACETATE	C4H8O2	141-78-6	Neutrodine Unisorb	2260	9.73 kPa	88	77	400 ppm		400 ppm	400 ppm	400 ppm	400 ppm	200 mg/m ³	200 ppm	-	S
ETHYL ACRYLATE	C5H8O2	140-88-5	Neutrodine Unisorb	1674	3.9 kPa	100	99,5	-		5 ppm	5 ppm	5 ppm	-	-	5 ppm	-	S
ETHYL ALCOHOL	C2H6O	64-17-5	Neutrodine Unisorb	319	7.87 kPa	46	78	1000 ppm		1000 ppm	500 ppm	500 ppm	-	-	-	-	S
ETHYL ALDEHYDE	C2H4O	75-07-0	Neutrodine Unisorb	66	120 kPa	44	20	18 ppm		100 ppm	50 ppm	50 ppm	50 ppm	-	-	-	S
ETHYL BENZENE	C8H10	100-41-4	Neutrodine Unisorb	1740	1.28 kPa	106	136	100ppm		20 ppm	20ppm	20ppm	20ppm	100mg/m ³	100 ppm	-	S
ETHYL BROMIDE	C2H5Br	74-96-4	Neutrodine Unisorb	1080	62.5 kPa	113	38,5	-		200 ppm	-	-	-	-	-	-	S
ETHYL CHLORIDE	C2H5Cl	75-00-3	Neutrodine Unisorb	48	160 kPa	65	12	-		100 ppm	40 ppm	-	-	-	-	-	S
ETHYL CYANOACRYLATE	C6H7NO2	7085-85-0	Neutrodine Unisorb	1248	0.27 kPa	125	66	-		-	-	-	-	-	0,3 ppm	-	-
ETHYL ETHANOATE	C4H8O2	141-78-6	Neutrodine Unisorb	2260	9.73 kPa	88	77	400 ppm		400 ppm	400 ppm	400 ppm	400 ppm	200 mg/m ³	200 ppm	-	S
ETHYL ETHER	C4H10O	60-29-7	Neutrodine Unisorb	822	71.7 kPa	74	35	-		100 ppm	400 ppm	400 ppm	400 ppm	300 mg/m ³	-	-	S
ETHYL FORMATE	C3H6O2	109-94-4	Neutrodine Unisorb	1032	32.3 kPa	74	35	100 ppm		100 ppm	100 ppm	100 ppm	-	-	100 ppm	-	S
ETHYL KETONE	C5H10O	96-22-0	Neutrodine Unisorb	1506	4.72 kPa	86	102	200 ppm		200 ppm	-	-	-	700 mg/m ³	200 ppm	5 ppm	S
ETHYL METHYL KETONE	C4H8O	78-93-3	Neutrodine Unisorb	984	12.6 kPa	72	80	200 ppm		200 ppm	200 ppm	200 ppm	200 ppm	300 mg/m ³	200 ppm	200 ppm	S
ETHYL NITRILE	C2H3N	75-05-8	Neutrodine Unisorb	*	11.9 kPa	41	82	20 ppm		40 ppm	20 ppm	20 ppm	-	30 mg/m ³	-	40 ppm	S
ETHYL OXIDE	C4H10O	60-29-7	Neutrodine Unisorb	822	71.7 kPa	74	35	-		100 ppm	400 ppm	400 ppm	400 ppm	300 mg/m ³	-	-	S

* More accurate estimation through eValiquet service

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
ETHYLAMINE	C2H7N	75-04-7	Neutrodine Unisorb	114	116 kPa	45	17	10 ppm	5 ppm	5 ppm	5 ppm	-	9 mg/m ³	-	-	S
ETHYLAMYL KETONE	C8H16O	106-68-3	Neutrodine Unisorb	1116	0.286 kPa	128,21	157	-	-	-	-	-	-	-	-	S
ETHYLEN CHLORHYDRIN	C2H5OCl	107-07-3	Neutrodine Unisorb	1440	4.45 kPa (50°C)	81	129	-	-	1 ppm	1 ppm	-	-	-	-	-
ETHYLENE ALCOHOL	C2H6O2	107-21-1	Neutrodine Unisorb	1116	0.010 kPa	66	198	-	20 ppm	10 ppm	10 ppm	-	20 mg/m ³	20 ppm	-	-
ETHYLENE BROMIDE	C2H4Br2	106-93-4	Neutrodine Unisorb	4680	1.55 kPa	188	131	-	-	-	-	-	-	0,5 ppm	-	S
ETHYLENE CHLORIDE	C2H4Cl2	107-06-2	Neutrodine Unisorb	1056	10.6 kPa	98	84	1 ppm	10 ppm	-	-	10ppm	7 mg/m ³	5 ppm	-	S
ETHYLENE CHLOROXYDRIN	C2H5OCl	107-07-3	Neutrodine Unisorb	1440	4.45 kPa (50°C)	81	129	-	-	1 ppm	1 ppm	-	-	-	-	-
ETHYLENE DIAMINE (SOLUTION)	C2H8N2	107-15-3	Neutrodine Unisorb	834	1.42 kPa	60	117	10 ppm	10 ppm	-	-	-	4 mg/m ³	10 ppm	-	-
ETHYLENE DIBROMIDE	C2H4Br2	106-93-4	Neutrodine Unisorb	4680	1.55 kPa	188	131	-	-	-	-	-	-	0,5 ppm	-	S
ETHYLENE DICHLORIDE	C2H4Cl2	107-06-2	Neutrodine Unisorb	1056	10.6 kPa	98	84	1 ppm	10 ppm	-	-	10ppm	7 mg/m ³	5 ppm	-	S
ETHYLENE GLYCOL	C2H6O2	107-21-1	Neutrodine Unisorb	1116	0.010 kPa	66	198	-	20 ppm	10 ppm	10 ppm	-	20 mg/m ³	20 ppm	-	S
ETHYLENE GLYCOL MONO ETHYL ETHER	C4H10O2	110-80-5	Neutrodine Unisorb	3182	0.71 kPa	90	135	0,5 ppm	2 ppm	2 ppm	2 ppm	5 ppm	18 mg/m ³	10 ppm	-	S
ETHYLENE GLYCOL MONO ETHYL ETHER ACETATE	C6H12O3	111-15-9	Neutrodine Unisorb	1518	0.24 kPa	132	157	0,5 ppm	2 ppm	2 ppm	2 ppm	5 ppm	30 mg/m ³	10 ppm	-	S
ETHYLENE TRICHLORIDE	C2HCl3	79-01-6	Neutrodine Unisorb	1956	9.91 kPa	130	86	25 ppm	75 ppm	6 ppm	-	25 ppm	30 mg/m ³	100 ppm	-	-
ETHYLENEDIAMINE	C2H8N2	107-15-3	Neutrodine Unisorb	834	1.42 kPa	60,1	118	10ppm	10 ppm	-	-	-	4 mg/m ³	10 ppm	-	-
ETHYLIDENE CHLORIDE	C2H4Cl2	75-34-3	Neutrodine Unisorb	456	30.5 kPa	98	57	100 ppm	100 ppm	100 ppm	100 ppm	100 ppm	-	-	100 ppm	S
ETHYNE	C2H2	74-86-2	Neutrodine Unisorb	48	4400 kPa	26	-84	2500 ppm	-	-	-	-	-	-	-	-

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapour pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
FERRIC OXIDE	Fe2O3	1309-37-1	PF + HEPA or Neutrodine Unisorb	-	-	159,7	-	5 mg/m ³ (total particulate)	-	-	-	-	-	5 mg/m ³	-	
FORMALDEHYDE	CH2O	50-00-0	Neutrodine Unisorb	258	220 kPa (100°C)	30	-	0,016 ppm	0,5 ppm	0,3 ppm	0,3 ppm	0,1 ppm	-	-	-	F
FORMALDEHYDE SOLUTION	CH2O	50-00-0	Neutrodine Unisorb	258	220 kPa (100°C)	30	-	0,016 ppm	0,5 ppm	0,3 ppm	0,3 ppm	0,1 ppm	-	-	-	F
FORMALIN (AS FORMALDEHYDE)	CH2O	50-00-0	Neutrodine Unisorb	258	220 kPa (100°C)	30	-	0,016 ppm	0,5 ppm	0,3 ppm	0,3 ppm	0,1 ppm	-	-	-	F
FORMAMIDE	CH2O2	75-12-7	Neutrodine Unisorb	1482	-	46	101	10 ppm	20 ppm	-	-	-	-	-	-	
FORMIC ACID	CH2O2	64-18-6	Neutrodine Unisorb	1482	5.75 kPa	46	101	5 ppm	5 ppm	5 ppm	5 ppm	-	10 mg/m ³	-	-	S
FORMIC ALDEHYDE	CH2O	50-00-0	Neutrodine Unisorb	258	220 kPa (100°C)	30	-	0,016 ppm	0,5 ppm	0,3 ppm	0,3 ppm	0,1 ppm	-	-	-	F
FORMONITRILE	HCN	74-90-8	Neutrodine Unisorb	270	82.7 kPa	27	26	-	2 ppm	-	1,9 ppm	5 ppm	-	-	-	
FURFURYL ALCOHOL	C5H6O2	98-00-0	Neutrodine Unisorb	1674	0.097 kPa	98	170	10 ppm	10 ppm	-	-	5 ppm	40 mg/m ³	-	-	S
FURYL CARBINOL	C5H6O2	98-00-0	Neutrodine Unisorb	1674	0.097 kPa	98	170	10 ppm	10 ppm	-	-	5 ppm	40 mg/m ³	-	-	S
GASOLINE 60	gasoline 60	8006-61-9	Neutrodine Unisorb	990	< 39 kPa	-	-	15 ppm LOQ	-	-	-	-	-	-	-	S
GLACIAL ACETIC ACID (PURE COMPOUND)	C2H4O2	64-19-7	Neutrodine Unisorb	5100	2.07 kPa	60	118	10 ppm	-	10 ppm	10 ppm	-	10 mg/m ³	-	-	S
GLUCOSE	C6H12O6	01-10-5996	PF + HEPA or Neutrodine Unisorb	-	-	180,16	-	-	-	-	-	-	-	-	-	
GLUTARALDEHYDE	C5H8O2	111-30-8	Neutrodine Unisorb	384	2.27 kPa	100	187	-	0,1 ppm	0,05 ppm	0,05 ppm	-	-	0,05 ppm	-	S
GLYCEROL , MIST	C3H8O3	56-81-5	HEPA	-	0.39 Pa (20°C)	92,09	182	-	10 mg/m ³	-	50 mg/m ³ inhalable aerosol	-	-	-	-	
GLYCID	C3H6O2	556-52-5	Neutrodine Unisorb	2190	0.12 kPa	74	167	25 ppm	25 ppm	-	-	-	-	-	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapour pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
GLYCIDOL	C3H6O2	556-52-5	Neutrodine Unisorb	2190	0.12 kPa	74	167	25 ppm	25 ppm	-	-	-	-	-	-	-
GLYCOL	C2H6O2	107-21-1	Neutrodine Unisorb	1116	0.010 kPa	66	198	-	20 ppm	10 ppm	10 ppm	-	20 mg/m ³	20 ppm	-	S
GRAPHITE (SYNTHETIC)	C	7440-44-0	PF + HEPA or Neutrodine Unisorb	-	-	12,01	3825	-	2 mg/m ³ respirable aerosol	-	4 mg/m ³ inhalable aerosol	-	-	-	-	-
HALOTHANE	C2HBrClF3	151-67-7	Neutrodine Unisorb	780	32.5 kPa	197	50	2 ppm	-	5 ppm	5 ppm	-	-	10 ppm	-	-
HEPTAN-4-ONE	C7H14O	123-19-3	Neutrodine Unisorb	1530	0.164 kPa	114	173	50 ppm	50 ppm	-	-	-	-	-	-	S
HEPTANE	C7H16	142-82-5	Neutrodine Unisorb	1482	6.09 kPa	100	99	85 ppm	400 ppm	500 ppm	500 ppm	-	500 mg/m ³	500 ppm	-	-
HEXANE	C6H14	110-54-3	Neutrodine Unisorb	1913	20.2 kPa	86	69	50 ppm	20 ppm	50 ppm	50 ppm	50 ppm	100 mg/m ³	20 ppm	-	S
HEXANOL	C6H14O	111-27-3	Neutrodine Unisorb	2755	0,1 kPa	102,2	157	-	-	50 ppm	-	-	-	-	-	-
HEXONE	C6H12O	108-10-1	Neutrodine Unisorb	1572	2.64 kPa	100	116	50 ppm	20 ppm	20 ppm	20 ppm	50ppm	-	50 ppm	-	S
HEXYL ALCOHOL	C6H14O	111-27-3	Neutrodine Unisorb	2755	0,1 kPa	102,2	157	-	-	50 ppm	-	-	-	-	-	-
HYDRAZINE	N2H4	302-01-2	Neutrodine Unisorb	768	1.3 kPa	32,05	-	-	0,1 ppm	0,017 ppm	-	-	0,06 mg/m ³	0,02 ppm	-	-
HYDROBROMIC ACID	HBr	10035-10-6	Neutrodine Unisorb	1626	0.15 kPa	81	-66	-	-	6,7 mg/m ³	2 ppm	-	10mg/m ³	3 ppm	-	A
HYDROCHLORIC ACID	HCl aq. sol.	7647-01-0	Neutrodine Unisorb	2184	4103 kPa	37	120	-	-	2 ppm	2 ppm	-	-	1 ppm	5 ppm	A
HYDROCYANIC ACID	HCN	74-90-8	Neutrodine Unisorb	270	82.7 kPa	27	26	-	2 ppm	-	1,9 ppm	5 ppm	-	-	-	-
HYDROFLUORIC ACID	HF aq. sol.	7664-39-3	Neutrodine Unisorb	540	104 kPa	20	112	3 ppm	1,8 ppm	1 ppm	1 ppm	0,5 ppm	-	1,8 ppm	1,8 ppm	-
HYDROGEN BROMIDE	HBr	10035-10-6	Neutrodine Unisorb	1626	0.15 kPa	81	-66	-	-	6,7 mg/m ³	2 ppm	-	10mg/m ³	3 ppm	-	A

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapour pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
HYDROGEN CHLORIDE	HCl	7647-01-0	Neutrodine Unisorb	2184	4103 kPa	37	-85	-	-	2 ppm	2 ppm	-	-	1 ppm	5 ppm	A
HYDROGEN CYANIDE	HCN	74-90-8	Neutrodine Unisorb	270	82.7 kPa	27	26	-	2 ppm	-	1,9 ppm	5 ppm	-	-	-	-
HYDROGEN DIOXIDE	H2O2	7722-84-1	Neutrodine Unisorb	1854	1.32 kPa	34	158	1 ppm	1 ppm	-	0,5 ppm	-	1,5 mg/m ³	1 ppm	-	-
HYDROGEN FLUORIDE	HF	7664-39-3	Neutrodine Unisorb	540	104 kPa	20	20	3 ppm	1,8 ppm	1 ppm	1 ppm	0,5 ppm	-	1,8 ppm	1,8 ppm	-
HYDROGEN NITRATE	HNO3	7697-37-2	Neutrodine Unisorb	1608	6.39 kPa	63	120	2 ppm	-	-	-	2 ppm	-	1 ppm	-	A
HYDROGEN PEROXIDE	H2O2	7722-84-1	Neutrodine Unisorb	1854	1.32 kPa	34	158	1 ppm	1 ppm	-	0,5 ppm	-	1,5 mg/m ³	1 ppm	-	-
HYDROGEN SULFATE	H2SO4	7664-93-9	Neutrodine Unisorb	1674	1.3 Pa	98	296	1 mg/m ³	0,05 mg/m ³ Thoracic fraction	0,1 mg/m ³ inhalable aerosol	0,1 mg/m ³ inhalable aerosol	1 ppm	1 mg/m ³	-	0,05 mg/m ³	-
HYDROGEN SULFIDE	H2S	04-06-7783	Neutrodine Unisorb	810	1783 kPa	34	-60	-	5 ppm	5 ppm	5 ppm	5 ppm	-	5 ppm	5 ppm	-
(subject to risk assessment)																
HYDROQUINONE	C6H6O2	123-31-9	Neutrodine Unisorb	1878	1.3 Pa	110,11	285	-	2 mg/m ³	-	-	-	1 mg/m ³	0,5 mg/m ³	-	-
HYDROXYBENZENE	C6H6O	108-95-2	Neutrodine Unisorb	1290	0.055 kPa	94	182	5 ppm	2 ppm	2 ppm	-	-	10 mg/m ³	2 ppm	-	S
HYDROXYCELLULOSE	(C6H10O5) _n	9004-34-6	PF + HEPA or Neutrodine Unisorb	-	-	160,000-560,000	-	10 mg/m ³ total dust	10 mg/m ³ inhalable aerosol	-	-	-	10 mg/m ³	10 mg/m ³ inhalable aerosol	-	-
HYDROXYCYCLOHEXANE	C6H12O	108-93-0	Neutrodine Unisorb	1830	0.1 kPa	100	162	50 ppm	50 ppm	50 ppm	50 ppm	25 ppm	100 mg/m ³	50 ppm	-	S
HYPOCHLOROUS ACID	HClO	7790-92-3	Neutrodine Unisorb	1608	-	52,5	100	-	-	-	-	-	-	-	-	-
IODINE	I2	7553-56-2	Neutrodine Unisorb	2700	10 kPa (9°C)	254	185	-	-	-	-	0,1 ppm	-	-	-	-
IODOFORM	CHI3	75-47-8	PF + HEPA or Neutrodine Unisorb	-	-	393,73	210	0,6 ppm	0,6 ppm	-	-	-	10 mg/m ³	-	-	-

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapour pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
IPA	C3H8O	67-63-0	Neutrodine Unisorb	1421	6.02 kPa	60	83	400 ppm	-	200 ppm	200 ppm	400 ppm	350 mg/m ³	-	-	S
IRON OXIDE DUST AND FUME (AS FE)	Fe2O3	1309-37-1	PF + HEPA or Neutrodine Unisorb	-	-	159,7	-	5 mg/m ³ (total particulate)	-	-	-	-	-	5 mg/m ³	-	
ISOAMYL ACETATE	C7H14O2	123-92-2	Neutrodine Unisorb	1812	0.728 kPa	130	142	100 ppm	50 ppm	50 ppm	50 ppm	100 ppm	-	50 ppm	-	S
ISOAMYL ALCOHOL	C5H12O	71-41-0	Neutrodine Unisorb	3416	0.259 kPa	88	132	-	-	-	20 ppm	-	100 mg/m ³	-	-	S
ISOAMYL ALCOHOL (PRIMARY)	C5H12O	123-51-3	Neutrodine Unisorb	1542	0.315 kPa	88	132	100 ppm	100 ppm	-	20 ppm	100 ppm	-	100 ppm	-	S
ISOBUTANE	C4H10	75-28-5	Neutrodine Unisorb	60	348 kPa	58	-12	-	-	1000 ppm	1000 ppm	-	-	-	-	
ISOBUTANOL	C4H10O	78-83-1	Neutrodine Unisorb	1542	1.39 kPa	74	108	50 ppm	50 ppm	100 ppm	100 ppm	50 ppm	-	50 ppm	-	S
ISOBUTENYL METHYL KETONE	C6H10O	141-79-7	Neutrodine Unisorb	2178	1.47 kPa	98	130	10 ppm	15 ppm	-	5 ppm	-	60 mg/m ³	15 ppm	-	S
ISOBUTYL ACETATE	C6H12O2	110-19-0	Neutrodine Unisorb	1740	2.39 kPa	116	117	150 ppm	150 ppm	62 ppm	100 ppm	150 ppm	-	150 ppm	-	S
ISOBUTYL ALCOHOL	C4H10O	78-83-1	Neutrodine Unisorb	1542	1.39 kPa	74	108	50 ppm	50 ppm	100 ppm	100 ppm	50 ppm	-	50 ppm	-	S
ISOBUTYL CARBINOL	C5H12O	123-51-3	Neutrodine Unisorb	1542	0.315 kPa	88	132	100 ppm	100 ppm	-	20 ppm	100 ppm	-	100 ppm	-	S
ISOBUTYL METHYL CARBINOL	C6H14O	108-11-2	Neutrodine Unisorb	1572	0.39 kPa	102	132	25 ppm	25 ppm	20 ppm	20 ppm	-	-	25 ppm	-	S
ISOBUTYRONE	C7H14O	565-80-0	Neutrodine Unisorb	1458	6.87 kPa	114	124	-	-	-	-	-	-	-	-	
ISOFLURANE	C3H2F5ClO	26675-46-7	Neutrodine Unisorb	780	34.9 kPa (22°C)	184,5	48,5	-	-	-	-	-	-	50 ppm	-	
ISO-NITROPROPANE	C3H7NO2	79-46-9	Neutrodine Unisorb	1044	2.3 kPa	89	120	-	-	0,05 ppm	-	-	30 mg/m ³	5 ppm	-	S
ISOOCTANE	C8H18	540-84-1	Neutrodine Unisorb	1488	6.5 kPa	114	99	-	-	-	-	-	-	-	-	S
ISOOCTANOL	C8H18O	104-76-7	Neutrodine Unisorb	1500	0.019 kPa	130	190	-	-	10 ppm	10 ppm	-	-	-	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
ISOOCTYLALCOHOL	C8H18O	104-76-7	Neutrodine Unisorb	1500	0.019 kPa	130	190	-	-	10 ppm	10 ppm	-	-	-	-	
ISOPENTANE	C5H12	78-78-4	Neutrodine Unisorb	678	91.7 kPa	72	28	120 ppm	-	1000 ppm	1000 ppm	-	500 mg/m ³	-	1000 ppm	S
ISOPENTYL ACETATE	C7H14O2	123-92-2	Neutrodine Unisorb	1812	0.728 kPa	130	142	100 ppm	50 ppm	50 ppm	50 ppm	100 ppm	-	50 ppm	-	S
ISOPHORONE	C9H14O	78-59-1	Neutrodine Unisorb	1674	0.039 kPa	138	215	4 ppm	-	2 ppm	2 ppm	-	-	-	-	
ISOPRENE	C5H8	78-79-5	Neutrodine Unisorb	768	73.4 kPa	68	34	-	-	3 ppm	3 ppm	-	-	-	-	S
ISOPROPANOL	C3H8O	67-63-0	Neutrodine Unisorb	1421	6.02 kPa	60	83	400 ppm	-	200 ppm	200 ppm	400 ppm	350 mg/m ³	-	-	S
ISOPROPYL ACETATE	C5H10O2	108-21-4	Neutrodine Unisorb	1572	5.59 kPa	102	88	-	250 ppm	-	100 ppm	250 ppm	-	-	-	S
ISOPROPYL ALCOHOL	C3H8O	67-63-0	Neutrodine Unisorb	1421	6.02 kPa	60	83	400 ppm	-	200 ppm	200 ppm	400 ppm	350 mg/m ³	-	-	S
ISOPROPYL BENZENE	C9H12	98-82-8	Neutrodine Unisorb	1776	0.61 kPa	120	152	50 ppm	20 ppm	10 ppm	10 ppm	-	-	25 ppm	20 ppm	S
ISOPROPYL ETHER	C6H14O	108-20-3	Neutrodine Unisorb	2486	8.35 kPa	102	69	-	250 ppm	200 ppm	200 ppm	-	-	250 ppm	-	
ISOPROPYL GLYCIDYL ETHER	C6H12O2	4016-14-2	Neutrodine Unisorb	1788	1.2 kPa	116	127	-	50 ppm	-	-	-	-	50 ppm	-	
ISOPROPYLAMINE	C3H9N	75-31-0	Neutrodine Unisorb	234	78 kPa	59	34	-	5 ppm	5 ppm	5 ppm	-	12 mg/m ³	-	-	
ISOPROPYLCARBINOL	C4H10O	78-83-1	Neutrodine Unisorb	1542	1.39 kPa	74	108	50 ppm	50 ppm	100 ppm	100 ppm	50 ppm	-	50 ppm	-	S
ISOPROPYLIDENEACETONE	C6H10O	141-79-7	Neutrodine Unisorb	2178	1.47 kPa	98	130	10 ppm	15 ppm	-	5 ppm	-	60 mg3/m3	15 ppm	-	S
KORAX	C3H6ClNO2	600-25-9	Neutrodine Unisorb	2250	0.79 kPa	123	171	2 ppm	2 ppm	-	-	-	-	-	-	S
LIMONENE	C10H16	5989-54-8	Neutrodine Unisorb	1668	0.4 kPa (14.4°C)	136	178	-	-	-	-	-	-	-	-	
LITHIUM HYDRIDE	LiH	7580-67-8	HEPA	-	-	7,95	-	0,025 mg/m ³	0,025 mg/m ³	0,025 mg/m ³ inhalable aerosol	-	-	0,025mg/m ³	0,025mg/m ³	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
MAGNESIA FUME	MgO	1309-48-4	PF + HEPA or Neutrodine Unisorb	-	-	40,31	3600	-	10 mg/m ³ respirable aerosol	-	4 inhalable aerosol mg/m ³	-	10 mg/m ³	"10 mg/m ³ inhalable aerosol 4 mg/m ³ respirable aerosol"	-	
MAGNESITE	MgCO ₃	546-93-0	PF + HEPA or Neutrodine Unisorb	-	-	84	-	10 mg/m ³ total dust	10 mg/m ³ respirable aerosol	-	-	-	-	-	-	
MAGNESIUM OXIDE FUME	MgO	1309-48-4	PF + HEPA or Neutrodine Unisorb	-	-	40,31	3600	-	10 mg/m ³ respirable aerosol	-	4 inhalable aerosol mg/m ³	-	10 mg/m ³	"10 mg/m ³ inhalable aerosol 4 mg/m ³ respirable aerosol"	-	
MANGANESE COMPOUNDS (AS MN)	Mn	7439-96-5	PF + HEPA or Neutrodine Unisorb	-	-	55	1962	1 mg/m ³	-	0,5 mg/m ³ inhalable aerosol	0,2 mg/m ³ inhalable aerosol	1 mg/m ³	-	-	-	
MANGANESE OXIDE	MnO ₂	1344-43-0	PF + HEPA or Neutrodine Unisorb	-	-	86,94	-	-	1 mg/m ³	-	-	-	-	-	-	
M-CRESOL	C ₇ H ₈ O	108-39-4	Neutrodine Unisorb	1578	0.019 kPa	108,14	203	2,3 ppm	-	-	-	-	-	5 ppm	-	
MEK	C ₄ H ₈ O	78-93-3	Neutrodine Unisorb	984	12.6 kPa	72	80	200 ppm	200 ppm	200 ppm	200 ppm	200 ppm	300 mg/m ³	200 ppm	200 ppm	S
MERCAPTO-2 ETHANOL	C ₂ H ₆ SO	60-24-2	Neutrodine Unisorb	1170	0.1 kPa	78	157	-	-	-	-	-	-	-	-	
MERCURY	Hg	7439-97-6	No filtration	-	0.16 Pa	200,59	356	-	0,02 mg/m ³	0,02 mg/m ³	0,02 mg/m ³ inhalable aerosol	0,05 mg/m ³	-	0,025 mg/m ³	0,02 mg/m ³	
MESITYL OXIDE	C ₆ H ₁₀ O	141-79-7	Neutrodine Unisorb	2178	1.47 kPa	98	130	10 ppm	15 ppm	-	5 ppm	-	60 mg ₃ /m ₃	15 ppm	-	S
MESITYLENE	C ₉ H ₁₂	108-67-8	Neutrodine Unisorb	1776	16.6 kPa	120	152	-	20 ppm	20 ppm	20 ppm	-	-	25 ppm	-	S
METALLIC MERCURY	Hg	7439-97-6	No filtration	-	0.16 Pa	200,59	356,73	0,05 mg/m ³	-	0,02 mg/m ³	0,02 mg/m ³ inhalable aerosol	-	0,02 mg/m ³	-	-	
METHACRYLIC ACID	C ₄ H ₆ O ₂	79-41-4	Neutrodine Unisorb	2280	0.703 kPa (50°C)	86	161	20 ppm	20 ppm	50 ppm	50 ppm	2 ppm	3 mg/m ³	20 ppm	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
METHANAL	CH2O	50-00-0	Neutrodine Unisorb	258	220 kPa (100°C)	30	-	0,016 ppm	0,5 ppm	0,3 ppm	0,3 ppm	0,1 ppm	-	-	-	F
METHANE TETRABROMIDE	CBr4	558-13-4	Neutrodine Unisorb	2700	5.3 kPa	332	189,5	0,1 ppm	0,1 ppm	-	-	-	1,5 mg/m ³	-	-	
METHANOIC ACID	CH2O2	64-18-6	Neutrodine Unisorb	1482	5.75 kPa	46	101	5 ppm	5 ppm	5 ppm	5 ppm	-	10 mg/m ³	-	-	S
METHANOL	CH4O	67-56-1	Neutrodine Unisorb	*	16.9 kPa	32	65	200 ppm	200 ppm	200 ppm	200 ppm	200 ppm	25 mg/m ³	-	-	S
METHOXYCARBONYLETHYLENE	C4H5O2	96-33-3	Neutrodine Unisorb	924	5.3 kPa	86	80,5	10 ppm	5 ppm	5 ppm	2 ppm	2 ppm	20 mg/m ³	-	-	
METHYL ACETATE	C3H6O2	79-20-9	Neutrodine Unisorb	1246	28.8 kPa	74	58	200 ppm	200 ppm	200 ppm	100 ppm	200 ppm	200 mg/m ³	200 ppm	-	S
METHYL ACETONE	C4H8O	78-93-3	Neutrodine Unisorb	984	12.6 kPa	72	80	200 ppm	200 ppm	200 ppm	200 ppm	200 ppm	300 mg/m ³	200 ppm	200 ppm	S
METHYL ACETYLENE	C3H4	74-99-7	Neutrodine Unisorb	24	145 kPa (-25°C)	40	-23	1000 ppm	1000 ppm	-	-	-	-	-	-	
METHYL ALCOHOL	CH4O	67-56-1	Neutrodine Unisorb	*	16.9 kPa	32	65	200 ppm	200 ppm	200 ppm	200 ppm	200 ppm	25 mg/m ³	-	-	S
METHYL ALDEHYDE	CH2O	50-00-0	Neutrodine Unisorb	258	220 kPa (100°C)	30	-	0,016 ppm	0,5 ppm	0,3 ppm	0,3 ppm	0,1 ppm	-	-	-	F
METHYL BENZENE	C7H8	108-88-3	Neutrodine Unisorb	2784	3.79 kPa	92	110	100 ppm	20 ppm	50 ppm	50 ppm	20 ppm	50 mg/m ³	50 ppm	-	S
METHYL BUTYL KETONE	C6H12O	591-78-6	Neutrodine Unisorb	1488	1.54 kPa	100	127	1 ppm	5 ppm	5 ppm	5 ppm	5 ppm	20 mg/m ³	-	-	S
METHYL CELLOSOLVE	C3H8O2	109-86-4	Neutrodine Unisorb	1632	0.79 kPa	76	125	0,1 ppm	1 ppm	1 ppm	1 ppm	5 ppm	15 mg/m ³	-	-	S
METHYL CHLORIDE	CH3Cl	74-87-3	Neutrodine Unisorb	18	575 kPa	51	-24	-	50 ppm	50 ppm	50 ppm	-	60 mg/m ³	-	-	S
METHYL CHLOROFORM	C2H3Cl3	71-55-6	Neutrodine Unisorb	1080	16.5 kPa	133	74	-	100 ppm	200 ppm	200 ppm	200 ppm	900 mg/m ³	-	-	
METHYL CYANIDE	C2H3N	75-05-8	Neutrodine Unisorb	*	11.9 kPa	41	82	20 ppm	40 ppm	20 ppm	20 ppm	-	30 mg/m ³	-	40 ppm	S
METHYL CYANOACRYLATE	C5H5NO2	137-05-3	Neutrodine Unisorb	456	0.27 kPa	111	49	2 ppm	2 ppm	2 ppm	2 ppm	-	-	-	-	

* More accurate estimation through eValiquest service

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
METHYL CYCLOHEXANE	C7H14	108-87-2	Neutrodine Unisorb	2241	6.18 kPa	98	100	400 ppm	400 ppm	200 ppm	200 ppm	-	-	196 ppm	-	S
METHYL CYCLOHEXANOL	C7H14O	25639-42-3	Neutrodine Unisorb	1884	0.13 kPa	114	155	50 ppm	50 ppm	6 ppm	-	-	-	50 ppm	-	S
METHYL CYCLOHEXANONE	C7H12O	1331-22-2	Neutrodine Unisorb	1938	0.13 kPa	112	165	50 ppm	50 ppm	-	-	50 ppm	-	-	-	S
METHYL ETHER	C2H6O	115-10-6	Neutrodine Unisorb	18	273 kPa (0°C)	46	-25	-	1000 ppm	1000 ppm	1000 ppm	-	-	400 ppm	-	S
METHYL ETHYL KETONE	C4H8O	78-93-3	Neutrodine Unisorb	984	12.6 kPa	72	80	200 ppm	200 ppm	200 ppm	200 ppm	200 ppm	300 mg/m ³	200 ppm	200 ppm	S
METHYL FORMATE	C2H4O2	107-31-3	Neutrodine Unisorb	108	78.1 kPa	60	32	100 ppm	100 ppm	50 ppm	50 ppm	-	-	100 ppm	-	S
METHYL ISOBUTENYL KETONE	C6H10O	141-79-7	Neutrodine Unisorb	2178	1.47 kPa	98	130	10 ppm	15 ppm	-	5 ppm	-	60 mg3/m3	15 ppm	-	S
METHYL ISOBUTYL KETONE	C6H12O	108-10-1	Neutrodine Unisorb	1572	2.64 kPa	100	116	50 ppm	20 ppm	20 ppm	20 ppm	50ppm	-	50 ppm	-	S
METHYL METACRYLATE	C5H8O2	80-62-6	Neutrodine Unisorb	1674	5.1 kPa	100	101	100 ppm	50 ppm	50 ppm	50 ppm	8,3 mg/m ³	100 mg/m ³	50 ppm	50 ppm	S
METHYL PHENOL ALL ISOMERS	C7H8O	1319-77-3	Neutrodine Unisorb	1578	0.018 kPa	108	191	2,3 ppm	5 ppm	-	-	5 ppm	10 mg/m ³	5 ppm	-	S
METHYL PROPENOATE	C4H5O2	96-33-3	Neutrodine Unisorb	924	5.3 kPa	86	80,5	10 ppm	5 ppm	5 ppm	2 ppm	2 ppm	20 mg/m ³	-	-	S
METHYL PROPYL KETONE	C5H10O	107-87-9	Neutrodine Unisorb	1500	4.97 kPa	86	102	150 ppm	200 ppm	-	-	-	-	200 ppm	-	S
METHYL STYRENE	C9H10	25013-15-4	Neutrodine Unisorb	1740	0.13 kPa	118	152	100 ppm	50 ppm	100 ppm	100 ppm	-	-	100 ppm	-	S
METHYL TERT-BUTYL ETHER	C5H12O	1634-04-4	Neutrodine Unisorb	1869	27 kPa	88,2	55,05	-	50 ppm	50 ppm	50 ppm	-	-	25 ppm	50 ppm	S
METHYL-2-PROPANE	C4H10	75-28-5	Neutrodine Unisorb	60	349 kPa	58	-12	-	-	1000 ppm	1000 ppm	-	-	-	-	S
METHYL-2-PROPANOL-2	C4H10O	75-65-0	Neutrodine Unisorb	1170	5.42 kPa	74	83	100 ppm	100 ppm	20 ppm	20 ppm	50 ppm	-	100 ppm	-	S

* More accurate estimation through eValiquest service

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h		France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
METHYL-3-BUTANOL-1	C5H12O	123-51-3	Neutrodine Unisorb	3416	0.259 kPa	88	132	-		-	-	20 ppm	-	100 mg/m ³	-	-	S
METHYLACRYLATE	C4H5O2	96-33-3	Neutrodine Unisorb	924	5.3 kPa	86	80,5	10 ppm		5 ppm	5 ppm	2 ppm	2 ppm	20 mg/m ³	-	-	
METHYLAL	C3H8O2	109-87-5	Neutrodine Unisorb	744	53.1 kPa	76	43	1000 ppm		1000 ppm	1000 ppm	1000 ppm	-	-	1000 ppm	-	S
METHYLAMINE	CH5N	74-89-5	Neutrodine Unisorb	36	353 kPa	31	-7	10 ppm		-	10 ppm	10 ppm	-	5 mg/m ³	-	-	
METHYLAMYL ALCOHOL	C6H14O	108-11-2	Neutrodine Unisorb	1572	0.39 kPa	102	132	25 ppm		25 ppm	20 ppm	20 ppm	-	-	25 ppm	-	S
METHYLENE CHLORIDE	CH2Cl2	75-09-2	Neutrodine Unisorb	604	58.2 kPa	85	40	-		50 ppm	50 ppm	50 ppm	100 ppm	200 mg/m ³	-	100 ppm	
METHYLENE CHLOROBROMIDE	CH2BrCl	74-97-5	Neutrodine Unisorb	1992	19.5 kPa	129	68	200 ppm		200 ppm	-	-	-	-	-	-	
METHYLENE DICHLORIDE	CH2Cl2	75-09-2	Neutrodine Unisorb	604	58.2 kPa	85	40	-		50 ppm	50 ppm	50 ppm	100 ppm	200 mg/m ³	-	100 ppm	
METHYLENE OXIDE	CH2O	50-00-0	Neutrodine Unisorb	258	220 kPa (100°C)	30	-	0,016 ppm		0,5 ppm	0,3 ppm	0,3 ppm	0,1 ppm	-	-	-	F
METHYLETHYL CARBINOL	C4H10O	78-92-2	Neutrodine Unisorb	1392	2.32 kPa	74	99,5	100 ppm		100 ppm	-	-	103 ppm	-	100 ppm	-	S
METHYL-N-AMYL KETONE	C7H14O	110-43-0	Neutrodine Unisorb	1620	-	114	151	100 ppm		50 ppm	238 mg/m ³	-	-	-	50 ppm	-	S
MIBC	C6H14O	108-11-2	Neutrodine Unisorb	1572	0.39 kPa	102	132	25 ppm		25 ppm	20 ppm	20 ppm	-	-	25 ppm	-	S
MIBK	C6H12O	108-10-1	Neutrodine Unisorb	1572	2.64 kPa	100	116	50 ppm		20 ppm	20 ppm	20 ppm	50ppm	-	50 ppm	-	S
MICA (CONTAINING LESS THAN 1% QUARTZ)	K(Mg, Fe)3AlSi3O10(F, OH)2	12001-26-2	PF + HEPA or Neutrodine Unisorb	-	-	797	-	3 mg/m ³ respirable fraction		-	-	-	-	"2 mg/m ³ inhalable fraction 1.5 mg/m ³ respirable fraction"	"10mg/m ³ inhalable fraction 0.8 mg/m ³ respirable fraction"	-	
MIRBANE OIL	C6H5NO2	98-95-3	Neutrodine Unisorb	1284	0.04 kPa	123	210	1 ppm		0,2 ppm	1 mg/m ³	0,1 ppm	1 ppm	2 mg/m ³	0,2 ppm	0,2 ppm	
MTBE	C5H12O	1634-04-4	Neutrodine Unisorb	1869	27 kPa	88,2	55,05	-		50 ppm	50 ppm	50 ppm	-	-	25 ppm	50 ppm	

BY ALPHABETICAL ORDER

M

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
MURIATIC ACID	HCl aq. sol.	7647-01-0	Neutrodine Unisorb	2184	4103 kPa	37	120	-	-	2 ppm	2 ppm	-	-	1 ppm	5 ppm	A
MUSCOVITE	K(Mg, Fe)3AlSi3O10(F, OH)2	12001-26-2	PF + HEPA or Neutrodine Unisorb	-	-	797	-	3 mg/m ³ respirable fraction	-	-	-	-	"2 mg/m ³ inhalable fraction 1.5 mg/m ³ respirable fraction"	"10 mg/m ³ inhalable fraction 0.8 mg/m ³ respirable fraction"	-	-
N,N-DIMETHYLETHYLAMINE	C4H11N	598-56-1	Neutrodine Unisorb	420	-	73,14	36	-	5 ppm	-	2 ppm	-	-	-	-	S
N-AMYL ACETATE	C7H14O2	123-92-2	Neutrodine Unisorb	1614	0.728 kPa	130	148	100 ppm	50 ppm	50 ppm	50 ppm	100 ppm	-	50 ppm	-	S
N-AMYL ACETATE	C7H14O2	628-63-7	Neutrodine Unisorb	1614	0.6 kPa	130,18	142	100 ppm	50 ppm	50 ppm	50 ppm	100 ppm	100 mg/m ³	-	-	-
NAPHTA 30/60	85% Nonane/15% trimethylbenzene	8052-41-3	Neutrodine Unisorb	768	-	-	220	350 mg/m ³	-	-	-	-	-	-	-	S
NAPHTHALENE	C10H8	91-20-3	PF + HEPA or Neutrodine Unisorb	-	7.2 Pa	128,2	-	10 ppm	10 ppm	0,1 ppm	-	-	50 mg/m ³	-	10 ppm	-
NAPHTHALIN	C10H8	91-20-3	PF + HEPA or Neutrodine Unisorb	-	7.2 Pa	128,2	-	10 ppm	10 ppm	0,1 ppm	-	-	50 mg/m ³	-	10 ppm	-
N-BUTANE	C4H10	106-97-8	Neutrodine Unisorb	60	242 kPa	58	-12	800 ppm	800 ppm	1000 ppm	1000 ppm	-	-	600 ppm	-	S
N-BUTANETHIOL	C4H10S	109-79-5	Neutrodine Unisorb	156	6.07 kPa	90	97	-	0,5 ppm	0,5 ppm	0,5 ppm	-	2 mg/m ³	-	-	S
N-BUTANOL	C4H10O	71-36-3	Neutrodine Unisorb	1760	0.86 kPa	74	117,5	-	-	100 ppm	100 ppm	25 ppm	100 mg/m ³	-	-	S
N-BUTYL ACETATE	C6H12O2	123-86-4	Neutrodine Unisorb	3808	1.66 kPa	116	127	150 ppm	150 ppm	62 ppm	100 ppm	150 ppm	200 mg/m ³	150 ppm	-	S
N-BUTYL AMINE	C4H9NH2	109-73-9	Neutrodine Unisorb	456	12.2 kPa	73	78	-	-	-	2 ppm	-	-	-	-	S
N-BUTYL CHLORIDE	C4H9Cl	109-69-3	Neutrodine Unisorb	1308	13.7 kPa	92	78,5	-	-	25 ppm	-	-	-	-	-	-

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
N-BUTYL MERCAPTAN	C4H10S	109-79-5	Neutrodine Unisorb	156	6.07 kPa	90	97	-	0,5 ppm	0,5 ppm	0,5 ppm	-	2 mg/m ³	-	-	S
N-ETHYLETHANAMINE	C4H11N	109-89-7	Neutrodine Unisorb	378	30.1 kPa	73	56	10 ppm	5 ppm	5 ppm	5 ppm	-	-	5 ppm	-	S
N-HEXANE	C6H14	110-54-3	Neutrodine Unisorb	1913	20.2 kPa	86	69	50 ppm	20 ppm	50 ppm	50 ppm	50 ppm	100 mg/m ³	20 ppm	-	S
N-HEXANOL	C6H14O	111-27-3	Neutrodine Unisorb	2755	0,1 kPa	102,2	157	-	-	50 ppm	-	-	-	-	-	-
NICKEL METAL AND OTHER COMPOUNDS (AS NI)	Ni	7440-02-0	PF + HEPA or Neutrodine Unisorb	-	-	58,69	2732	0,015 mg/m ³	1 mg/m ³	-	-	-	1 mg/m ³	-	-	-
NINHYDRIN (POWDER)	C9H4O3 .H2O	485-47-2	HEPA	-	-	178,14	240	-	-	-	-	-	-	-	-	-
NITRIC ACID	HNO3	7697-37-2	Neutrodine Unisorb	1608	6.39 kPa	63	120	2 ppm	-	-	-	2 ppm	-	1 ppm	-	A
NITRO BENZENE	C6H5NO2	98-95-3	Neutrodine Unisorb	1284	0.04 kPa	123	210	1 ppm	0,2 ppm	1 mg/m ³	0,1 ppm	1 ppm	2 mg/m ³	0,2 ppm	0,2 ppm	-
NITROETHANE	C2H5NO2	79-24-3	Neutrodine Unisorb	1170	2.79 kPa	75	114	100 ppm	100 ppm	100 ppm	10 ppm	-	300 mg/m ³	-	20 ppm	S
NITROGEN DIOXIDE	NO2	10102-44-0	No filtration	-	101 kPa	46	-	-	-	-	0,5ppm	-	5mg/m ³	3 ppm	0,5 ppm	-
NITROGLYCERINE	C3H5N3O9	55-63-0	No filtration	-	0.03 Pa	227,1	-	-	0,1 ppm	0,01 ppm	0,01 ppm	-	-	-	-	-
NITROMETHANE	CH3NO2	75-52-5	Neutrodine Unisorb	972	4.79 kPa	61	101	-	100 ppm	-	-	-	50 mg/m ³	-	-	S
NITROPROPANE 2	C3H7NO2	79-46-9	Neutrodine Unisorb	1044	2.3 kPa	89	120	-	-	0,55 ppm	-	-	30 mg/m ³	5 ppm	-	S
N-OCTANE	C8H18	111-65-9	Neutrodine Unisorb	1188	1.86 kPa	114	126	75 ppm	300 ppm	500 ppm	500 ppm	-	500 mg/m ³	210 ppm	-	S
NONANE ALL ISOMERS	C9H20	111-84-2	Neutrodine Unisorb	1380	0.57 kPa	128	151	200 ppm	200 ppm	-	-	-	500 mg/m ³	222 ppm	-	S
N-PENTANE	C5H12	109-66-0	Neutrodine Unisorb	920	68.3 kPa	72	36	120 ppm	1000 ppm	1000 ppm	1000 ppm	-	-	600 ppm	-	S
ORTHOPHOSPHORIC ACID	H3PO4	7664-38-2	Neutrodine Unisorb	1674	0.004 kPa	98	276	1 mg/m ³	0,2 ppm	2 mg/m ³ inhalable aerosol	2 mg/m ³ inhalable aerosol	1 mg/m ³	1 ppm	1 mg/m ³	1 mg/m ³	-
OSMIUM TETROXIDE (AS OS)	OsO4	20816-12-0	No filtration	-	0.93 kPa	-	-	0,0002 ppm	0,0002 ppm	-	-	-	-	0,0002 ppm	-	-

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
O-TOLUIDINE	C7H9N	119-93-7	Neutrodine Unisorb	465	0.039 kPa	107,15	199	-	-	-	-	-	-	-	-	-
OXALIC ACID	C2H2O4	144-62-7	PF + HEPA or Neutrodine Unisorb	-	0.13 Pa	90,03	-	1 mg/m ³	1 mg/m ³	1 mg/m ³ inhalable aerosol	-	-	1 mg/m ³	1 mg/m ³	-	-
O-XYLENE	C8H10	95-47-6	Neutrodine Unisorb	1920	0.88 kPa	106	138	100 ppm	50 ppm	100 ppm	100 ppm	100 ppm	50 mg/m ³	-	50 ppm	S
OZONE	O3	10028-15-6	No filtration	-	101 kPa	48	-	-	0,1 ppm	-	-	-	-	-	-	-
PARAQUAT	C12H14N2	4685-14-7	HEPA	-	-	186,26	-	0,1 mg/m ³	0,1 mg/m ³	-	-	-	0,5 mg/m ³ respirable aerosol	-	-	-
PARATHION	C10H14NO5PS	56-38-2	HEPA	-	0.005 kPa	291,26	-	0,05 mg/m ³	0,1 mg/m ³	0,1 mg/m ³ inhalable aerosol	0,1 mg/m ³ inhalable aerosol	-	0,05 mg/m ³	-	-	-
P-CRESOL	C7H8O	106-44-5	Neutrodine Unisorb	1578	0.017 kPa	108,14	202	2,3 ppm	-	-	-	-	-	5 ppm	-	-
P-DICHLOROBENZENE	C6H4Cl2	106-46-7	Neutrodine Unisorb	3060	0.003 kPa	147	173	-	0,75 ppm	1 ppm	-	-	30mg/m ³	25 ppm	-	S
PENTACHLOROETHANE	C2H5Cl5	76-01-7	Neutrodine Unisorb	3672	0.39 kPa	200	161	-	-	-	5 ppm	-	-	-	-	-
PENTANOL 1	C5H12O	71-41-0	Neutrodine Unisorb	3416	0.259 kPa	88	138	-	-	-	20 ppm	-	100 mg/m ³	-	-	S
PENTYL ACETATE	C7H14O2	628-63-7	Neutrodine Unisorb	1614	0.6 kPa	130,18	142	100 ppm	50 ppm	50 ppm	50 ppm	100 ppm	100 mg/m ³	-	-	-
PERCHLORIC ACID	HClO4	7601-90-3	Neutrodine Unisorb	1674	-	100	203	-	-	-	-	-	-	-	-	-
PERCHLOROETHYLENE	C2Cl4	127-18-4	Neutrodine Unisorb	2796	2.42 kPa	166	121	-	20 ppm	20 ppm	-	50 ppm	200 mg/m ³	50 ppm	-	-
PETROLEUM ETHER 30/60	ether de petrole 30/60	8032-32-4	Neutrodine Unisorb	768	5.3 kPa	75	-	350 mg/m ³	-	-	-	-	-	-	-	S
PHENOL	C6H6O	108-95-2	Neutrodine Unisorb	1290	0.055 kPa	94	182	5 ppm	2 ppm	2 ppm	-	-	10 mg/m ³	2 ppm	-	S
PHENOLPHTALEINE	C20H14O4	77-09-8	HEPA	-	-	318,32	-	-	-	-	-	-	-	-	-	-
PHENYL AMINE	C6H5NH2	62-53-3	Neutrodine Unisorb	1944	0.09 kPa	93	184	-	2 ppm	2 ppm	2 ppm	-	3 mg/m ³	-	-	-
PHENYL CHLORIDE	C6H5Cl	108-90-7	Neutrodine Unisorb	2364	1.6 kPa	113	133	-	5 ppm	10 ppm	10 ppm	10 ppm	50 mg/m ³	1 ppm	-	S

BY ALPHABETICAL ORDER
N-P



Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h		France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
PHENYL ETHER	C12H10O	101-84-8	Neutrodine Unisorb	2304	< 1 hPa	170	259	1 ppm		1 ppm	1 ppm	1 ppm	-	7mg/m ³	1 ppm	-	
PHENYL GLYCIDYL ETHER	C9H10O2	122-60-1	Neutrodine Unisorb	2214	1.3 Pa	150	245	-		1 ppm	-	-	-	-	1 ppm	-	
PHENYL HYDROXIDE	C6H6O	108-95-2	Neutrodine Unisorb	1290	0.055 kPa	94	182	5 ppm		2 ppm	2 ppm	-	-	10 mg/m ³	2 ppm	-	S
PHENYLETHANE	C8H10	100-41-4	Neutrodine Unisorb	1740	1.28 kPa	106	136	100ppm		20 ppm	20ppm	20ppm	20ppm	100mg/m ³	100 ppm	-	S
PHENYLETHYLENE	C8H8	100-42-5	Neutrodine Unisorb	1260	0.81 kPa	104	146	50ppm		50 ppm	20 ppm	20 ppm	50ppm	50mg/m ³	100 ppm	-	S
PHENYLMETHANE	C7H8	108-88-3	Neutrodine Unisorb	2784	3.79 kPa	92	110	100 ppm		20 ppm	50 ppm	50 ppm	20 ppm	50 mg/m ³	50 ppm	-	S
PHOSPHORIC ACID	H3PO4	7664-38-2	Neutrodine Unisorb	1674	0.004 kPa	98	276	1 mg/m ³		0,2 ppm	2 mg/m ³ inhalable aerosol	2 mg/m ³ inhalable aerosol	1 mg/m ³	1 ppm	1 mg/m ³	1 mg/m ³	
PLATINUM	Pt	04-06-7440	PF + HEPA or Neutrodine Unisorb	-	-	195,8	4300	1 mg/m ³		1 mg/m ³	1 mg/m ³ inhalable aerosol	-	-	-	-	-	
POTASSIUM HYDRATE	KOH	1310-58-3	PF + HEPA or Neutrodine Unisorb	-	0.13 kPa (1044°C)	56,11	-	-		-	-	-	-	-	-	-	
POTASSIUM HYDROXIDE	KOH	1310-58-3	PF + HEPA or Neutrodine Unisorb	-	0.13 kPa (1044°C)	56,11	-	-		-	-	-	-	-	-	-	
PROPANAL	C3H6O	123-38-6	Neutrodine Unisorb	597	42.2 kPa	58,08	46	-		-	-	-	-	-	-	-	
PROPANE	C3H8	74-98-6	Neutrodine Unisorb	12	939 kPa	44	-42	1000 ppm		-	1000 ppm	1000 ppm	-	-	-	-	S
PROPANOL-1	C3H8O	71-23-8	Neutrodine Unisorb	834	2.76 kPa	60	97	200 ppm		200 ppm	-	-	-	200 mg/m ³	-	-	S
PROPANAMINE	C3H9NO	107-10-8	Neutrodine Unisorb	384	42.1 kPa	75	180	-		-	-	-	-	-	-	-	
PROPARGYL ALCOHOL	C3H4O	107-19-7	Neutrodine Unisorb	750	1.59 kPa	56	113	1 ppm		1 ppm	2 ppm	2 ppm	-	-	1 ppm	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
PROPENENITRILE	C3H3N	107-13-1	Neutrodine Unisorb	486	11.3 kPa	53	77	1 ppm	2 ppm	1,2 ppm	-	2 ppm	1 mg/m ³	2 ppm	-	
PROPYL ACETATE	C5H10O2	109-60-4	Neutrodine Unisorb	3396	4.49 kPa	102	102	200 ppm	200 ppm	-	100 ppm	200 ppm	200 mg/m ³	200 ppm	-	S
PROPYL ALCOHOL	C3H8O	71-23-8	Neutrodine Unisorb	834	2.76 kPa	60	97	200 ppm	200 ppm	-	-	-	200 mg/m ³	-	-	S
PROPYL MERCAPTAN	C3H8S	107-03-9	Neutrodine Unisorb	78	20.6 kPa	76,2	67	0,3 ppm	-	-	-	-	-	-	-	
PROPYLACETATE	C5H10O2	109-60-4	Neutrodine Unisorb	3396	4.49 kPa	102,13	102	200 ppm	200 ppm	-	100 ppm	200 ppm	200 mg/m ³	200 ppm	-	S
PROPYLAMINE	C3H9N	107-10-8	Neutrodine Unisorb	384	42.1 kPa	59	48	-	-	-	-	-	-	-	-	
PROPYLBENZENE	C9H12	108-67-8	Neutrodine Unisorb	1776	16.6 kPa	120	152	-	20 ppm	20 ppm	20 ppm	-	-	25 ppm	-	S
PROPYLENE ALDEHYDE	C4H6O	4170-30-3	Neutrodine Unisorb	990	4.92 kPa	70	102	-	2 ppm	-	-	-	-	-	-	S
PROPYLENE CHLORIDE	C3H6Cl2	78-87-5	Neutrodine Unisorb	1350	6.62 kPa	113	97	-	75 ppm	-	-	1 ppm	350 mg/m ³	-	-	S
PROPYLENE DICHLORIDE	C3H6Cl2	78-87-5	Neutrodine Unisorb	1350	6.62 kPa	113	97	-	75 ppm	-	-	1 ppm	350 mg/m ³	-	-	S
PROPYLENE OXIDE	C3H6O	75-56-9	Neutrodine Unisorb	192	59.3 kPa	58	34	-	20 ppm	2 ppm	2 ppm	-	5 mg/m ³	5 ppm	-	S
PROPYNE	C3H4	74-99-7	Neutrodine Unisorb	24	145 kPa (-25°C)	40	-23	1000 ppm	1000 ppm	-	-	-	-	-	-	
PRUSSIC ACID	HCN	74-90-8	Neutrodine Unisorb	270	82.7 kPa	27	26	-	2 ppm	-	1,9 ppm	5 ppm	-	-	-	
P-TOLUIDINE	C7H9N	106-49-0	Neutrodine Unisorb	465	1.74 kPa	107,2	200	-	-	-	-	-	-	-	-	
PYRIDINE	C5H5N	110-86-1	Neutrodine Unisorb	960	2.13 kPa	79	115	5 ppm	5 ppm	-	-	-	4 mg/m ³	5 ppm	-	S
PYROCELLULOSE	(C6H10O5) _n	9004-34-6	PF + HEPA or Neutrodine Unisorb	-	-	160,000-560,000	-	10 mg/m ³ total dust	10 mg/m ³ inhalable aerosol	-	-	-	10 mg/m ³	10 mg/m ³ inhalable aerosol	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h		France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
QUARTZ	SiO ₂	14808-60-7	PF + HEPA or Neutrodine Unisorb	-	-	60,1	2230	30/(%silica+2) mg/m ³ total dust		0,1 mg/m ³ respirable aerosol	-	-	-	1 mg/m ³ respirable fraction	-	-	
RED IRON OXIDE	Fe ₂ O ₃	1309-37-1	PF + HEPA or Neutrodine Unisorb	-	-	159,7	-	5 mg/m ³ (total particulate)		-	-	-	-	-	5 mg/m ³	-	
SACCAROSE	C ₁₂ H ₂₂ O ₁₁	57-50-1	PF + HEPA or Neutrodine Unisorb	-	-	342,3	-	"10 mg/m ³ total dust 5 mg/m ³ respirable fraction"		10 mg/m ³	-	-	-	-	-	-	
SBA	C ₄ H ₁₀ O	78-92-2	Neutrodine Unisorb	1392	2.32 kPa	74	99,5	100 ppm		100 ppm	-	-	104 ppm	-	100 ppm	-	S
SEC-AMYL ACETATE	C ₇ H ₁₄ O ₂	123-92-2	Neutrodine Unisorb	1830	0.728 kPa	130	123	100 ppm		50 ppm	50 ppm	50 ppm	100 ppm	-	50 ppm	-	S
SEC-BUTYL AMINE	C ₄ H ₉ NH ₂	13952-84-6	Neutrodine Unisorb	420	23 kPa	73	63	-		-	-	2 ppm	-	-	-	-	S
SILICA GEL	SiO ₂	7631-86-9	HEPA	-	-	60,1	2230	-		-	4 mg/m ³ inhalable aerosol	4 mg/m ³ inhalable aerosol	-	2 mg/m ³ inhalable fraction	-	-	
SILICA, AMORPHOUS	SiO ₂	7631-86-9	HEPA	-	-	60,1	2230	-		-	4 mg/m ³ inhalable aerosol	4 mg/m ³ inhalable aerosol	-	2 mg/m ³ inhalable fraction	-	-	
SILICON	Si	7440-21-3	HEPA	-	-	28,1	2355	-		10 mg/m ³ respirable aerosol	-	-	-	-	-	-	
SILVER (DUST)	Ag	7440-22-4	HEPA	-	0.34 Pa (961 °C)	107,87	2000	-		-	0,01 mg/m ³ inhalable aerosol	0,01 mg/m ³ inhalable aerosol	-	-	-	-	
SODIUM BISULPHITE	HNaO ₃ S	7631-90-5	PF + HEPA or Neutrodine Unisorb	-	-	104,06	-	5 mg/m ³		5 mg/m ³	-	-	-	-	5 mg/m ³	-	
SODIUM HYDRATE	NaOH	1310-73-2	PF + HEPA or Neutrodine Unisorb	-	-	40	1390	-		2 mg/m ³	-	-	-	-	-	-	
SODIUM HYDROXIDE	NaOH	1310-73-2	PF + HEPA or Neutrodine Unisorb	-	-	40	1390	-		2 mg/m ³	-	-	-	-	-	-	

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
STODDARD SOLVENT	85% No-nane/15% trimethylbenzene	8052-41-3	Neutrodine Unisorb	768	-	-	220	350 mg/m ³	-	-	-	-	-	-	-	S
STRYCHNINE	C ₂₁ H ₂₂ N ₂ O ₂	57-24-9	HEPA	-	-	334,42	-	0,15 mg/m ³	0,15 mg/m ³	-	-	-	-	-	-	-
STYRENE	C ₈ H ₈	100-42-5	Neutrodine Unisorb	1260	0.81 kPa	104	146	50ppm	50 ppm	20 ppm	20 ppm	50ppm	50mg/m ³	100 ppm	-	S
SULFURIC ACID	H ₂ SO ₄	7664-93-9	Neutrodine Unisorb	1674	1.3 Pa	98	296	1 mg/m ³	0,05 mg/m ³ thoracic fraction	0,1 mg/m ³ inhalable aerosol	0,1 mg/m ³ inhalable aerosol	1 ppm	1 mg/m ³	-	0,05 mg/m ³	-
TERT-BUTYL ACETATE	C ₆ H ₁₂ O ₂	540-88-5	Neutrodine Unisorb	1884	-	116	96	200 ppm	200 ppm	42 ppm	50 ppm	-	-	-	-	S
TERT-BUTYL ALCOHOL	C ₄ H ₁₀ O	75-65-0	Neutrodine Unisorb	1170	5.42 kPa	74	83	100 ppm	100 ppm	20 ppm	20 ppm	50 ppm	-	100 ppm	-	S
TERT-BUTYL CHLORIDE	C ₄ H ₉ Cl	507-20-0	Neutrodine Unisorb	1128	21.2 kPa	96	68	-	-	-	-	-	-	-	-	-
TERT-BUTYL METHYL ETHER	C ₅ H ₁₂ O	1634-04-4	Neutrodine Unisorb	1869	27 kPa	88,2	55,05	-	50 ppm	50 ppm	50 ppm	-	-	25 ppm	50 ppm	-
TETRABROMOMETHANE	CBr ₄	558-13-4	Neutrodine Unisorb	2700	5.3 kPa	332	189,5	0,1 ppm	0,1 ppm	-	-	-	1,5 mg/m ³	-	-	-
TETRACHLOROETHYLENE	C ₂ Cl ₄	127-18-4	Neutrodine Unisorb	2796	2.42 kPa	166	121	-	20 ppm	20 ppm	-	50 ppm	200 mg/m ³	50 ppm	-	-
TETRACHLOROMETHANE	CCl ₄	56-23-5	Neutrodine Unisorb	2790	15.2 kPa	154	77	-	2 ppm	0,5 ppm	0,5 ppm	-	15 mg/m ³	-	-	-
TETRAHYDROFURAN	C ₄ H ₈ O	109-99-9	Neutrodine Unisorb	1044	21.6 kPa	72	65	200 ppm	50 ppm	50 ppm	50 ppm	50 ppm	300 mg/m ³	50 ppm	-	S
TFA	C ₂ HF ₃ O ₂	76-05-1	Neutrodine Unisorb	1080	15.1kPa	114	72	-	-	-	-	-	-	-	-	-
THF	C ₄ H ₈ O	109-99-9	Neutrodine Unisorb	1044	21.6 kPa	72	65	200 ppm	50 ppm	50 ppm	50 ppm	50 ppm	300 mg/m ³	50 ppm	-	S
TIN (INORGANIC COMPOUNDS, AS SN)	Sn	7440-31-5	PF + HEPA or Neutrodine Unisorb	-	-	118,69	2260	-	-	-	-	-	-	-	-	-
TIN(IV) OXIDE (AS SN)	O ₂ Sn	18282-10-5	PF + HEPA or Neutrodine Unisorb	-	-	150,69	-	-	-	-	-	-	-	-	-	-

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
TITANIUM DIOXIDE	TiO2	13463-67-7	HEPA	-	-	79,9	2900	-	11 mg/m ³ inhalable aerosol	-	-	-	8 mg/m ³ inhalable fraction	"10 mg/m ³ inhalable aerosol 4 mg/m ³ respirable aerosol"	-	
TMA	C3H9N	75-50-3	Neutrodine Unisorb	48	215 kPa	59	-4	10 ppm	-	-	2 ppm	-	-	-	-	S
TOLUENE	C7H8	108-88-3	Neutrodine Unisorb	2784	3.79 kPa	92	110	100 ppm	20 ppm	50 ppm	50 ppm	20 ppm	50 mg/m ³	50 ppm	-	S
TOLUOL	C7H8	108-88-3	Neutrodine Unisorb	2784	3.79 kPa	92	110	100 ppm	20 ppm	50 ppm	50 ppm	20 ppm	50 mg/m ³	50 ppm	-	S
TRIBROMOMETHANE	CHBr3	75-25-2	Neutrodine Unisorb	900	0.726 kPa	253	149,5	0,5 ppm	0,5 ppm	-	-	1 ppm	-	-	-	S
TRICHLOROACETIC ACID	C2HCl3O2	76-03-9	Neutrodine Unisorb	2700	0.101 kPa (50°C)	163	198	1 ppm	1 ppm	-	0,2 ppm	-	-	-	-	
TRICHLOROETHANE-1,1,2	C2H3Cl3	79-00-5	Neutrodine Unisorb	1860	3.1 kPa	132	114	10 ppm	-	10 ppm	10 ppm	10 ppm	-	-	-	
TRICHLOROETHANOIC ACID	C2HCl3O2	76-03-9	Neutrodine Unisorb	2700	0.101 kPa (50°C)	163,39	197,55	1 ppm	1 ppm	-	0,2 ppm	-	-	-	-	
TRICHLOROETHENE	C2HCl3	79-01-6	Neutrodine Unisorb	1956	9.91 kPa	130	86	25 ppm	75 ppm	6 ppm	-	25 ppm	30 mg/m ³	100 ppm	-	
TRICHLOROETHYLENE	C2HCl3	79-01-6	Neutrodine Unisorb	1956	9.91 kPa	130	86	25 ppm	75 ppm	6 ppm	-	25 ppm	30 mg/m ³	100 ppm	-	
TRICHLOROMETHANE	CHCl3	67-66-3	Neutrodine Unisorb	780	26.2 kPa	119	61	-	2 ppm	0,5 ppm	0,5 ppm	3 ppm	20 mg/m ³	-	-	S
TRIETHYLAMINE	C6H15N	121-44-8	Neutrodine Unisorb	246	7.7 kPa	101	90	-	1 ppm	1 ppm	1 ppm	-	-	2 ppm	-	S
TRIFLUOROACETIC ACID	C2HF3O2	76-05-1	Neutrodine Unisorb	1080	15.1 kPa	114	72	-	-	-	-	-	-	-	-	
TRIMETHYL CARBINOL	C4H10O	75-65-0	Neutrodine Unisorb	1170	5.42 kPa	74	83	100ppm	100 ppm	20 ppm	20 ppm	50 ppm	-	100 ppm	-	S
TRIMETHYL METHANE	C4H10	75-28-5	Neutrodine Unisorb	60	350 kPa	58	-12	-	-	1000 ppm	1000 ppm	-	-	-	-	
TRIMETHYL PENTANE-2,2,4	C8H18	540-84-1	Neutrodine Unisorb	1488	6.5 kPa	114	99	-	-	-	-	-	-	-	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h		France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
TRIMETHYLAMINE	C3H9N	75-50-3	Neutrodine Unisorb	48	215 kPa	59	-4	10 ppm		-	-	2 ppm	-	-	-	-	S
TRIMETHYLBENZENE	C9H12	108-67-8	Neutrodine Unisorb	1776	16.6 kPa	120	152	-		20 ppm	20 ppm	20 ppm	-	-	25 ppm	-	S
TRINITROGLYCERINE	C3H5N3O9	55-63-0	No filtration	-	0.001 kPa (100°C)	227,1	-	-		0,1 ppm	0,01 ppm	0,01 ppm	-	-	-	-	
TUNGSTEN insoluble	W	7440-33-7	PF + HEPA or Neutrodine Unisorb	-	-	183,84	5555	5 mg/m³		-	-	-	-	5 mg/m³	-	-	
TURPENTINE OIL	C10H16	8006-64-2	Neutrodine Unisorb	1422	0.53 kPa	-	160	100 ppm		100 ppm	-	5 ppm	50 ppm	300 mg/m³	100 ppm	-	
UREA	CH4N2O	57-13-6	PF + HEPA or Neutrodine Unisorb	-	-	60,06	-	-		-	-	-	-	-	-	-	
VC	C2H3Cl	75-01-4	Neutrodine Unisorb	48	355 kPa	61	-14	-		1 ppm	3 ppm	-	2 ppm	10 mg/m³	-	-	S
VERT DE BROMOCRESOL	C21H14Br4O5S	76-60-8	PF + HEPA or Neutrodine Unisorb	-	-	698,01	-	-		-	-	-	-	-	-	-	
VINYL ACETATE	C4H6O2	108-05-4	Neutrodine Unisorb	924	15.4 kPa	86	73	-		5 ppm	5 ppm	-	-	10mg/m³	10 ppm	-	S
VINYL BROMIDE	C2H3Br	593-60-2	Neutrodine Unisorb	48	141 kPa	107	16	-		-	-	-	-	-	-	-	
VINYL CARBINOL	C3H6O	107-18-6	Neutrodine Unisorb	678	3.14 kPa	58	97	2 ppm		0,2 ppm	2 ppm	-	-	2 mg/m³	2 ppm	-	S
VINYL CHLORIDE	C2H3Cl	75-01-4	Neutrodine Unisorb	48	355 kPa	61	-14	-		1 ppm	3 ppm	-	2 ppm	10 mg/m³	-	-	S
VINYL CYANIDE	C3H3N	107-13-1	Neutrodine Unisorb	486	11.3 kPa	53	77	1 ppm		2 ppm	1,2 ppm	-	2 ppm	1 mg/m³	2 ppm	-	
VINYL ETHYLENE	C4H6	106-99-0	Neutrodine Unisorb	60	120 kPa (0°C)	54	-4,5	0,19 ppm		-	2 ppm	-	-	5 mg/m³	10 ppm	-	S

Chemical name	Formula	CAS number	Suitable filter	Retention capacity (g)	Vapor pressure	MM (g/mol)	Boiling point (°C)	NIOSH 8h	France 8h	AGS 8h	DFG 8h	Japan 8h	China 8h	UK 8h	European union 8h	Detection Sensor
VINYL TOLUENE	C9H10	25013-15-4	Neutrodine Unisorb	1740	0.13 kPa	118	170	100 ppm	50 ppm	100 ppm	100 ppm	-	-	100 ppm	-	S
VINYL TRICHLORIDE	C2H3Cl3	79-00-5	Neutrodine Unisorb	1860	3.1 kPa	132	114	10 ppm	-	10 ppm	10 ppm	10 ppm	-	-	-	
VINYLBENZENE	C8H8	100-42-5	Neutrodine Unisorb	1260	0.81 kPa	104	146	50ppm	50 ppm	20 ppm	20 ppm	50ppm	50mg/m ³	100 ppm	-	S
WHITE SPIRIT	85% Nonane / 15% trimethylbenzene	8052-41-3	Neutrodine Unisorb	768	-	-	220	350 mg/m ³	-	-	-	-	-	-	-	S
XYLENE	C8H10	1330-20-7	Neutrodine Unisorb	3360	1.15 kPa	106	138	100 ppm	50 ppm	100 ppm	100 ppm	100 ppm	50 mg/m ³	50 ppm	-	S
XYLENE (ISOMERS)	C8H10	1330-20-7	Neutrodine Unisorb	3360	1.15 kPa	106	138	100 ppm	50 ppm	100 ppm	100 ppm	100 ppm	50 mg/m ³	50 ppm	-	S
ZINC OXIDE	ZnO	1314-13-2	PF + HEPA or Neutrodine Unisorb	-	-	81,38	-	5 mg/m ³	5 mg/m ³	-	-	-	-	5 mg/m ³	-	

List of chemicals substances by formula

Formula	Chemical name
Ag	SILVER (DUST)
Al	ALUMINIUM
Al ₂ O ₃	ALUMINA
Al ₂ O ₃	ALUMINUM OXIDE
Al ₂ O ₃	ALUMINUM TRIOXIDE
Al ₂ Si ₂ O ₅ (OH) ₄	CLAY
As	ARSENIC (INORGANIC COMPOUNDS, AS AS)
BF ₃	BORON TRIFLUORIDE
B ₂ O ₃	BORON OXIDE
BaCl ₂ .2H ₂ O	BARIUM CHLORIDE
Be	BERYLLIUM COMPOUNDS (AS BE)
Br ₂	BROMINE
C	CARBON BLACK
C	GRAPHITE (SYNTHETIC)
CBr ₄	CARBON BROMIDE
CBr ₄	CARBON TETRABROMIDE
CBr ₄	METHANE TETRABROMIDE
CBr ₄	TETRABROMOMETHANE
CCl ₄	CARBON TETRACHLORIDE
CCl ₄	TETRACHLOROMETHANE
CHBr ₃	TRIBROMOMETHANE
CHBr ₃	BROMOFORM
CHCl ₃	CHLOROFORM
CHCl ₃	TRICHLOROMETHANE

Formula	Chemical name
CHI ₃	IODOFORM
CH ₂ BrCl	BROMOCHLOROMETHANE
CH ₂ BrCl	CHLOROBROMOMETHANE
CH ₂ BrCl	METHYLENE CHLOROBROMIDE
CH ₂ Cl ₂	DICHLOROMETHANE
CH ₂ Cl ₂	METHYLENE CHLORIDE
CH ₂ Cl ₂	METHYLENE DICHLORIDE
CH ₂ O	FORMALDEHYDE
CH ₂ O	FORMALDEHYDE SOLUTION
CH ₂ O	FORMALIN (AS FORMALDEHYDE)
CH ₂ O	FORMIC ALDEHYDE
CH ₂ O	METHANAL
CH ₂ O	METHYL ALDEHYDE
CH ₂ O	METHYLENE OXIDE
CH ₂ O ₂	FORMAMIDE
CH ₂ O ₂	FORMIC ACID
CH ₂ O ₂	METHANOIC ACID
CH ₃ Cl	CHLOROMETHANE
CH ₃ Cl	METHYL CHLORIDE
CH ₃ NO ₂	NITROMETHANE
CH ₄ N ₂ O	UREA
CH ₄ O	METHANOL
CH ₄ O	METHYL ALCOHOL
CH ₅ N	AMINOMETHANE

Formula	Chemical name
CH ₅ N	METHYLAMINE
C ₂ Cl ₄	PERCHLOROETHYLENE
C ₂ Cl ₄	TETRACHLOROETHYLENE
C ₂ HBrClF ₃	HALOTHANE
C ₂ HCl ₃	ETHYLENE TRICHLORIDE
C ₂ HCl ₃	TRICHLOROETHENE
C ₂ HCl ₃	TRICHLOROETHYLENE
C ₂ HCl ₃ O ₂	TRICHLOROACETIC ACID
C ₂ HCl ₃ O ₂	TRICHLOROETHANOIC ACID
C ₂ HF ₃ O ₂	TFA
C ₂ HF ₃ O ₂	TRIFLUOROACETIC ACID
C ₂ H ₂	ACETYLENE
C ₂ H ₂	ETHYNE
C ₂ H ₂ Br ₄	1,1,2,2-TETRABROMOETHANE
C ₂ H ₂ Br ₄	ACETYLENE TETRABROMIDE
C ₂ H ₂ Cl ₂	1,2-DICHLOROETHYLENE
C ₂ H ₂ Cl ₂	ACETYLENE DICHLORIDE
C ₂ H ₂ Cl ₄	ACETYLENE TETRACHLORIDE
C ₂ H ₂ Cl ₄	1,1,2,2-TETRACHLOROETHANE
C ₂ H ₂ O ₄	ETHANEDIOIC ACID
C ₂ H ₂ O ₄	OXALIC ACID
C ₂ H ₃ Br	BROMOETHENE
C ₂ H ₃ Br	BROMOETHYLENE
C ₂ H ₃ Br	VINYL BROMIDE
C ₂ H ₃ Cl	CHLOROETHENE
C ₂ H ₃ Cl	CHLOROETHYLENE
C ₂ H ₃ Cl	VC
C ₂ H ₃ Cl	VINYL CHLORIDE

Formula	Chemical name
C ₂ H ₃ Cl ₃	1,1,1-TRICHLOROETHANE
C ₂ H ₃ Cl ₃	CHLOROTHENE
C ₂ H ₃ Cl ₃	METHYL CHLOROFORM
C ₂ H ₃ Cl ₃	TRICHLOROETHANE-1,1,2
C ₂ H ₃ Cl ₃	VINYL TRICHLORIDE
C ₂ H ₃ N	ACETONITRILE
C ₂ H ₃ N	CYANOMETHANE
C ₂ H ₃ N	ETHYL NITRILE
C ₂ H ₃ N	METHYL CYANIDE
C ₂ H ₃ OCl	2-CHLOROACETALDEHYDE
C ₂ H ₃ OCl	2-CHLOROETHANAL
C ₂ H ₄ Br ₂	1,2-DIBROMOETHANE
C ₂ H ₄ Br ₂	ETHYLENE BROMIDE
C ₂ H ₄ Br ₂	ETHYLENE DIBROMIDE
C ₂ H ₄ Cl ₂	1,1-DICHLOROETHANE
C ₂ H ₄ Cl ₂	1,2-DICHLOROETHANE
C ₂ H ₄ Cl ₂	ETHYLENE CHLORIDE
C ₂ H ₄ Cl ₂	ETHYLENE DICHLORIDE
C ₂ H ₄ Cl ₂	ETHYLIDENE CHLORIDE
C ₂ H ₄ O	ACETALDEHYDE
C ₂ H ₄ O	ETHANAL
C ₂ H ₄ O	ETHYL ALDEHYDE
C ₂ H ₄ O ₂	ACETIC ACID
C ₂ H ₄ O ₂	ETHANOIC ACID
C ₂ H ₄ O ₂	GLACIAL ACETIC ACID (PURE COMPOUND)
C ₂ H ₄ O ₂	METHYL FORMATE
C ₂ H ₅ Br	BROMOETHANE

Formula	Chemical name
C2H5Br	ETHYL BROMIDE
C2H5Cl	CHLOROETHANE
C2H5Cl	ETHYL CHLORIDE
C2H5Cl5	PENTACHLOROETHANE
C2H5NO2	NITROETHANE
C2H5OCl	ETHYLEN CHLORHYDRIN
C2H5OCl	ETHYLENE CHLOROXYDRIN
C2H5OCl	2-CHLOROETHANOL
C2H5OCl	2-CHLOROETHYL ALCOHOL
C2H6O	ABSOLUTE ALCOHOL
C2H6O	ALCOHOL
C2H6O	DIMETHYL ETHER
C2H6O	ETHANOL
C2H6O	ETHYL ALCOHOL
C2H6O	METHYL ETHER
C2H6O2	1,2-ETHANEDIOL
C2H6O2	5-METHYL-3-HEPTANONE
C2H6O2	ETHYLENE ALCOHOL
C2H6O2	ETHYLENE GLYCOL
C2H6O2	GLYCOL
C2H6SO	DIMETHYL SULFOXIDE
C2H6SO	DMSO
C2H6SO	MERCAPTO-2 ETHANOL
C2H7N	AMINOETHANE
C2H7N	DIMETHYL AMINE
C2H7N	DMA
C2H7N	ETHYLAMINE
C2H7NO	2-AMINOETHANOL

Formula	Chemical name
C2H7NO	BETA-AMINOETHYL ALCOHOL
C2H7NO	ETHANOLAMINE
C2H8N2	ETHYLENE DIAMINE (SOLUTION)
C2H8N2	ETHYLENEDIAMINE
C3H2F5ClO	ISOFLURANE
C3H3N	2-PROPENENITRILE
C3H3N	ACRYLONITRILE
C3H3N	PROPENENITRILE
C3H3N	VINYL CYANIDE
C3H4	ALLYLENE
C3H4	METHYL ACETYLENE
C3H4	PROPYNE
C3H4Cl2	1,3-DICHLOROPROPENE
C3H4Cl2	1,3-DICHLOROPROPYLENE
C3H4O	2-PROPENAL
C3H4O	2-PROPYN-1-OL
C3H4O	2-PROPYNYL ALCOHOL
C3H4O	ACROLEIN
C3H4O	ACRYLIC ALDEHYDE
C3H4O	ALLYL ALDEHYDE
C3H4O	PROPARGYL ALCOHOL
C3H4O2	2-PROPENOIC ACID
C3H4O2	ACROLEIC ACID
C3H4O2	ACRYLIC ACID
C3H5Cl	ALLYL CHLORIDE
C3H5Cl	3-CHLORO-1-PROPENE
C3H5ClO	1-CHLORO-2,3-EPOXYPROPANE
C3H5ClO	EPICHLORHYDRINE

Formula	Chemical name
C3H5ClO	2-CHLOROPROPYLENE OXIDE
C3H5NO	2-PROPENAMIDE
C3H5NO	ACRYLAMIDE
C3H5N3O9	NITROGLYCERINE
C3H5N3O9	TRINITROGLYCERINE
C3H6ClNO2	CHLORO-1-NITROPROPANE 1
C3H6ClNO2	KORAX
C3H6Cl2	DICHLOROPROPANE 1, 2
C3H6Cl2	PROPYLENE CHLORIDE
C3H6Cl2	PROPYLENE DICHLORIDE
C3H6O	2-PROPANONE
C3H6O	2-PROPEN-1-OL
C3H6O	2-PROPENOL
C3H6O	ACETONE
C3H6O	ALLYL ALCOHOL
C3H6O	ALLYLIC ALCOHOL
C3H6O	DIMETHYL KETONE
C3H6O	PROPANAL
C3H6O	PROPYLENE OXIDE
C3H6O	VINYL CARBINOL
C3H6O2	1,3-DIOXOLANE
C3H6O2	EPOXY-2,3-PROPANOL-1
C3H6O2	ETHYL FORMATE
C3H6O2	GLYCIDE
C3H6O2	GLYCIDOL
C3H6O2	METHYL ACETATE
C3H7NO	DIMETHYLFORMAMIDE
C3H7NO	DMF

Formula	Chemical name
C3H7NO2	ISO-NITROPROPANE
C3H7NO2	NITROPROPANE 2
C3H8	DIMETHYLMETHANE
C3H8	PROPANE
C3H8O	2-PROPANOL
C3H8O	DIMETHYL CARBINOL
C3H8O	IPA
C3H8O	ISOPROPANOL
C3H8O	ISOPROPYL ALCOHOL
C3H8O	PROPANOL-1
C3H8O	PROPYL ALCOHOL
C3H8O	1-PROPANOL
C3H8O2	DIMETHOXYMETHANE
C3H8O2	METHYL CELLOSOLVE
C3H8O2	METHYLAL
C3H8O3	GLYCEROL, MIST
C3H8S	1-PROPANETHIOL
C3H8S	PROPYL MERCAPTAN
C3H9N	2-AMINOPROPANE
C3H9N	ISOPROPYLAMINE
C3H9N	PROPYLAMINE
C3H9N	TMA
C3H9N	TRIMETHYLAMINE
C3H9N	2-PROPYLAMINE
C3H9NO	1-AMINOPROPANE
C3H9NO	2-AMINO 1-PROPANOL
C3H9NO	3-AMINO-1-PROPANOL
C3H9NO	PROPANAMINE

Formula	Chemical name
C4H4Cl	BETA-CHLOROPRENE
C4H4Cl	CHLOROBUTADIENE
C4H5Cl	CHLOROPRENE
C4H5O2	METHOXYCARBONYLETHYLENE
C4H5O2	METHYL PROPENOATE
C4H5O2	METHYLACRYLATE
C4H6	1,3-BUTADIENE
C4H6	DIVINYL
C4H6	ERYTHRENE
C4H6	VINYL ETHYLENE
C4H6O	2-BUTENAL
C4H6O	BETA-METHYL ACROLEIN
C4H6O	CROTONALDEHYDE
C4H6O	PROPYLENE ALDEHYDE
C4H6O2	METHACRYLIC ACID
C4H6O2	VINYL ACETATE
C4H6O3	ACETIC ANHYDRE
C4H6O3	ACETIC OXIDE
C4H8O	2-BUTANONE
C4H8O	DIETHYLENE OXIDE
C4H8O	ETHYL METHYL KETONE
C4H8O	MEK
C4H8O	METHYL ACETONE
C4H8O	METHYL ETHYL KETONE
C4H8O	TETRAHYDROFURAN
C4H8O	THF
C4H8O2	1, 4-DIOXANE
C4H8O2	BUTANOIC ACID

Formula	Chemical name
C4H8O2	BUTYRIC ACID
C4H8O2	DIETHYLENE DIOXIDE
C4H8O2	ETHYL ACETATE
C4H8O2	ETHYL ETHANOATE
C4H8OC12	2, 2'-DICHLORODIETHYL ETHER
C4H9Cl	1-CHLORO BUTANE
C4H9Cl	N-BUTYL CHLORIDE
C4H9Cl	TERT-BUTYL CHLORIDE
C4H9NH2	1-AMINOBUTANE
C4H10	ISOBUTANE
C4H10	METHYL-2-PROPANE
C4H10	N-BUTANE
C4H10	TRIMETHYL METHANE
C4H10O	2-BUTANOL
C4H10O	2-METHYL-1-PROPANOL
C4H10O	BUTYL ALCOHOL
C4H10O	BUTYL ALCOHOL SEC
C4H10O	BUTYL ALCOHOL TER
C4H10O	BUTYLENE HYDRATE
C4H10O	DIETHYL ETHER
C4H10O	DIETHYL OXIDE
C4H10O	ETHER
C4H10O	ETHYL ETHER
C4H10O	ETHYL OXIDE
C4H10O	ISOBUTANOL
C4H10O	ISOBUTYL ALCOHOL
C4H10O	ISOPROPYL CARBINOL
C4H10O	METHYL-2-PROPANOL-2

Formula	Chemical name
C4H10O	METHYLETHYL CARBINOL
C4H10O	N-BUTANOL
C4H10O	TERT-BUTYL ALCOHOL
C4H10O	TRIMETHYL CARBINOL
C4H10O	1-BUTANOL
C4H10O	SBA
C4H10O2	2-ETHOXYETHANOL
C4H10O2	CELLOSOLVE®
C4H10O2	ETHYLENE GLYCOL MONO ETHYL ETHER
C4H10S	1-MERCAPTOBUTANE
C4H10S	N-BUTANETHIOL
C4H10S	N-BUTYL MERCAPTAN
C4H10S	1-BUTANETHIOL
C4H11N	DIETHAMINE
C4H11N	DIETHYLAMINE
C4H11N	N,N-DIMETHYLETHYLAMINE
C4H11N	N-BUTYL AMINE
C4H11N	N-ETHYLETHANAMINE
C4H11N	SEC-BUTYL AMINE
C4H11N	2-AMINO BUTANE
C4H11NO2	DEA
C4H11NO2	DIETHANOLAMINE
C4H13N3	DIETHYLENE TRIAMINE
C5H5N	AZINE
C5H5N	PYRIDINE
C5H5NO2	METHYL CYANOACRYLATE
C5H6	1,3-CYCLOPENTADIENE
C5H6N2	2-AMINO PYRIDINE

Formula	Chemical name
C5H6O2	2-FURYL METHANOL
C5H6O2	2-HYDROXYMETHYLFURAN
C5H6O2	FURFURYL ALCOHOL
C5H6O2	FURYL CARBINOL
C5H8	2-METHYL-1,3-BUTADIENE
C5H8	ISOPRENE
C5H8O2	ETHYL ACRYLATE
C5H8O2	GLUTARALDEHYDE
C5H8O2	METHYL METACRYLATE
C5H9NO	1-METHYL-2-PYRROLIDINONE
C5H10	CYCLOPENTANE
C5H10O	3-PENTANONE
C5H10O	DIETHYL KETONE
C5H10O	2-PENTANONE
C5H10O	DIMETHYLACETONE
C5H10O	ETHYL KETONE
C5H10O	METHYL PROPYL KETONE
C5H10O2	ISOPROPYL ACETATE
C5H10O2	2-PROPYL ACETATE
C5H10O2	PROPYL ACETATE
C5H10O2	PROPYLACETATE
C5H12	2-METHYLBUTANE
C5H12	ISOPENTANE
C5H12	N-PENTANE
C5H12O	3-METHYL-1-BUTANOL
C5H12O	AMYL ALCOHOL N
C5H12O	2-METHOXY-2-METHYLPROPANE
C5H12O	BUTYL CARBINOL

Formula	Chemical name
C5H12O	ISOAMYL ALCOHOL
C5H12O	ISOAMYL ALCOHOL (PRIMARY)
C5H12O	ISOBUTYL CARBINOL
C5H12O	METHYL-3-BUTANOL-1
C5H12O	METHYL TERT-BUTYL ETHER
C5H12O	MTBE
C5H12O	PENTANOL 1
C5H12O	TERT-BUTYL METHYL ETHER
C6H4Cl2	1,2-DICHLOROBENZENE
C6H4Cl2	P-DICHLOROBENZENE
C6H5Cl	BENZENE CHLORIDE
C6H5Cl	CHLOROBENZENE
C6H5Cl	PHENYL CHLORIDE
C6H5NO2	ESSENCE OF MIRBANE
C6H5NO2	MIRBANE OIL
C6H5NO2	NITRO BENZENE
C6H6	BENZENE
C6H6O	HYDROXYBENZENE
C6H6O	PHENOL
C6H6O	PHENYL HYDROXIDE
C6H6O2	HYDROQUINONE
C6H7N	AMINO-BENZENE
C6H7N	ANILINE
C6H7N	BENZENAMINE
C6H7N	PHENYL AMINE
C6H7NO2	ETHYL CYANOACRYLATE
C6H10	CYCLOHEXENE
C6H10O	3-METHYL-3-PENTEN-2-ONE

Formula	Chemical name
C6H10O	CYCLOHEXANONE
C6H10O	CYCLOHEXYL KETONE
C6H10O	ISOBUTENYL METHYL KETONE
C6H10O	ISOPROPYLIDENEACETONE
C6H10O	MESITYL OXIDE
C6H10O	METHYL ISOBUTENYL KETONE
C6H10O2	ALLYLGLYCIDYLETHER
[C6H10O5] _n	CELLULOSE
[C6H10O5] _n	HYDROXYCELLULOSE
[C6H10O5] _n	PYROCELLULOSE
C6H12	CYCLOHEXANE
C6H12O	2-HEXANONE
C6H12O	4-METHYL 2-PENTANONE
C6H12O	BUTYL VINYL ETHER
C6H12O	BVE
C6H12O	CYCLOHEXANOL
C6H12O	CYCLOHEXYL ALCOHOL
C6H12O	MIBK
C6H12O	HEXONE
C6H12O	HYDROXYCYCLOHEXANE
C6H12O	METHYL BUTYL KETONE
C6H12O	METHYL ISOBUTYL KETONE
C6H12O2	1,2-EPOXY-3-ISOPROPOXYPROPANE
C6H12O2	2-METHYLPROPYL ACETATE
C6H12O2	2-METHYLPROPYL ESTER OF ACETIC ACID
C6H12O2	BETA-METHYLPROPYL ETHANOATE
C6H12O2	DIACETONE
C6H12O2	DIACETONE ALCOHOL

Formula	Chemical name
C6H12O2	ISOBUTYL ACETATE
C6H12O2	ISOPROPYL GLYCIDYL ETHER
C6H12O2	N-BUTYL ACETATE
C6H12O2	TERT-BUTYL ACETATE
C6H12O3	2-ETHOXY ACETATE
C6H12O3	CELLOSOLVE "ACETATE"
C6H12O3	ETHYLENE GLYCOL MONO ETHYL ETHER ACETATE
C6H12O6	GLUCOSE
C6H13N	AMINOCYCLOHEXANE
C6H13N	CYCLOHEXYLAMINE
C6H14	HEXANE
C6H14	N-HEXANE
C6H14O	1-hexanol
C6H14O	DIISOPROPYL ETHER
C6H14O	DIISOPROPYL OXIDE
C6H14O	HEXANOL
C6H14O	HEXYL ALCOHOL
C6H14O	ISOBUTYL METHYL CARBINOL
C6H14O	ISOPROPYL ETHER
C6H14O	METHYLAMYL ALCOHOL
C6H14O	MIBC
C6H14O	N-HEXANOL
C6H14O2	2-BUTOXYETHANOL
C6H14O2	3-METHOXY-3-METHYL-1-BUTANOL
C6H14O2	BUTYL CELLOSOLVE®
C6H14O2	BUTYL GLYCOL

Formula	Chemical name
C6H15N	DIISOPROPYLAMINE
C6H15N	TRIETHYLAMINE
C6H15NO	DIETHYLAMINOETHANOL-2
C7H7Cl	A-CHLOROTOLUENE
C7H7Cl	BENZYL CHLORIDE
C7H7Cl	CHLOROTOLUENE (ORTHO)
C7H8	METHYL BENZENE
C7H8	PHENYLMETHANE
C7H8	TOLUENE
C7H8	TOLUOL
C7H8O	3-CRESOL
C7H8O	3-HYDROXYTOLUENE
C7H8O	3-METHYL PHENOL
C7H8O	4-CRESOL
C7H8O	4-HYDROXYTOLUENE
C7H8O	BENZYL ALCOHOL
C7H8O	CRESOL ALL ISOMERS
C7H8O	M-CRESOL
C7H8O	METHYL PHENOL ALL ISOMERS
C7H8O	P-CRESOL
C7H9N	4-AMINOTOLUENE
C7H9N	4-METHYLANILINE
C7H9N	O-TOLUIDINE
C7H9N	P-TOLUIDINE
C7H12O	METHYL CYCLOHEXANONE
C7H12O2	BUTYL ACRYLATE
C7H14	METHYL CYCLOHEXANE
C7H14O	2,4-DIMETHYL-3-PENTANONE

Formula	Chemical name
C7H14O	2-HEPTANONE
C7H14O	DIISOPROPYL KETONE
C7H14O	DIPROPYL KETONE
C7H14O	HEPTAN-4-ONE
C7H14O	ISOBUTYRONE
C7H14O	METHYL CYCLOHEXANOL
C7H14O	METHYL-N-AMYL KETONE
C7H14O2	BUTYL GLYCIDYL ETHER
C7H14O2	ISOAMYL ACETATE
C7H14O2	ISOPENTYL ACETATE
C7H14O2	N-AMYL ACETATE
C7H14O2	N-AMYL ACETATE
C7H14O2	PENTYL ACETATE
C7H14O2	SEC-AMYL ACETATE
C7H14O3	BUTYL LACTATE
C7H16	2, 4-DIMETHYL PENTANE
C7H16	DIPROPYLMETHANE
C7H16	HEPTANE
C8H8	CINAMENE
C8H8	PHENYLETHYLENE
C8H8	STYRENE
C8H8	VINYLBENZENE
C8H10	DIMETHYL BENZENE (AND ISOMERS)
C8H10	ETHYL BENZENE
C8H10	O-XYLENE
C8H10	PHENYLETHANE
C8H10	XYLENE
C8H14ClN5	ATRAZINE

Formula	Chemical name
C8H16N2	1,1'-BIPHENYL-4,4'-DIAMINE
C8H16N2	4,4'-BIANILINE
C8H16N2	4,4'-BIPHENYLDIAMINE
C8H16N2	4,4'-DIAMINOBIIPHENYL
C8H16N2	BENZINE 35 80
C8H16O	3-OCTANONE
C8H16O	EAK
C8H16O	ETHYLAMYL KETONE
C8H18	ISOOCTANE
C8H18	N-OCTANE
C8H10	XYLENE (ISOMERS)
C8H18	TRIMETHYL PENTANE-2,2,4
C8H18O	2-ETHYL-1-HEXANOL
C8H18O	BUTYL ETHER
C8H18O	DIBUTYL ETHER
C8H18O	ISOOCTANOL
C8H18O	ISOOCTYLALCOHOL
C8H18O3	DIETHYLENE GLYCOL MONOBUTYL ETHER
C9H4O3 .H2O	NINHYDRIN (POWDER)
C9H8O4	ACETYSALICYCLIC ACID
C9H8O4	ASPIRIN
C9H10	VINYL TOLUENE
C9H10	METHYL STYRENE
C9H10Cl2N2O	DIURON
C9H10O2	PHENYL GLYCIDYL ETHER
C9H12	2-PHENYL PROPANE
C9H12	CUMENE
C9H12	CUMOL

Formula	Chemical name
C9H12	ISOPROPYL BENZENE
C9H12	ISOPROPYL BENZENE
C9H12	MESITYLENE
C9H12	PROPYLBENZENE
C9H12	TRIMETHYLBENZENE
C9H14O	ISOPHORONE
C9H18O	2,6-DIMETHYL-4-HEPTANONE
C9H18O	DIISOBUTYL KETONE
C9H20	NONANE ALL ISOMERS
C10H8	NAPHTHALENE
C10H8	NAPHTHALIN
C10H10	1,3-DIVINYLBENZENE
C10H12	BICYCLOPENTADIENE
C10H12	DICYCLOPENTADIENE
C10H14NO5PS	PARATHION
C10H16N2O8	EDTA
C10H16	LIMONENE
C10H16	TURPENTINE OIL
C10H22	DECANE
C11H16	4-TERT-BUTYL TOLUENE
C12H4Cl4O2	DIOXIN
C12H4Cl4O2	DIOXINE
C12H10O	DIPHENYL OXIDE
C12H10O	PHENYL ETHER
C12H14N2	PARAQUAT
C12H22O11	SACCAROSE
C18H14O2	BUTYL METACRYLATE
C20H14O4	PHENOLPHTALEINE

Formula	Chemical name
C21H14Br4O5S	VERT DE BROMOCRESOL
C21H20N3Br	BET
C21H20N3Br	ETHIDIUM BROMIDE
C21H22N2O2	STRYCHNINE
C27H39O5S	BLEU DE THYMOL
CO	CARBON MONOXIDE
CO2	CARBON DIOXIDE
CS2	CARBON DISULFIDE
CaCO3	CALCIUM CARBONATE
Ca(OH)2	CALCIUM HYDRATE
Ca(OH)2	CALCIUM HYDROXIDE
CaO	CALCIUM OXIDE
Ca(SO4). 2H2O	CALCIUM SULFATE
Cd	CADMIUM DUST (AS CD)
Cd	CADMIUM FUME (AS CD)
ClO2	CHLORINE DIOXIDE
ClO2	CHLORINE OXIDE
Cl2	CHLORINE
CrO3	CHROMIC ACID
CrO3	CHROMIC OXIDE
CrO3	CHROMIUM(VI) OXIDE (1:3)
CuO	COPPER(II) OXIDE FUME
Cu	COPPER (DUSTS AND MISTS, AS CU)
Fe2O3	FERRIC OXIDE
Fe2O3	IRON OXIDE DUST AND FUME (AS FE)
Fe2O3	RED IRON OXIDE
HBr	ANHYDROUS HYDROGEN BROMIDE
HBr	AQUEOUS HYDROGEN BROMIDE (I.E.

Formula	Chemical name
HBr	HYDROBROMIC ACID
HBr	HYDROGEN BROMIDE
HCN	FORMONITRILE
HCN	HYDROCYANIC ACID
HCN	HYDROGEN CYANIDE
HCN	PRUSSIC ACID
HClO	HYPOCHLOROUS ACID
HClO4	PERCHLORIC ACID
HCl+HNO3	AQUA REGIA
HCl	HYDROGEN CHLORIDE
HCl aq. sol.	HYDROCHLORIC ACID
HCl aq. sol.	MURIATIC ACID
HCl aq. sol.	AQUEOUS HYDROGEN CHLORIDE (I.E.
HF	HYDROGEN FLUORIDE
HF aq. sol.	HYDROFLUORIC ACID
HNO3	AQUA FORTIS
HNO3	HYDROGEN NITRATE
HNO3	NITRIC ACID
HNaO3S	SODIUM BISULPHITE
H2O2	HYDROGEN DIOXIDE
H2O2	HYDROGEN PEROXIDE
H2S	HYDROGEN SULFIDE
H2SO4	HYDROGEN SULFATE
H2SO4	SULFURIC ACID
H3PO4	ORTHOPHOSPHORIC ACID
H3PO4	PHOSPHORIC ACID

Formula	Chemical name
Hg	MERCURY
Hg	METALLIC MERCURY
I2	IODINE
KOH	CAUSTIC POTASH
KOH	POTASSIUM HYDRATE
KOH	POTASSIUM HYDROXIDE
K(Mg, Fe)3AlSi3O10(F, OH)2	MICA (CONTAINING LESS THAN 1% QUARTZ)
K(Mg, Fe)3AlSi3O10(F, OH)2	MUSCOVITE
K(Mg, Fe)3AlSi3O10(F, OH)2	BIOTITE
LiH	LITHIUM HYDRIDE
MgCO3	MAGNESITE
MgO	MAGNESIA FUME
MgO	MAGNESIUM OXIDE FUME
Mn	MANGANESE COMPOUNDS (AS MN)
MnO2	MANGANESE OXIDE
NH3	AMMONIA
NH4Cl	AMMONIUM CHLORIDE
NH4Cl	AMMONIUM CHLORIDE FUME
NH4OH	AMMONIUM HYDROXYDE SOL
N2H4	HYDRAZINE
N2H4	DIAMINE
N2O	DINITROGEN MONOXIDE
NaCl	CHLORURE DE SODIUM
NaOH	CAUSTIC SODA

Formula	Chemical name
NaOH	SODIUM HYDRATE
NaOH	SODIUM HYDROXIDE
Na2B4O7 • 10H2O	BORAX
NO2	NITROGEN DIOXIDE
Ni	NICKEL METAL AND OTHER COMPOUNDS (AS NI)
O2Sn	TIN(IV) OXIDE (AS SN)
O3	OZONE
OsO4	OSMIUM TETROXIDE (AS OS)
Pt	PLATINUM
Si	SILICON
SiO2	QUARTZ
SiO2	SILICA GEL
SiO2	SILICA, AMORPHOUS
Sn	TIN (INORGANIC COMPOUNDS, AS SN)
TiO2	TITANIUM DIOXIDE
W	TUNGSTEN insoluble
ZnO	ZINC OXIDE
85% Nonane/15% trimethylbenzene	NAPHTA 30/60
85% Nonane/15% trimethylbenzene	STODDARD SOLVENT

Formula	Chemical name
85% Nonane/15% trimethylbenzene	WHITE SPIRIT

List of chemicals substances by CAS number

CAS number	Chemical name
50-00-0	FORMALDEHYDE
50-00-0	FORMALDEHYDE SOLUTION
50-00-0	FORMALIN (AS FORMALDEHYDE)
50-00-0	FORMIC ALDEHYDE
50-00-0	METHANAL
50-00-0	METHYL ALDEHYDE
50-00-0	METHYLENE OXIDE
50-78-2	ACETYLSALICYCLIC ACID
50-78-2	ASPIRIN
55-63-0	NITROGLYCERINE
55-63-0	TRINITROGLYCERINE
56-23-5	CARBON TETRACHLORIDE
56-23-5	TETRACHLOROMETHANE
56-38-2	PARATHION
56-81-5	GLYCEROL , MIST
57-13-6	UREA
57-24-9	STRYCHNINE
57-50-1	SACCAROSE
60-00-4	EDTA
60-24-2	MERCAPTO-2 ETHANOL
60-29-7	DIETHYL ETHER
60-29-7	DIETHYL OXIDE
60-29-7	ETHER
60-29-7	ETHYL ETHER

CAS number	Chemical name
60-29-7	ETHYL OXIDE
62-53-3	AMINO-BENZENE
62-53-3	ANILINE
62-53-3	BENZENAMINE
62-53-3	PHENYL AMINE
64-17-5	ABSOLUTE ALCOHOL
64-17-5	ALCOHOL
64-17-5	ETHANOL
64-17-5	ETHYL ALCOHOL
64-18-6	FORMIC ACID
64-18-6	METHANOIC ACID
64-19-7	ACETIC ACID
64-19-7	ETHANOIC ACID
64-19-7	GLACIAL ACETIC ACID (PURE COMPOUND)
67-56-1	METHANOL
67-56-1	METHYL ALCOHOL
67-63-0	2-PROPANOL
67-63-0	DIMETHYL CARBINOL
67-63-0	IPA
67-63-0	ISOPROPANOL
67-63-0	ISOPROPYL ALCOHOL
67-64-1	2-PROPANONE
67-64-1	ACETONE
67-64-1	DIMETHYL KETONE

CAS number	Chemical name
67-66-3	CHLOROFORM
67-66-3	TRICHLOROMETHANE
67-68-5	DIMETHYL SULFOXIDE
67-68-5	DMSO
68-12-2	DIMETHYLFORMAMIDE
68-12-2	DMF
71-23-8	1-PROPANOL
71-23-8	PROPANOL-1
71-23-8	PROPYL ALCOHOL
71-36-3	1-BUTANOL
71-36-3	BUTYL ALCOHOL
71-36-3	N-BUTANOL
71-41-0	AMYL ALCOHOL N
71-41-0	BUTYL CARBINOL
71-41-0	ISOAMYL ALCOHOL
71-41-0	PENTANOL 1
71-43-2	BENZENE
71-55-6	1,1,1-TRICHLOROETHANE
71-55-6	CHLOROTHENE
71-55-6	METHYL CHLOROFORM
74-86-2	ACETYLENE
74-86-2	ETHENE
74-86-2	ETHYNE
74-87-3	CHLOROMETHANE
74-87-3	METHYL CHLORIDE
74-89-5	AMINOMETHANE
74-89-5	METHYLAMINE
74-90-8	FORMONITRILE

CAS number	Chemical name
74-90-8	HYDROCYANIC ACID
74-90-8	HYDROGEN CYANIDE
74-90-8	PRUSSIC ACID
74-96-4	BROMOETHANE
74-96-4	ETHYL BROMIDE
74-97-5	BROMOCHLOROMETHANE
74-97-5	CHLOROBROMOMETHANE
74-97-5	METHYLENE CHLOROBROMIDE
74-98-6	DIMETHYLMETHANE
74-98-6	PROPANE
74-99-7	ALLYLENE
74-99-7	METHYL ACETYLENE
74-99-7	PROPYNE
75-00-3	CHLOROETHANE
75-00-3	ETHYL CHLORIDE
75-01-4	CHLOROETHENE
75-01-4	CHLOROETHYLENE
75-01-4	VC
75-01-4	VINYL CHLORIDE
75-04-7	AMINOETHANE
75-04-7	ETHYLAMINE
75-05-8	ACETONITRILE
75-05-8	CYANOMETHANE
75-05-8	ETHYL NITRILE
75-05-8	METHYL CYANIDE
75-07-0	ACETALDEHYDE
75-07-0	ETHANAL
75-07-0	ETHYL ALDEHYDE

CAS number	Chemical name
75-09-2	DICHLOROMETHANE
75-09-2	METHYLENE CHLORIDE
75-09-2	METHYLENE DICHLORIDE
75-12-7	FORMAMIDE
75-15-0	CARBON DISULFIDE
75-25-2	BROMOFORM
75-25-2	TRIBROMOMETHANE
75-28-5	ISOBUTANE
75-28-5	METHYL-2-PROPANE
75-28-5	TRIMETHYL METHANE
75-31-0	2-AMINOPROPANE
75-31-0	2-PROPYLAMINE
75-31-0	ISOPROPYLAMINE
75-34-3	1,1-DICHLOROETHANE
75-34-3	ETHYLIDENE CHLORIDE
75-47-8	IODOFORM
75-50-3	TMA
75-50-3	TRIMETHYLAMINE
75-52-5	NITROMETHANE
75-56-9	PROPYLENE OXIDE
75-65-0	BUTYL ALCOHOL TER
75-65-0	METHYL-2-PROPANOL-2
75-65-0	TERT-BUTYL ALCOHOL
75-65-0	TRIMETHYL CARBINOL
76-01-7	PENTACHLOROETHANE
76-03-9	TRICHLOROACETIC ACID
76-03-9	TRICHLOROETHANOIC ACID
76-05-1	TFA

CAS number	Chemical name
76-05-1	TRIFLUOROACETIC ACID
76-60-8	VERT DE BROMOCRESOL
76-61-9	BLEU DE THYMOL
77-09-8	PHENOLPHTALEINE
77-73-6	BICYCLOPENTADIENE
77-73-6	DICYCLOPENTADIENE
78-59-1	ISOPHORONE
78-78-4	2-METHYLBUTANE
78-78-4	ISOPENTANE
78-79-5	2-METHYL-1,3-BUTADIENE
78-79-5	ISOPRENE
78-83-1	2-METHYL-1-PROPANOL
78-83-1	ISOBUTANOL
78-83-1	ISOBUTYL ALCOHOL
78-83-1	ISOPROPYL CARBINOL
78-87-5	DICHLOROPROPANE 1, 2
78-87-5	PROPYLENE CHLORIDE
78-87-5	PROPYLENE DICHLORIDE
78-92-2	2-BUTANOL
78-92-2	BUTYL ALCOHOL SEC
78-92-2	BUTYLENE HYDRATE
78-92-2	METHYLETHYL CARBINOL
78-92-2	SBA
78-93-3	2-BUTANONE
78-93-3	ETHYL METHYL KETONE
78-93-3	MEK
78-93-3	METHYL ACETONE
78-93-3	METHYL ETHYL KETONE

CAS number	Chemical name
79-00-5	TRICHLOROETHANE-1,1,2
79-00-5	VINYL TRICHLORIDE
79-01-6	ETHYLENE TRICHLORIDE
79-01-6	TRICHLOROETHENE
79-01-6	TRICHLOROETHYLENE
79-06-1	2-PROPENAMIDE
79-06-1	ACRYLAMIDE
79-10-7	2-PROPENOIC ACID
79-10-7	ACROLEIC ACID
79-10-7	ACRYLIC ACID
79-20-9	METHYL ACETATE
79-24-3	NITROETHANE
79-27-6	1,1,2,2-TETRABROMOETHANE
79-27-6	ACETYLENE TETRABROMIDE
79-34-5	1,1,2,2-TETRACHLOROETHANE
79-34-5	ACETYLENE TETRACHLORIDE
79-41-4	METHACRYLIC ACID
79-46-9	ISO-NITROPROPANE
79-46-9	NITROPROPANE 2
80-62-6	METHYL METACRYLATE
91-20-3	NAPHTHALENE
91-20-3	NAPHTHALIN
92-87-5	1,1'-BIPHENYL-4,4'-DIAMINE
92-87-5	4,4'-BIANILINE
92-87-5	4,4'-BIPHENYLDIAMINE
92-87-5	4,4'-DIAMINOBIIPHENYL
92-87-5	BENZINE 35 80
95-47-6	DIMETHYL BENZENE (AND ISOMERS)

CAS number	Chemical name
95-47-6	O-XYLENE
95-50-1	1,2-DICHLOROBENZENE
96-22-0	3-PENTANONE
96-22-0	DIETHYL KETONE
96-22-0	DIMETHYLACETONE
96-22-0	ETHYL KETONE
96-33-3	METHOXYCARBONYLETHYLENE
96-33-3	METHYL PROPENOATE
96-33-3	METHYLACRYLATE
97-88-1	BUTYL METACRYLATE
98-00-0	2-FURYL METHANOL
98-00-0	2-HYDROXYMETHYLFURAN
98-00-0	FURFURYL ALCOHOL
98-00-0	FURYL CARBINOL
98-51-1	4-TERT-BUTYL TOLUENE
98-82-8	2-PHENYL PROPANE
98-82-8	CUMENE
98-82-8	CUMOL
98-82-8	ISOPROPYL BENZENE
98-82-8	ISOPROPYL BENZENE
98-95-3	ESSENCE OF MIRBANE
98-95-3	MIRBANE OIL
98-95-3	NITRO BENZENE
100-37-8	DIETHYLAMINOETHANOL-2
100-41-4	ETHYL BENZENE
100-41-4	PHENYLETHANE
100-42-5	CINAMENE
100-42-5	PHENYLETHYLENE

CAS number	Chemical name
100-42-5	STYRENE
100-42-5	VINYLBENZENE
100-44-7	A-CHLOROTOLUENE
100-44-7	BENZYL CHLORIDE
100-51-6	BENZYL ALCOHOL
101-84-8	DIPHENYL OXIDE
101-84-8	PHENYL ETHER
104-76-7	2-ETHYL-1-HEXANOL
104-76-7	ISOOCTANOL
104-76-7	ISOOCTYLALCOHOL
106-43-4	CHLOROTOLUENE (ORTHO)
106-44-5	4-CRESOL
106-44-5	4-HYDROXYTOLUENE
106-44-5	P-CRESOL
106-46-7	P-DICHLOROBENZENE
106-49-0	4-AMINOTOLUENE
106-49-0	4-METHYLANILINE
106-49-0	P-TOLUIDINE
106-68-3	3-OCTANONE
106-68-3	EAK
106-68-3	ETHYLAMYL KETONE
106-89-8	1-CHLORO-2,3-EPOXYPROPANE
106-89-8	2-CHLOROPROPYLENE OXIDE
106-89-8	EPICHLORHYDRINE
106-92-3	ALLYLGLYCIDYLETHER
106-93-4	1,2-DIBROMOETHANE
106-93-4	ETHYLENE BROMIDE
106-93-4	ETHYLENE DIBROMIDE

CAS number	Chemical name
106-97-8	N-BUTANE
106-99-0	1,3-BUTADIENE
106-99-0	DIVINYL
106-99-0	ERYTHRENE
106-99-0	VINYL ETHYLENE
107-02-8	2-PROPENAL
107-02-8	ACROLEIN
107-02-8	ACRYLIC ALDEHYDE
107-02-8	ALLYL ALDEHYDE
107-03-9	1-PROPANETHIOL
107-03-9	PROPYL MERCAPTAN
107-05-1	3-CHLORO-1-PROPENE
107-05-1	ALLYL CHLORIDE
107-06-2	1,2-DICHLOROETHANE
107-06-2	ETHYLENE CHLORIDE
107-06-2	ETHYLENE DICHLORIDE
107-07-3	2-CHLOROETHANOL
107-07-3	2-CHLOROETHYL ALCOHOL
107-07-3	ETHYLEN CHLORHYDRIN
107-07-3	ETHYLENE CHLOROXYDRIN
107-10-8	1-AMINOPROPANE
107-10-8	PROPANAMINE
107-10-8	PROPYLAMINE
107-13-1	2-PROPENENITRILE
107-13-1	ACRYLONITRILE
107-13-1	PROPENENITRILE
107-13-1	VINYL CYANIDE
107-15-3	ETHYLENE DIAMINE (SOLUTION)

CAS number	Chemical name
107-15-3	ETHYLENEDIAMINE
107-18-6	2-PROPEN-1-OL
107-18-6	2-PROPENOL
107-18-6	ALLYL ALCOHOL
107-18-6	ALLYLIC ALCOHOL
107-18-6	VINYL CARBINOL
107-19-7	2-PROPYN-1-OL
107-19-7	2-PROPYNYL ALCOHOL
107-19-7	PROPARGYL ALCOHOL
107-20-0	2-CHLOROACETALDEHYDE
107-20-0	2-CHLOROETHANAL
107-21-1	1,2-ETHANEDIOL
107-21-1	ETHYLENE ALCOHOL
107-21-1	ETHYLENE GLYCOL
107-21-1	GLYCOL
107-31-3	METHYL FORMATE
107-87-9	2-PENTANONE
107-87-9	METHYL PROPYL KETONE
107-92-6	BUTANOIC ACID
107-92-6	BUTYRIC ACID
108-05-4	VINYL ACETATE
108-08-7	2, 4-DIMETHYL PENTANE
108-10-1	4-METHYL 2-PENTANONE
108-10-1	HEXONE
108-10-1	METHYL ISOBUTYL KETONE
108-10-1	MIBK
108-11-2	ISOBUTYL METHYL CARBINOL
108-11-2	METHYLAMYL ALCOHOL

CAS number	Chemical name
108-11-2	MIBC
108-18-9	DIISOPROPYLAMINE
108-20-3	DIISOPROPYL ETHER
108-20-3	DIISOPROPYL OXIDE
108-20-3	ISOPROPYL ETHER
108-21-4	2-PROPYL ACETATE
108-21-4	ISOPROPYL ACETATE
108-24-7	ACETIC ANHYDRE
108-24-7	ACETIC OXIDE
108-39-4	3-CRESOL
108-39-4	3-HYDROXYTOLUENE
108-39-4	3-METHYL PHENOL
108-39-4	M-CRESOL
108-67-8	MESITYLENE
108-67-8	PROPYLBENZENE
108-67-8	TRIMETHYLBENZENE
108-83-8	2,6-DIMETHYL-4-HEPTANONE
108-83-8	DIISOBUTYL KETONE
108-87-2	METHYL CYCLOHEXANE
108-88-3	METHYL BENZENE
108-88-3	PHENYLMETHANE
108-88-3	TOLUENE
108-88-3	TOLUOL
108-90-7	BENZENE CHLORIDE
108-90-7	CHLOROBENZENE
108-90-7	PHENYL CHLORIDE
108-91-8	AMINOCYCLOHEXANE
108-91-8	CYCLOHEXYLAMINE

CAS number	Chemical name
108-93-0	CYCLOHEXANOL
108-93-0	CYCLOHEXYL ALCOHOL
108-93-0	HYDROXYCYCLOHEXANE
108-94-1	CYCLOHEXANONE
108-94-1	CYCLOHEXYL KETONE
108-95-2	HYDROXYBENZENE
108-95-2	PHENOL
108-95-2	PHENYL HYDROXIDE
109-60-4	PROPYL ACETATE
109-60-4	PROPYLACETATE
109-66-0	N-PENTANE
109-69-3	1-CHLORO BUTANE
109-69-3	N-BUTYL CHLORIDE
109-73-9	1-AMINOBUTANE
109-73-9	N-BUTYL AMINE
109-79-5	1-BUTANETHIOL
109-79-5	1-BUTHANETHIOL
109-79-5	1-MERCAPTOBUTANE
109-79-5	N-BUTANETHIOL
109-79-5	N-BUTYL MERCAPTAN
109-86-4	METHYL CELLOSOLVE
109-87-5	DIMETHOXYMETHANE
109-87-5	METHYLAL
109-89-7	DIETHAMINE
109-89-7	DIETHYLAMINE
109-89-7	N-ETHYLETHANAMINE
109-94-4	ETHYL FORMATE
109-99-9	DIETHYLENE OXIDE

CAS number	Chemical name
109-99-9	TETRAHYDROFURAN
109-99-9	THF
110-19-0	2-METHYLPROPYL ACETATE
110-19-0	2-METHYLPROPYL ESTER OF ACETIC ACID
110-19-0	BETA-METHYLPROPYL ETHANOATE
110-19-0	ISOBUTYL ACETATE
110-43-0	2-HEPTANONE
110-43-0	METHYL-N-AMYL KETONE
110-54-3	HEXANE
110-54-3	N-HEXANE
110-80-5	2-ETHOXYETHANOL
110-80-5	CELLOSOLVE®
110-80-5	ETHYLENE GLYCOL MONO ETHYL ETHER
110-82-7	CYCLOHEXANE
110-83-8	CYCLOHEXENE
110-86-1	AZINE
110-86-1	PYRIDINE
111-15-9	2-ETHOXY ACETATE
111-15-9	CELLOSOLVE "ACETATE"
111-15-9	ETHYLENE GLYCOL MONO ETHYL ETHER ACETATE
111-27-3	1-hexanol
111-27-3	HEXANOL
111-27-3	HEXYL ALCOHOL
111-27-3	N-HEXANOL
111-30-8	GLUTARALDEHYDE
111-34-2	BUTYL VINYL ETHER

CAS number	Chemical name
111-34-2	BVE
111-40-0	DIETHYLENE TRIAMINE
111-42-2	DEA
111-42-2	DIETHANOLAMINE
111-44-4	2, 2'-DICHLORODIETHYL ETHER
111-65-9	N-OCTANE
111-76-2	2-BUTOXYETHANOL
111-76-2	BUTYL CELLOSOLVE®
111-76-2	BUTYL GLYCOL
111-84-2	NONANE ALL ISOMERS
112-34-5	DIETHYLENE GLYCOL MONOBUTYL ETHER
115-10-6	DIMETHYL ETHER
115-10-6	METHYL ETHER
119-93-7	O-TOLUIDINE
121-44-8	TRIETHYLAMINE
122-60-1	PHENYL GLYCIDYL ETHER
123-19-3	DIPROPYL KETONE
123-19-3	HEPTAN-4-ONE
123-31-9	HYDROQUINONE
123-38-6	PROPANAL
123-42-2	DIACETONE
123-42-2	DIACETONE ALCOHOL
123-51-3	3-METHYL-1-BUTANOL
123-51-3	ISOAMYL ALCOHOL (PRIMARY)
123-51-3	ISOBUTYL CARBINOL
123-86-4	N-BUTYL ACETATE
123-91-1	1, 4-DIOXANE
123-91-1	DIETHYLENE DIOXIDE

CAS number	Chemical name
123-92-2	ISOAMYL ACETATE
123-92-2	ISOPENTYL ACETATE
123-92-2	N-AMYL ACETATE
123-92-2	SEC-AMYL ACETATE
124-18-5	DECANE
124-38-9	CARBON DIOXIDE
124-40-3	DIMETHYL AMINE
124-40-3	DMA
126-99-8	BETA-CHLOROPRENE
126-99-8	CHLOROBUTADIENE
126-99-8	CHLOROPRENE
127-18-4	PERCHLOROETHYLENE
127-18-4	TETRACHLOROETHYLENE
137-05-3	METHYL CYANOACRYLATE
138-22-7	BUTYL LACTATE
140-88-5	ETHYL ACRYLATE
141-32-2	BUTYL ACRYLATE
141-43-5	2-AMINOETHANOL
141-43-5	BETA-AMINOETHYL ALCOHOL
141-43-5	ETHANOLAMINE
141-78-6	ETHYL ACETATE
141-78-6	ETHYL ETHANOATE
141-79-7	ISOBUTENYL METHYL KETONE
141-79-7	ISOPROPYLIDENEACETONE
141-79-7	MESITYL OXIDE
141-79-7	METHYL ISOBUTENYL KETONE
142-82-5	DIPROPYLMETHANE
142-82-5	HEPTANE

CAS number	Chemical name
142-96-1	BUTYL ETHER
142-96-1	DIBUTYL ETHER
144-62-7	ETHANEDIOIC ACID
144-62-7	OXALIC ACID
151-67-7	HALOTHANE
156-87-6	3-AMINO-1-PROPANOL
287-92-3	CYCLOPENTANE
302-01-2	DIAMINE
302-01-2	HYDRAZINE
330-54-1	DIURON
485-47-2	NINHYDRIN (POWDER)
504-29-0	2-AMINO PYRIDINE
507-20-0	TERT-BUTYL CHLORIDE
540-59-0	1,2-DICHLOROETHYLENE
540-59-0	ACETYLENE DICHLORIDE
540-84-1	ISOOCTANE
540-84-1	TRIMETHYL PENTANE-2,2,4
540-88-5	TERT-BUTYL ACETATE
541-85-5	5-METHYL-3-HEPTANONE
542-75-6	1,3-DICHLOROPROPENE
542-75-6	1,3-DICHLOROPROPYLENE
542-92-7	1,3-CYCLOPENTADIENE
546-93-0	MAGNESITE
556-52-5	EPOXY-2,3-PROPANOL-1
556-52-5	GLYCIDE
556-52-5	GLYCIDOL
558-13-4	CARBON BROMIDE
558-13-4	CARBON TETRABROMIDE

CAS number	Chemical name
558-13-4	METHANE TETRABROMIDE
558-13-4	TETRABROMOMETHANE
565-62-8	3-METHYL-3-PENTEN-2-ONE
565-80-0	2,4-DIMETHYL-3-PENTANONE
565-80-0	DIISOPROPYL KETONE
565-80-0	ISOBUTYRONE
591-78-6	2-HEXANONE
591-78-6	METHYL BUTYL KETONE
593-60-2	BROMOETHENE
593-60-2	BROMOETHYLENE
593-60-2	VINYL BROMIDE
598-56-1	N,N-DIMETHYLETHYLAMINE
600-25-9	CHLORO-1-NITROPROPANE 1
600-25-9	KORAX
628-63-7	N-AMYL ACETATE
628-63-7	PENTYL ACETATE
630-08-0	CARBON MONOXIDE
646-06-0	1,3-DIOXOLANE
872-50-4	1-METHYL-2-PYRROLIDINONE
1239-45-8	BET
1239-45-8	ETHIDIUM BROMIDE
1303-86-2	BORON OXIDE
1303-96-4	BORAX
1305-62-0	CALCIUM HYDRATE
1305-62-0	CALCIUM HYDROXIDE
1305-78-8	CALCIUM OXIDE
1309-37-1	FERRIC OXIDE
1309-37-1	IRON OXIDE DUST AND FUME (AS FE)

CAS number	Chemical name
1309-37-1	RED IRON OXIDE
1309-48-4	MAGNESIA FUME
1309-48-4	MAGNESIUM OXIDE FUME
1310-58-3	CAUSTIC POTASH
1310-58-3	POTASSIUM HYDRATE
1310-58-3	POTASSIUM HYDROXIDE
1310-73-2	CAUSTIC SODA
1310-73-2	SODIUM HYDRATE
1310-73-2	SODIUM HYDROXIDE
1314-13-2	ZINC OXIDE
1317-38-0	COPPER(II) OXIDE FUME
1317-65-3	CALCIUM CARBONATE
1319-77-3	CRESOL ALL ISOMERS
1319-77-3	METHYL PHENOL ALL ISOMERS
1321-74-0	1,3-DIVINYLBENZENE
1330-20-7	XYLENE
1330-20-7	XYLENE (ISOMERS)
1331-22-2	METHYL CYCLOHEXANONE
1332-21-4	ASBESTOS
1332-58-7	CLAY
1333-82-0	CHROMIC ACID
1333-82-0	CHROMIC OXIDE
1333-82-0	CHROMIUM(VI) OXIDE (1:3)
1333-86-4	CARBON BLACK
1344-28-1	ALUMINA
1344-28-1	ALUMINUM OXIDE
1344-28-1	ALUMINUM TRIOXIDE
1344-43-0	MANGANESE OXIDE

CAS number	Chemical name
1634-04-4	2-METHOXY-2-METHYLPROPANE
1634-04-4	METHYL TERT-BUTYL ETHER
1634-04-4	MTBE
1634-04-4	TERT-BUTYL METHYL ETHER
1746-01-6	DIOXIN
1746-01-6	DIOXINE
1912-24-9	ATRAZINE
2426-08-6	BUTYL GLYCIDYL ETHER
4016-14-2	1,2-EPOXY-3-ISOPROPOXYPROPANE
4016-14-2	ISOPROPYL GLYCIDYL ETHER
4170-30-3	2-BUTENAL
4170-30-3	BETA-METHYL ACROLEIN
4170-30-3	CROTONALDEHYDE
4170-30-3	PROPYLENE ALDEHYDE
4685-14-7	PARAQUAT
5989-54-8	LIMONENE
5996-10-1	GLUCOSE
7085-85-0	ETHYL CYANOACRYLATE
7429-90-5	ALUMINIUM
7439-96-5	MANGANESE COMPOUNDS (AS MN)
7439-97-6	MERCURY
7439-97-6	METALLIC MERCURY
7440-02-0	NICKEL METAL AND OTHER COMPOUNDS (AS NI)
7440-06-4	PLATINUM
7440-21-3	SILICON
7440-22-4	SILVER (DUST)

CAS number	Chemical name
7440-31-5	TIN (INORGANIC COMPOUNDS, AS SN)
7440-33-7	TUNGSTEN insoluble
7440-38-2	ARSENIC (INORGANIC COMPOUNDS, AS AS)
7440-41-7	BERYLLIUM COMPOUNDS (AS BE)
7440-43-9	CADMIUM DUST (AS CD)
7440-43-9	CADMIUM FUME (AS CD)
7440-44-0	GRAPHITE (SYNTHETIC)
7440-50-8	COPPER (DUSTS AND MISTS, AS CU)
7553-56-2	IODINE
7580-67-8	LITHIUM HYDRIDE
7601-90-3	PERCHLORIC ACID
7631-86-9	SILICA GEL
7631-86-9	SILICA, AMORPHOUS
7631-90-5	SODIUM BISULPHITE
7637-07-2	BORON TRIFLUORIDE
7647-01-0	AQUEOUS HYDROGEN CHLORIDE (I.E.
7647-01-0	HYDROCHLORIC ACID
7647-01-0	HYDROGEN CHLORIDE
7647-01-0	MURIATIC ACID
7647-14-5	CHLORURE DE SODIUM
7664-38-2	ORTHOPHOSPHORIC ACID
7664-38-2	PHOSPHORIC ACID
7664-39-3	HYDROFLUORIC ACID
7664-39-3	HYDROGEN FLUORIDE
7664-41-7	AMMONIA
7664-41-7	AMMONIUM HYDROXYDE SOL
7664-93-9	HYDROGEN SULFATE
7664-93-9	SULFURIC ACID

CAS number	Chemical name
7697-37-2	AQUA FORTIS
7697-37-2	HYDROGEN NITRATE
7697-37-2	NITRIC ACID
7722-84-1	HYDROGEN DIOXIDE
7722-84-1	HYDROGEN PEROXIDE
7726-95-6	BROMINE
7778-18-9	CALCIUM SULFATE
7782-50-5	CHLORINE
7783-06-4	HYDROGEN SULFIDE
7790-92-3	HYPOCHLOROUS ACID
8006-61-9	GASOLINE 60
8006-64-2	TURPENTINE OIL
8032-32-4	PETROLEUM ETHER 30/60
8052-41-3	NAPHTA 30/60
8052-41-3	STODDARD SOLVENT
8052-41-3	WHITE SPIRIT
9004-34-6	CELLULOSE
9004-34-6	HYDROXYCELLULOSE
9004-34-6	PYROCELLULOSE
10024-97-2	DINITROGEN MONOXIDE
10028-15-6	OZONE
10035-10-6	ANHYDROUS HYDROGEN BROMIDE
10035-10-6	AQUEOUS HYDROGEN BROMIDE (I.E.
10035-10-6	HYDROBROMIC ACID
10035-10-6	HYDROGEN BROMIDE
10049-04-4	CHLORINE DIOXIDE
10049-04-4	CHLORINE OXIDE
10102-44-0	NITROGEN DIOXIDE

CAS number	Chemical name
10326-38-9	BARIUM CHLORIDE
12001-26-2	BIOTITE
12001-26-2	MICA (CONTAINING LESS THAN 1% QUARTZ)
12001-26-2	MUSCOVITE
12125-02-9	AMMONIUM CHLORIDE
12125-02-9	AMMONIUM CHLORIDE FUME
13463-67-7	TITANIUM DIOXIDE
13952-84-6	2-AMINO BUTANE
13952-84-6	SEC-BUTYL AMINE
14808-60-7	QUARTZ
18282-10-5	TIN(IV) OXIDE (AS SN)
20816-12-0	OSMIUM TETROXIDE (AS OS)
25013-15-4	METHYL STYRENE
25013-15-4	VINYL TOLUENE
25639-42-3	METHYL CYCLOHEXANOL
26675-46-7	ISOFLURANE
35320-23-1	2-AMINO 1-PROPANOL
56539-66-3	3-METHOXY-3-METHYL-1-BUTANOL
-	AQUA REGIA
-	FORENE



This document is not contractual. Erlab reserves the right to modify the present document without prior notice. CHEMICAL LISTING_GFH_2021_EN

France
+33 (0) 2 32 09 55 80 | ventes@erlab.net

United States
+1 800-964-4434 | captairsales@erlab.com

China
+86 (0) 512 5781 4085 | sales.china@erlab.com.cn

Spain
+34 936 732 474 | export.south@erlab.net

Germany
0800 330 47 31 | export.north@erlab.net

United Kingdom
+44 (0) 1722 341 940 | export.north@erlab.net

Italy
+39 (0) 2 89 00 771 | export.south@erlab.net



www.erlab.com

ecosystem